



Jun 27, 2019

## Copy of Fluorescence analysis using CF imager-v2

DOI

[dx.doi.org/10.17504/protocols.io.4t9gwr6](https://dx.doi.org/10.17504/protocols.io.4t9gwr6)

Steven J Burgess<sup>1</sup>

<sup>1</sup>University of Illinois at Urbana-Champaign



Steven J Burgess

University of Illinois at Urbana-Champaign

### Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.4t9gwr6>

**Protocol Citation:** Steven J Burgess 2019. Copy of Fluorescence analysis using CF imager-v2. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.4t9gwr6>

**License:** This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** In development

We are still developing and optimizing this protocol

**Created:** June 27, 2019

**Last Modified:** June 27, 2019

**Protocol Integer ID:** 25185

**Keywords:** copy of fluorescence analysis, fluorescence analysis, fluorescence, using cf

## Troubleshooting

## Preparing plants

- 1 Dark adapt plants for at least 20 minutes prior to taking measurements.

### Note


This is done to ensure the photosynthetic electron transport chain is fully oxidized and reaction centres are open. In an ideal situation plants are allowed to dark adapt overnight prior to measurement.

A properly adapted, healthy plant should give a Fv/Fm value of ~0.8. This has been shown to be highly stable between species. Significant deviation from this (e.g. <0.7) either suggests incomplete adaptation or stressed plant material.

Prepare an extra plate with few leaf disks(or have extra plants/leaf disks) to adjust the focus (on step 5). During the focus adjustment, the light flashes inside the chamber, so you should not use the plates that are with your dark-adapted samples.

- 2 Turn on the cf imager and open the FluorImager software.



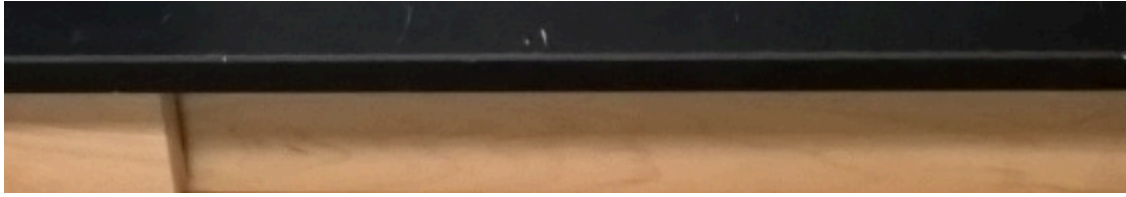

3

Start the FluorImager software



- 4 The surface of the leaf should be 140mm from the base of the imaging chamber, and can be adjusted by lowering or raising the plant under analysis. Position plant/leaf in the chamber



## Set focus

- 5 Before initiating this step, place the extra plate/plants inside the chamber.

Set the focus by adjusting the dial above the chamber, and lock in position by turning the screw on the side.

## Set exposure

- 6 Manually adjust the aperture as shown on the right to allow an optimal amount of light into the imager so as not to overexpose measurements.


After adjusting the focus and press the "map image" button, change the plate/plant to the one you are going to read. If you don't press "map image", the light will be still flashing when you open the chamber.




## 7 Isolate the plant or leaf of interest

Applying isolation background - it can be done more than once if the image is still with some background

right click  
on image

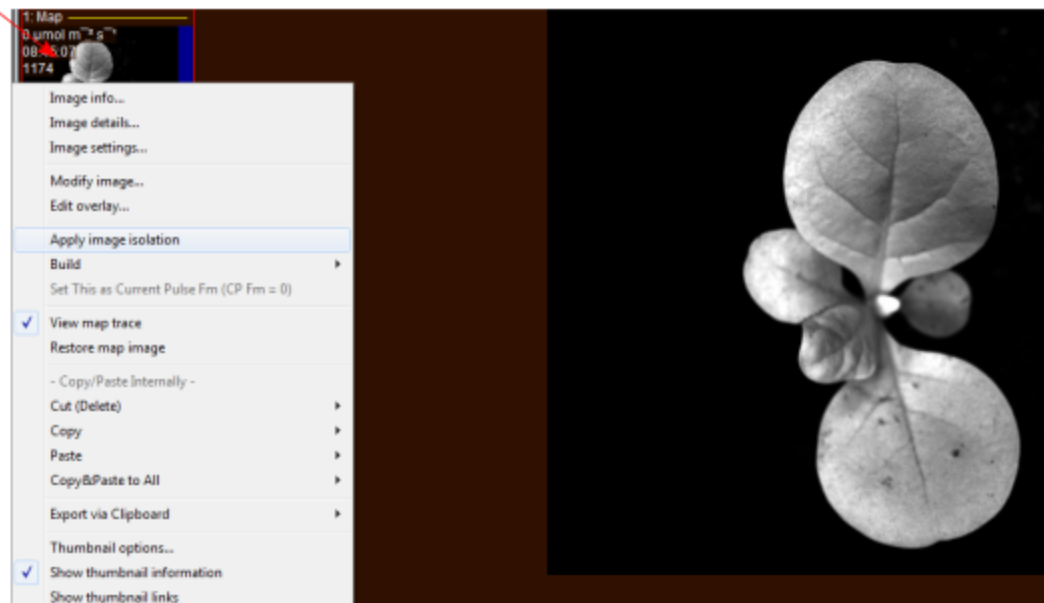


After applying isolation background areas will be masked out in blue as shown below.

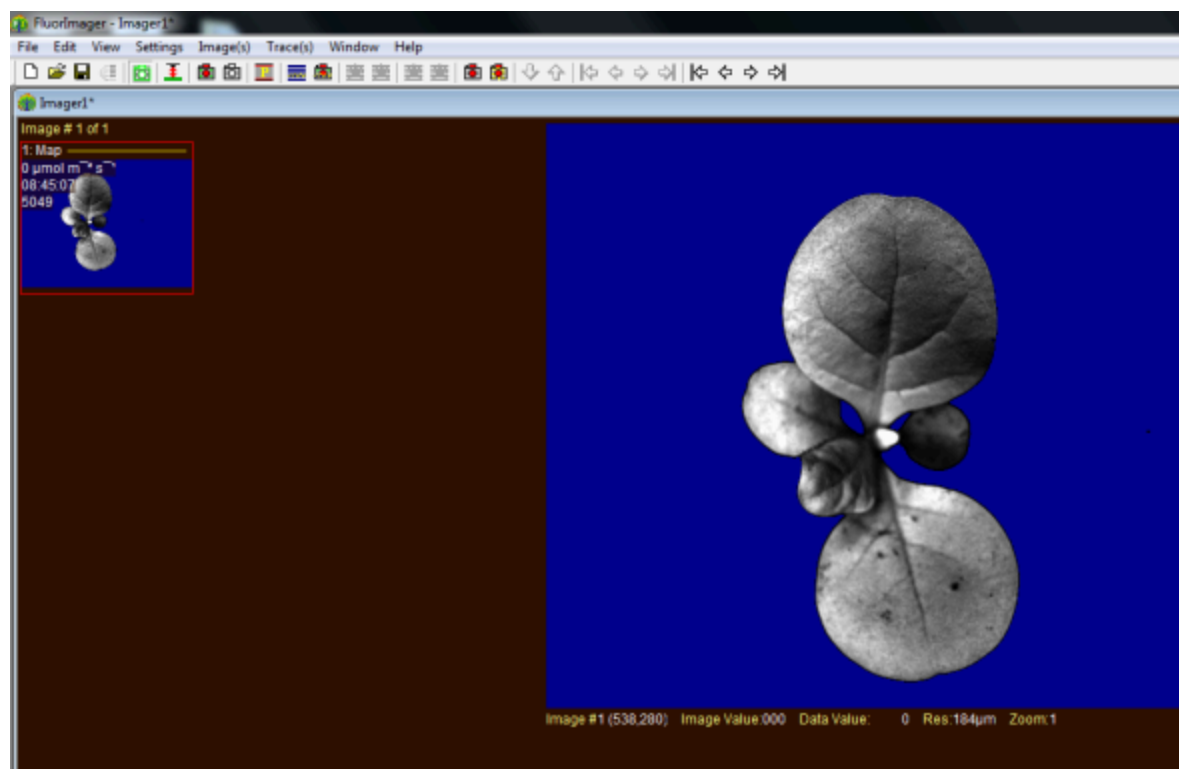


## 8 Remove residual noise from image

Sometimes the software picks up background noise as real signals, shown below as black dots on the blue background. It is advisable to mask these, otherwise they will be counted as a separate 'colony' during analysis and measurements will be recorded for each of these spots in the final data sheet. Noise can be masked by moving the cursor over the dot, pressing CTRL+left click simultaneously.

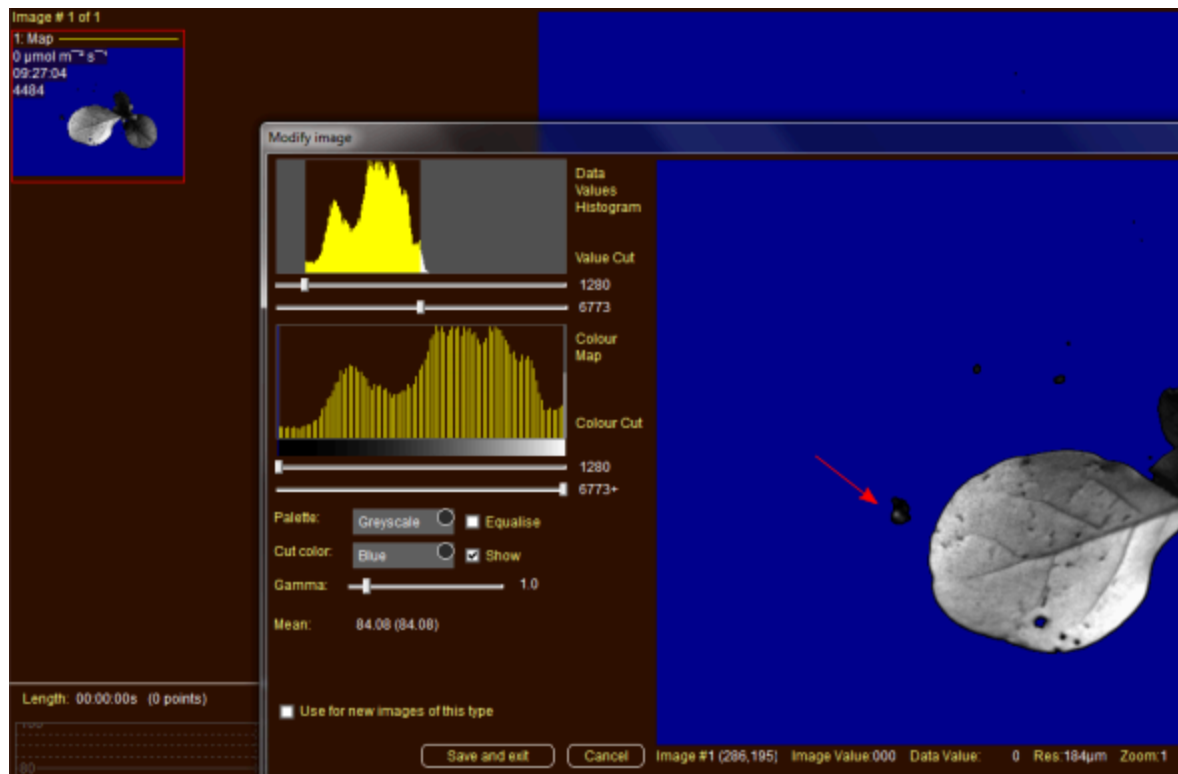

## 9 Load protocol



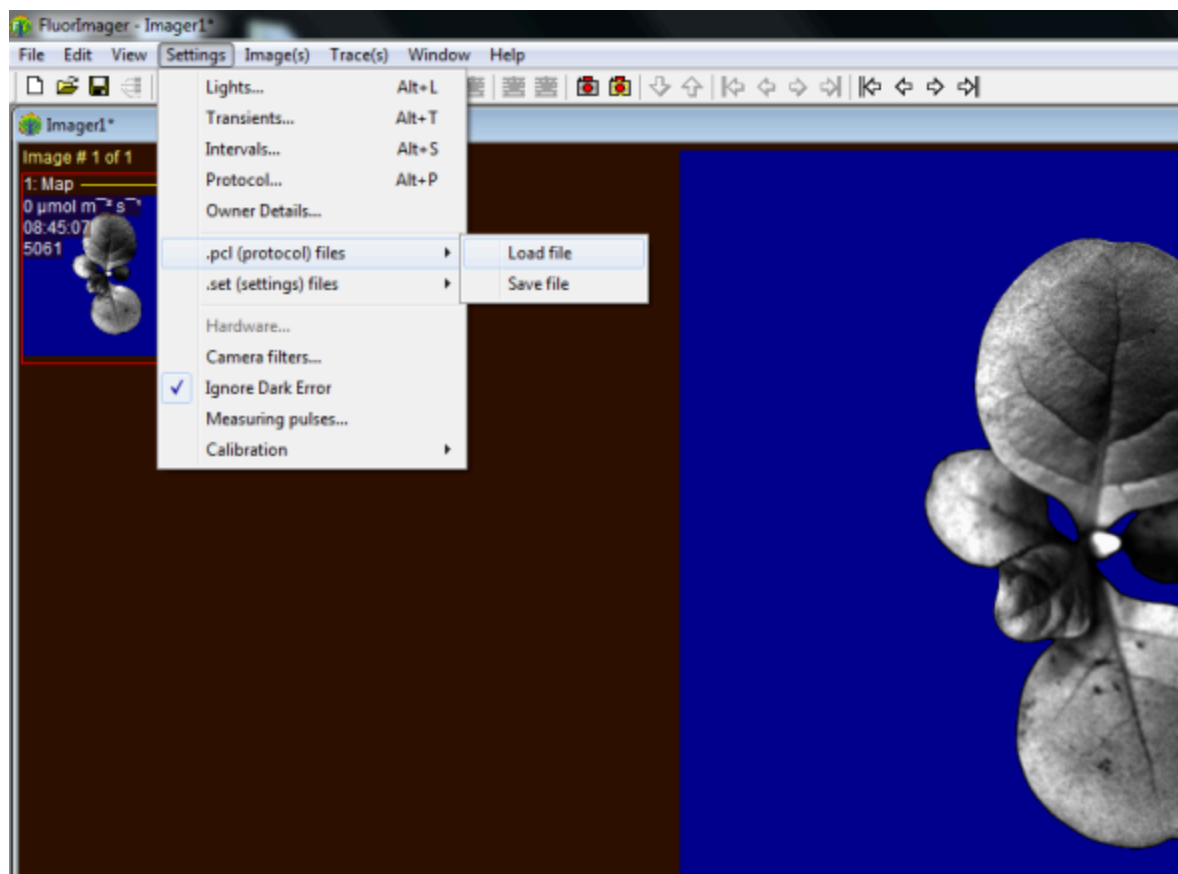

10 If the protocol needs change, press Settings > Protocol. Use the left menu to adjust the parameters you want to.

11



- 12 Once you are happy with the scheduled program click on the protocol icon in the toolbar (black arrow below) to start the run



- 13 Export data

After the reading is complete, save the images (File > Save As).

To obtain the files to process the images, go to File > Export to Folder. It will save in the folder that is opened.

