

Apr 15, 2019

DNA isolation from Formalin-Fixed, Paraffin-Embedded (FFPE) material

 [PLOS One](#)

DOI

dx.doi.org/10.17504/protocols.io.zjwf4pe

Hendrik F HFvE Van Essen¹

¹Amsterdam UMC



Hendrik F HFvE Van Essen

Amsterdam UMC

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.zjwf4pe

External link: <https://doi.org/10.1371/journal.pone.0223827>

Protocol Citation: Hendrik F HFvE Van Essen 2019. DNA isolation from Formalin-Fixed, Paraffin-Embedded (FFPE) material. [protocols.io](https://dx.doi.org/10.17504/protocols.io.zjwf4pe) <https://dx.doi.org/10.17504/protocols.io.zjwf4pe>

Manuscript citation:

Vincenten JPL, Essen HFv, Lissenberg-Witte BI, Bulkman NWJ, Krijgsman O, Sie D, Eijk PP, Smit EF, Ylstra B, Thunnissen E (2019) Clonality analysis of pulmonary tumors by genome-wide copy number profiling. PLoS ONE 14(10): e0223827. doi: [10.1371/journal.pone.0223827](https://doi.org/10.1371/journal.pone.0223827)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: March 27, 2019

Last Modified: April 15, 2019

Protocol Integer ID: 21846

Keywords: DNA, FFPE, QIAgen






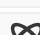

Abstract

This protocol describes the isolation of DNA from paraffin embedded tissue. The resulting DNA can be used for CGH array hybridizations.





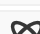

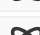
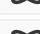
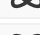
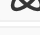
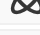



After removal of the paraffin, the tissue is lysed and the DNA is released from the cells. Then the DNA is isolated by using a silica-based column that binds the DNA.

Materials

STEP MATERIALS

-  Buffer ATL **Qiagen Catalog #19076**
-  Proteinase K **Thermo Fisher Scientific Catalog #EO0491**
-  Buffer AL, Lysis buffer **Qiagen Catalog #19076**
-  EtOH
-  Buffer AW1 **Qiagen Catalog #19081**
-  Buffer AW2 **Qiagen Catalog #19072**
-  Buffer AE **Qiagen Catalog #19077**

Protocol materials

-  Buffer ATL **Qiagen Catalog #19076**
-  Proteinase K **Thermo Fisher Scientific Catalog #EO0491**
-  Buffer AL, Lysis buffer **Qiagen Catalog #19076**
-  EtOH
-  Buffer AW1 **Qiagen Catalog #19081**
-  Buffer AW2 **Qiagen Catalog #19072**
-  Buffer AE **Qiagen Catalog #19077**
-  Buffer ATL **Qiagen Catalog #19076**
-  Proteinase K **Thermo Fisher Scientific Catalog #EO0491**
-  Buffer AL, Lysis buffer **Qiagen Catalog #19076**
-  EtOH
-  Buffer AW1 **Qiagen Catalog #19081**
-  Buffer AW2 **Qiagen Catalog #19072**
-  Buffer AE **Qiagen Catalog #19077**



H&E staining from paraffin blocks

- 1 Cut 3 μm (1x), 10 μm (2-3x, or more if needed depending on tumour area)
- 2 Mount all the section on slides with BSA 0.1%, and dry them o/n at 37°C (to a maximum of one week)
- 3 Place slides in xylene ⌚ 00:07:00 3 times
- 4 Hydrate by passing sequentially through 100%, 96%, 70% EtOH, and distilled-water
- 5 Stain slides with haematoxylin ⌚ 00:02:00
(if solution is not fresh staining time can be extended, max 4 min.)
- 6 Wash with running tap water ⌚ 00:05:00
- 7 Stain slides with eosin ⌚ 00:02:00
- 8 Wash quickly in water
- 9 Dehydrate by passing sequentially through 70%, 96%, 100% EtOH, xylene
- 10 Remove slides from xylene and use DePex to mount a covering glass.
- 11 Let the 3 μm H&E slide be judged by the pathologist for dissection.
Mark the tumor area and tumor percentage










DNA isolation

- 12 Place slides in xylene 3 \times 7 min. ⌚ 00:07:00 3 times











- 13 Hydrate by passing sequentially through 100%, 96%, 70% EtOH, and water
- 14 Stain slides with haematoxylin 1-2 min. ⌚ 00:02:00
(if solution is not fresh staining time can be extended, max 3 min.)
- 15 Wash with running tap water ⌚ 00:05:00
- 16 Place slides in distilled water
- 17 Dissect the tissue from the slides (when they are still slightly wet) with a scalpel or a needle
- 18 Place the material in a clean safelock eppendorf cup
- 19 Spin down in a centrifuge at full speed and get rid of the water layer ⌚ 00:05:00
- 20 Shake and invert samples (do not vortex).
- 21 Add 160 µl ATL buffer
⊗ Buffer ATL **Qiagen Catalog #19076**
- 22 Add 40 µl prot K (20 mg/mL)
⊗ Proteinase K **Thermo Fisher Scientific Catalog #EO0491**
- 23 vortex sample ⌚ 00:00:15
Make sure all tissue is in the liquid
- 24 Incubate o/n at 56°C in heat-block or waterbath and vortex regularly! ⌚ 08:00:00
🔥 56 °C



- Check if all tissue is digested
- 25 Vortex samples and spin down
- 26 Incubate at 98°C 10 min.  98 °C  00:10:00
- 27 Spin down at 16600 rcf. 1 min.  00:01:00
- 28 Add 200 µl AL buffer and mix very well by vortexing!
 Buffer AL, Lysis buffer **Qiagen Catalog #19076**
- 29 Add 200 µl Ethanol 100% and mix very well by vortexing!
 EtOH
- 30 Incubate at RT  00:05:00
- 31 Spin down the eppendorf cups at full speed 1 min.
- 32 Transfer up to 600 µl lysate to the QIAamp MinElute Column.
- 33 Spin down at 8.000 rpm  00:01:00
- 34 Place the QIA Column in a clean 2 ml collection tube, discard the flow-through make sure the column tip is clear of droplets
- 35 Add 500 µl AW1 buffer to the column
 Buffer AW1 **Qiagen Catalog #19081**
- 36 Spin down for at 8.000 rpm and place the column in a new tube  00:01:00



- 37 Add 500 μ l AW2 buffer (QIAamp micro-kit) to the column
 Buffer AW2 **Qiagen Catalog #19072**
- 38 Spin down at 8.000 rpm and place the column in a new tube  00:01:00
- 39 Spin down at full speed to dry the membrane  00:03:00
- 40 Place the QI column in a properly marked eppendorf cup, discard collection tube
- 41 Add 20 to 50 μ l (depending on sample size) of AE buffer to the column and incubate.
 00:05:00  20 °C room temperature
 Buffer AE **Qiagen Catalog #19077**
- 42 Spin down at full speed.  00:03:00
- 43 Throw away the column, close the eppendorf cup and store the DNA at  4 °C