

May 08, 2018 Version 4

🌍 Measuring spectral reflectance and transmittance (350-2500 nm) of large leaves using the Spectra Vista Corporation (SVC) DC-R/T Integrating Sphere V.4



DOI

dx.doi.org/10.17504/protocols.io.pzbdp2n

Etienne Laliberté¹, Raymond Soffer²

¹Université de Montréal; ²National Research Council

Canadian Airborne Biodiversity Observatory
Tech. support email: jocelyne.ayotte@umontreal.ca



Etienne Laliberté

Université de Montréal

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.pzbdp2n

External link: www.caboscience.org

Protocol Citation: Etienne Laliberté, Raymond Soffer 2018. Measuring spectral reflectance and transmittance (350-2500 nm) of large leaves using the Spectra Vista Corporation (SVC) DC-R/T Integrating Sphere. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.pzbdp2n>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: In development

We are still developing and optimizing this protocol.

Created: May 08, 2018

Last Modified: May 30, 2018

Protocol Integer ID: 12035



Abstract

Here we describe the standardised protocol used by the Canadian Airborne Biodiversity Observatory (CABO) to measure leaf spectral reflectance and transmittance, using the Spectra Vista Corporation (SVC) DC-R/T Integrating Sphere fitted to a portable full-range SVC HR-1024i field spectroradiometer. This standard version of our protocol describes the common case where an individual leaf is large enough to entirely cover the reflectance or transmission port of the integrating sphere. Briefly, six mature, healthy-looking and sunlit leaves from a canopy plant are selected for measurements of adaxial reflectance and transmittance. Leaf scans are referenced to a calibrated Spectralon® disk and corrected for stray light. Our leaf spectroscopy protocol builds from that of the Carnegie Airborne Observatory.

Guidelines

Handling Spectralon®

Do not touch Spectralon® (e.g. sphere interior, reference disks, plugs) with your fingers.

Do not use canned air to remove dust on the Spectralon® disk; canned air contains chemicals that can alter Spectralon®'s optical properties.

Do not attempt to clean Spectralon® in the field, other than **blowing surface dust only on the Spectralon® reference disk or sphere plugs** using the Canless Air Duster System; cleaning Spectralon® requires a special procedure that should only be done in the lab.

Never blow air inside of the integrating sphere, especially not when it is attached to the spectroradiometer, as this will blow dust inside the instrument.

Equipment

Spectra Vista Corporation HR-1024i full-range (350–2500 nm) field spectroradiometer

Spectra Vista Corporation 3-inch Spectralon® DC-R/T Sphere

Semi-rugged laptop or PDA running the SVC Scan software

Canless Air Duster System O₂ Hurricane (**never use canned air**) to remove dust from the surface of the Spectralon® reference disk

Plastic containers with lids to temporarily store leaf samples during measurements (optional)

Consumables

Nitrile gloves for handling leaves

Safety warnings

! The lamp of the integrating sphere can get **very hot** and should be handled from its slotted base to avoid burns.



Before start

Consult the user manual of the spectroradiometer and the integrating sphere to set up the instrument.

The instrument should be set up in the shade, sheltered as much as possible from the elements.

All canopy plants selected for measurements should have already been tagged, identified, and georeferenced before spectroscopy measurements start.

The spectroscopist should be positionned as close as possible to the sampled plants to minimise time from collection to measurement.

The spectroscopist should be in a comfortable position and have enough room around the instrument to spread leaf samples around without the risk of mixing up individual leaves during handling.

Six mature, fully-developped, healthy-looking leaves from the sunlit (>3 h per day of direct sunlight) portion of the canopy are selected for spectral measurements from the bulk leaf sample (often one of a few branches). Leaves should be collected from the uppermost surface of the branch (i.e. receiving the most direct sunlight),

Leaves should be stored in a sealed, clear plastic bag (breathe into the bag before sealing it) and **brought immediately to the spectroscopist for measurements.**

Instrument Set-Up

- 1 Install the integrating sphere onto the spectroradiometer.

Note

Follow the SVC integrating sphere manual p. 9-14.

- 2 Power the spectroradiometer and integrating sphere lamp on and **warm up for >15 min.**

 00:15:00 Spectroradiometer and lamp warm-up time

Note

Record the **time of day** when the spectroradiometer was started.

Review Protocol Summary Diagram

- 3 Review the **document** summarising the different sphere configurations (A–E), and the scans that need to be recorded in each configuration.



spectrostandard-4.pdf

Configuration A: Reflectance Mode, Reference

- 4 Position the lamp over the sphere **primary light entrance port**.

**Note**

Make sure lamp is secured in locked position.

Safety information

The lamp can get very hot. Grab it by the slotted heat shield.

- 5 Check lamp alignment.

Note

Use a thin piece of paper at the exit of the reflectance sample port (empty port) to ensure the light beam under-fills and is centered in the reflectance port. **If it is not, then proceed to lamp alignment** as described in the SVC integrating sphere user manual, p. 23-24.

- 6 Screw the tethered light trap on the **reflectance port** sample holder.

Note

The light trap can stay on the sample holder for the entire measurement session.

- 7 Screw the tethered light trap on the **transmission port** sample holder.

Note

The light trap can stay on the sample holder for all measurements made in reflectance mode (configurations A–C).

- 8 Place the tethered calibrated Spectralon® reflectance standard over the **reflectance port**.

Note


Place the standard over the reflectance port so that the light beam shines directly on its reflective surface (i.e facing into the sphere).

- 9 Position leaf #1 over the **transmission port** with its adaxial (upper) surface facing into the sphere.

Note

Position the leaf so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 10 Collect a '**Reference Scan**' in this configuration.

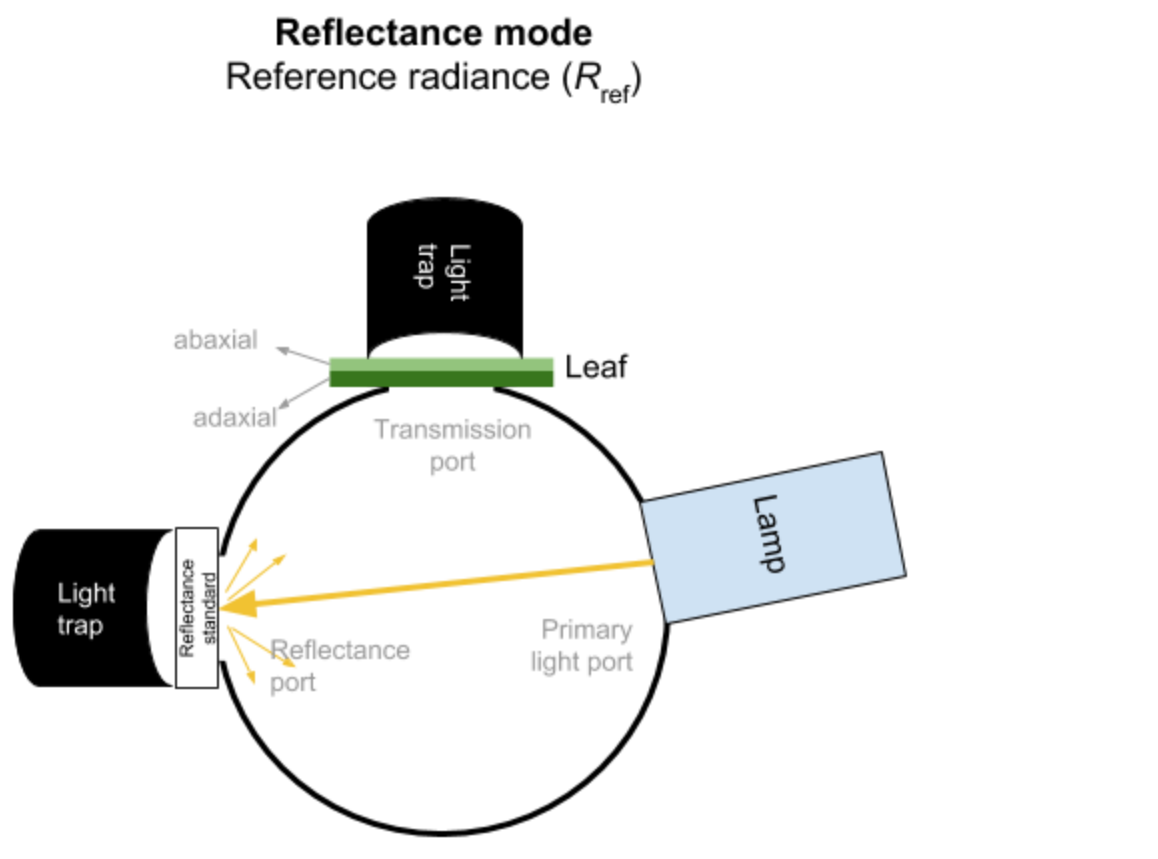

Note

This corresponds to the **reference radiance** in reflectance mode (R_{ref}). The reference data will be automatically saved in all successive target scan files until a new 'Reference Scan' is made.


- 11 Collect a '**Target Scan**' in this configuration and **save the file**.

Configuration B: Reflectance Mode, Stray light

- 12 Carefully remove leaf #1 from the transmission port sample holder.
- 13 Remove the tethered calibrated Spectralon® reflectance standard from the reflectance port.
- 14 Place the tethered calibrated Spectralon® reflectance standard over the **transmission port** sample holder.
- 15 Collect a '**Target Scan**' in this configuration and **save the file**.

Reflectance mode

Stray light radiance (R_{str})

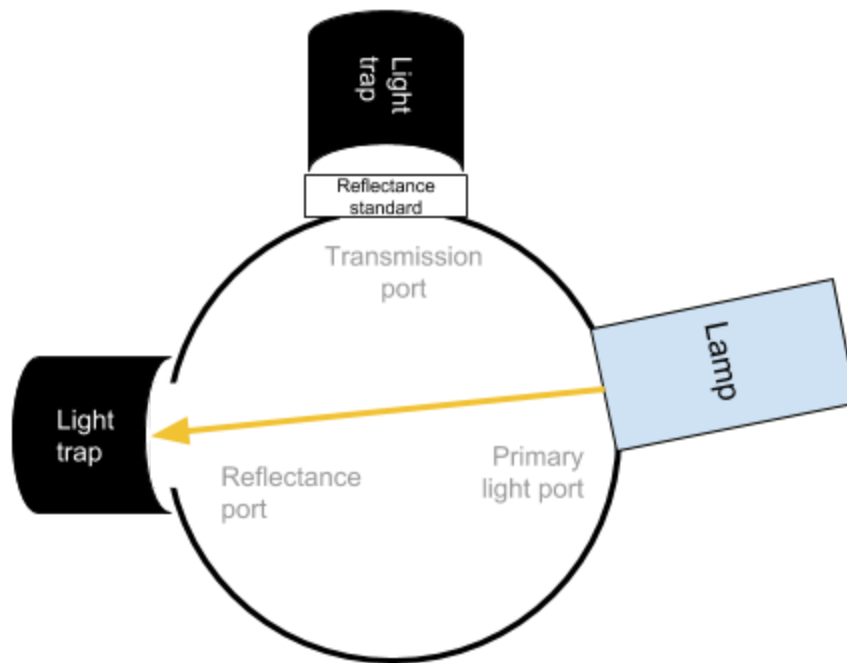

Note

This corresponds to the **stray light radiance** in reflectance mode (R_{str}).

Configuration C: Reflectance Mode, Target

- 16 Position leaf #1 over the **reflectance port** with its adaxial (upper) surface facing into the sphere.

Note


Position the leaf to target the same area measured for the reference radiance. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

Note

The light trap should remain on the reflectance port sample holder.

- 17 Collect a '**Target Scan**' for leaf #1 in this configuration and **save the file**.

Reflectance mode
Target radiance, leaf i ($R_{\text{tar},i}$)

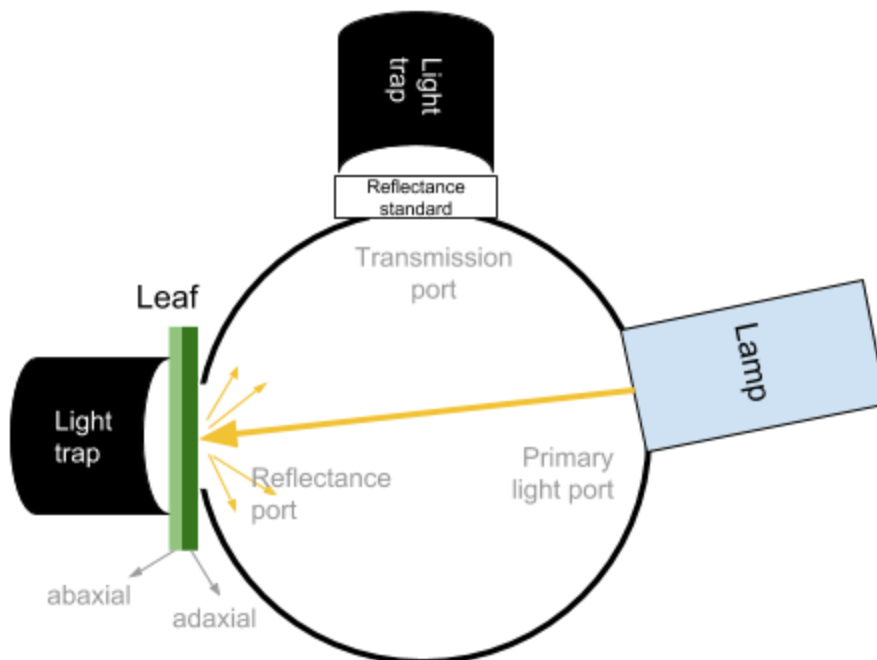

Note

This corresponds to the **target radiance** in reflectance mode for leaf #1 ($R_{\text{tar},1}$).

- 18 Carefully replace leaf #1 by leaf #2.

Note

Position the leaf to target the same area measured for the reference radiance. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 19 Collect a '**Target Scan**' for leaf #2 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in reflectance mode for leaf #2 ($R_{\text{tar},2}$).

- 20 Carefully replace leaf #2 by leaf #3.

Note

Position the leaf to target the same area measured for the reference radiance. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 21 Collect a '**Target Scan**' for leaf #3 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in reflectance mode for leaf #3 ($R_{\text{tar},3}$).

- 22 Carefully replace leaf #3 by leaf #4.

Note

Position the leaf to target the same area measured for the reference radiance. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 23 Collect a '**Target Scan**' for leaf #4 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in reflectance mode for leaf #4 ($R_{\text{tar},4}$).

- 24 Carefully replace leaf #4 by leaf #5.

Note

Position the leaf to target the same area measured for the reference radiance. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 25 Collect a '**Target Scan**' for leaf #5 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in reflectance mode for leaf #5 ($R_{\text{tar},5}$).

- 26 Carefully replace leaf #5 by leaf #6.

Note

Position the leaf to target the same area measured for the reference radiance. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 27 Collect a '**Target Scan**' for leaf #6 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in reflectance mode for leaf #6 ($R_{\text{tar},6}$).

Configuration A: Reflectance Mode, Reference

- 28 Remove leaf #6 from the reflectance sample port holder.
- 29 Place the tethered calibrated Spectralon® reflectance standard over the **reflectance port**.

Note


Place the standard over the reflectance port so that the light beam shines directly on its reflective surface (i.e facing into the sphere).

- 30 Position leaf #1 over the **transmission port** with its adaxial (upper) surface facing into the sphere.

Note

Position the leaf so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 31 Collect a '**Target Scan**' in this configuration and **save the file**.

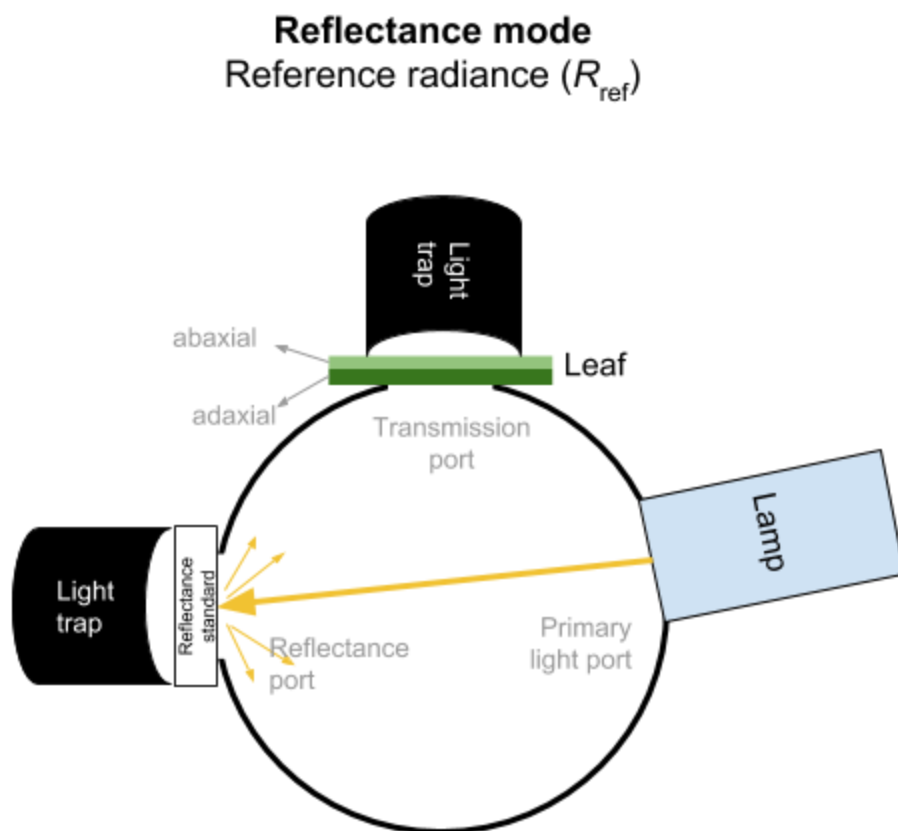

Note

This second reference radiance scan is only used to assess the stability of the system in reflectance mode.

Configuration D: Transmission Mode, Reference

- 32 Remove the tethered calibrated Spectralon® reflectance standard from the sphere transmission port.
- 33 Remove leaf #1 from the transmission port sample holder.
- 34 Position leaf #1 over the **reflectance port** so that its abaxial (lower) side is now facing the inside of the sphere.

Note

Position the leaf to target the same area measured for the reflectance radiance, with the exception that its abaxial surface now faces the inside of the sphere. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 35 Remove the light trap from the transmission port sample holder.



36 Position the lamp over the sphere **transmission port**.

Note

Make sure lamp is secured in locked position.

Safety information

The lamp can get very hot. Grab it by the slotted heat shield.

37 Install the Spectralon® plug over the **primary light port**.


Note

Ensure that the curved plug is placed the correct way to match the curvature of the sphere.

38 Collect a '**Reference Scan**' in this configuration.

Transmission mode

Reference radiance (T_{ref})

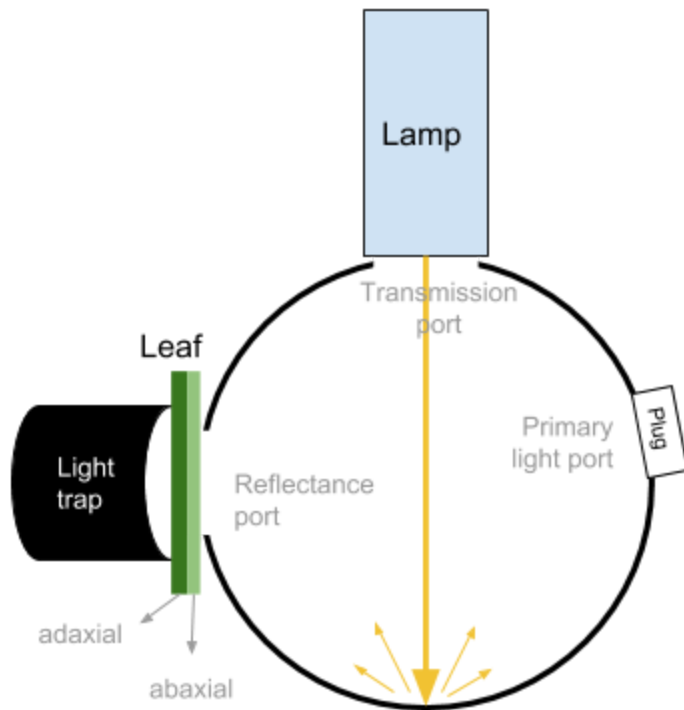

Note

This corresponds to the **reference radiance** in **transmission** mode (T_{ref}).

- 39 Collect a '**Target Scan**' in this configuration and **save the file**.

Configuration E: Transmission Mode, Target

- 40 Carefully remove leaf #1 from the reflectance port sample holder.

Note

The reflectance port should now be **empty** (but with **light trap on**).

- 41 Gently pull lamp away from the sphere.

Safety information

The lamp can get very hot. Grab it by the slotted heat shield.


- 42 Place leaf #1 over the **transmission port** with its **abaxial** (lower) surface facing into the sphere.

Note

Position the leaf to target the same area measured for the reflectance radiance, with the exception that its abaxial surface now faces into the sphere. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 43 Release the transmission sample holder and move lamp back to its locked position.

- 44 Collect a '**Target Scan**' for leaf #1 in this configuration and **save the file**.

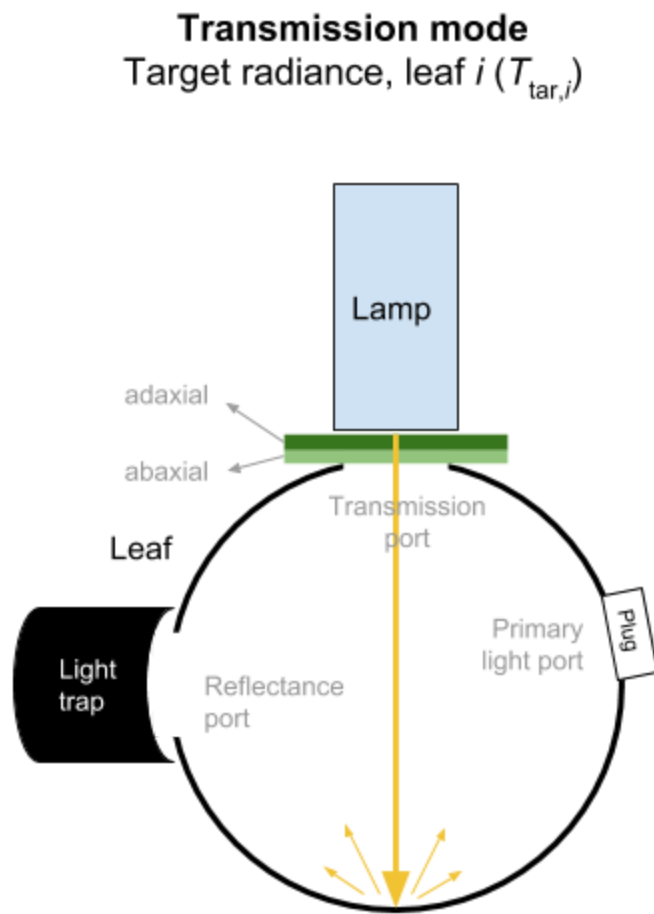

Note

This corresponds to the **target radiance** in **transmission** mode for leaf #1 ($T_{\text{tar},1}$).

- 45 Carefully replace leaf #1 by leaf #2.

Note

Position the leaf to target the same area measured for the reflectance radiance, with the exception that its abaxial surface now faces into the sphere. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 46 Collect a '**Target Scan**' for leaf #2 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in **transmission** mode for leaf #2 ($T_{\text{tar},2}$).

- 47 Carefully replace leaf #2 by leaf #3.

Note

Position the leaf to target the same area measured for the reflectance radiance, with the exception that its abaxial surface now faces the inside of the sphere. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 48 Collect a '**Target Scan**' for leaf #3 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in **transmission** mode for leaf #3 ($T_{\text{tar},3}$).

- 49 Carefully replace leaf #3 by leaf #4.

Note

Position the leaf to target the same area measured for the reflectance radiance, with the exception that its abaxial surface now faces the inside of the sphere. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 50 Collect a '**Target Scan**' for leaf #4 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in **transmission** mode for leaf #4 ($T_{\text{tar},4}$).

- 51 Carefully replace leaf #4 by leaf #5.

Note

Position the leaf to target the same area measured for the reflectance radiance, with the exception that its abaxial surface now faces the inside of the sphere. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 52 Collect a '**Target Scan**' for leaf #5 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in **transmission** mode for leaf #5 ($T_{\text{tar},5}$).

- 53 Carefully replace leaf #5 by leaf #6.

Note

Position the leaf to target the same area measured for the reflectance radiance, with the exception that its abaxial surface now faces the inside of the sphere. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 54 Collect a '**Target Scan**' for leaf #6 in this configuration and **save the file**.

Note

This corresponds to the **target radiance** in **transmission** mode for leaf #6 ($T_{\text{tar},6}$).

Configuration D: Transmission Mode, Target

- 55 Remove leaf #6 from the transmission sample port holder.

Configuration D: Transmission Mode, Reference


- 56 Position leaf #1 over the **reflectance port** so that its abaxial (lower) side is now facing the inside of the sphere.

Note

Position the leaf to target the same area measured for the reflectance radiance, with the exception that its abaxial surface now faces the inside of the sphere. Position it so that the amount of leaf and vein material over the port is roughly proportional to the area of leaf and vein found throughout the leaf, while avoiding the large midrib vein. Position the leaf so that it is approximately halfway between the mid-rib vein and the leaf margin, and halfway between the tip and the base of the leaf lamina.

- 57 Collect a '**Target Scan**' in this configuration and **save the file**.

Transmission mode Reference radiance (T_{ref})

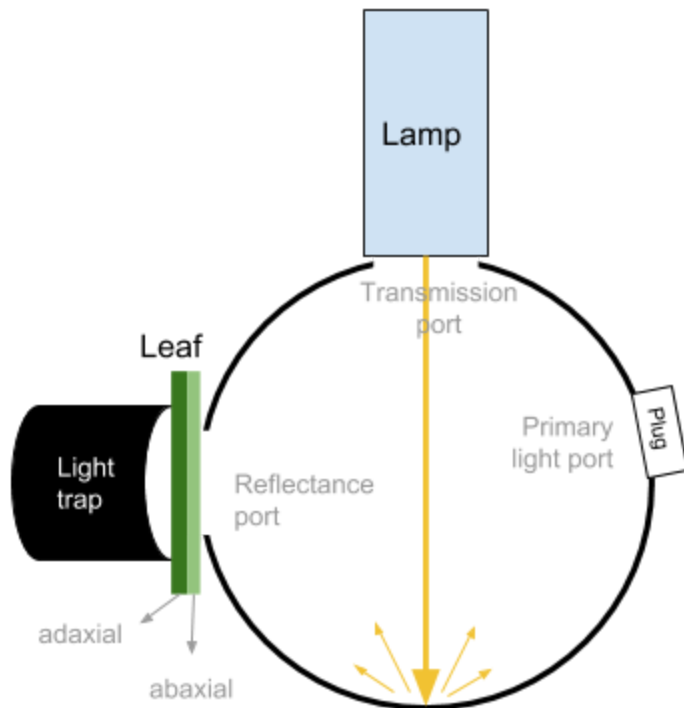

Note

This second reference radiance scan is only used to assess the stability of the system in transmission mode.

Calculating Leaf Reflectance (Adaxial Surface)

58 The equation for **adaxial reflectance** of leaf i , $\rho_{\text{leaf},i}$ is

$$\rho_{\text{leaf},i} = [(R_{\text{tar},i} - R_{\text{str}}) \div (R_{\text{ref}} - R_{\text{str}})] \times \rho_{\text{ref}}$$

where



ρ_{ref} is the absolute reflectance of the calibrated Spectralon® reflectance standard.

Calculating Leaf Transmittance (Adaxial Surface)

59 The equation for **adaxial transmittance** of leaf i , $\tau_{\text{leaf},i}$ is

$$\tau_{\text{leaf},i} = (T_{\text{tar},i} \div T_{\text{ref}})$$