

Jan 27, 2015

Version 1

Gibson Assembly® Protocol (E5510) V.1

DOI

<https://dx.doi.org/10.17504/protocols.io.cdms45>



New England Biolabs (NEB)

Tech. support phone: +1(800)632-7799 email: info@neb.com



New England Biolabs

New England Biolabs

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

[Create free account](#)

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.cdms45>

External link: <https://www.neb.com/protocols/2012/12/11/gibson-assembly-protocol-e5510>

Protocol Citation: New England Biolabs: Gibson Assembly® Protocol (E5510). **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.cdms45>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: July 10, 2014

Last Modified: March 28, 2018

Protocol Integer ID: 141

Keywords: cloning, Gibson cloning, Gibson HIFI, transformation protocol, transformed competent E.coli cells with gibson assembly, protocol for the gibson assembly, gibson assembly, e5510, protocol

Abstract

This is the protocol for the Gibson Assembly using the Gibson Assembly® Cloning Kit (E5510).

Guidelines

Optimal Quantities

NEB recommends a total of 0.02–0.5 pmols of DNA fragments when 1 or 2 fragments are being assembled into a vector and 0.2–1.0 pmoles of DNA fragments when 4–6 fragments are being assembled. Efficiency of assembly decreases as the number or length of fragments increases. To calculate the number of pmols of each fragment for optimal assembly, based on fragment length and weight, we recommend using NEB's online tool, [NEBioCalculator](#), or using the following formula:

$$\text{pmols} = (\text{weight in ng}) \times 1,000 / (\text{base pairs} \times 650 \text{ daltons})$$

50 ng of 5000 bp dsDNA is about 0.015 pmols.

50 ng of 500 bp dsDNA is about 0.15 pmols.

The mass of each fragment can be measured using the NanoDrop instrument, absorbance at 260 nm or estimated from agarose gel electrophoresis followed by ethidium bromide staining.

Optimized cloning efficiency is 50–100 ng of vectors with 2–3 fold of excess inserts. Use 5 times more of inserts if size is less than 200 bps. Total volume of unpurified PCR fragments in Gibson Assembly reaction should not exceed 20%.

OVERVIEW:


Gibson Assembly was developed by Dr. Daniel Gibson and his colleagues at the J. Craig Venter Institute and licensed to NEB by Synthetic Genomics, Inc. It allows for successful assembly of multiple DNA fragments, regardless of fragment length or end compatibility. It has been rapidly adopted by the synthetic biology community due to its ease-of-use, flexibility and suitability for large DNA constructs.

Gibson Assembly efficiently joins multiple overlapping DNA fragments in a single-tube isothermal reaction (1,2). The Gibson Assembly Master Mix includes three different enzymatic activities that perform in a single buffer:

- The exonuclease creates single-stranded 3' overhangs that facilitate the annealing of fragments that share complementarity at one end (overlap region).
- The proprietary DNA polymerase fills in gaps within each annealed fragment.
- The DNA ligase seals nicks in the assembled DNA.

The end result is a double-stranded fully sealed DNA molecule that can serve as template for PCR, RCA or a variety of other molecular biology applications, including direct transformation. The method has been successfully used by Gibson's group and others to assemble oligonucleotides, DNA with varied overlaps (15–80 bp) and fragments hundreds of kilobases long (1–2).

Overview of the Gibson Assembly Cloning Method

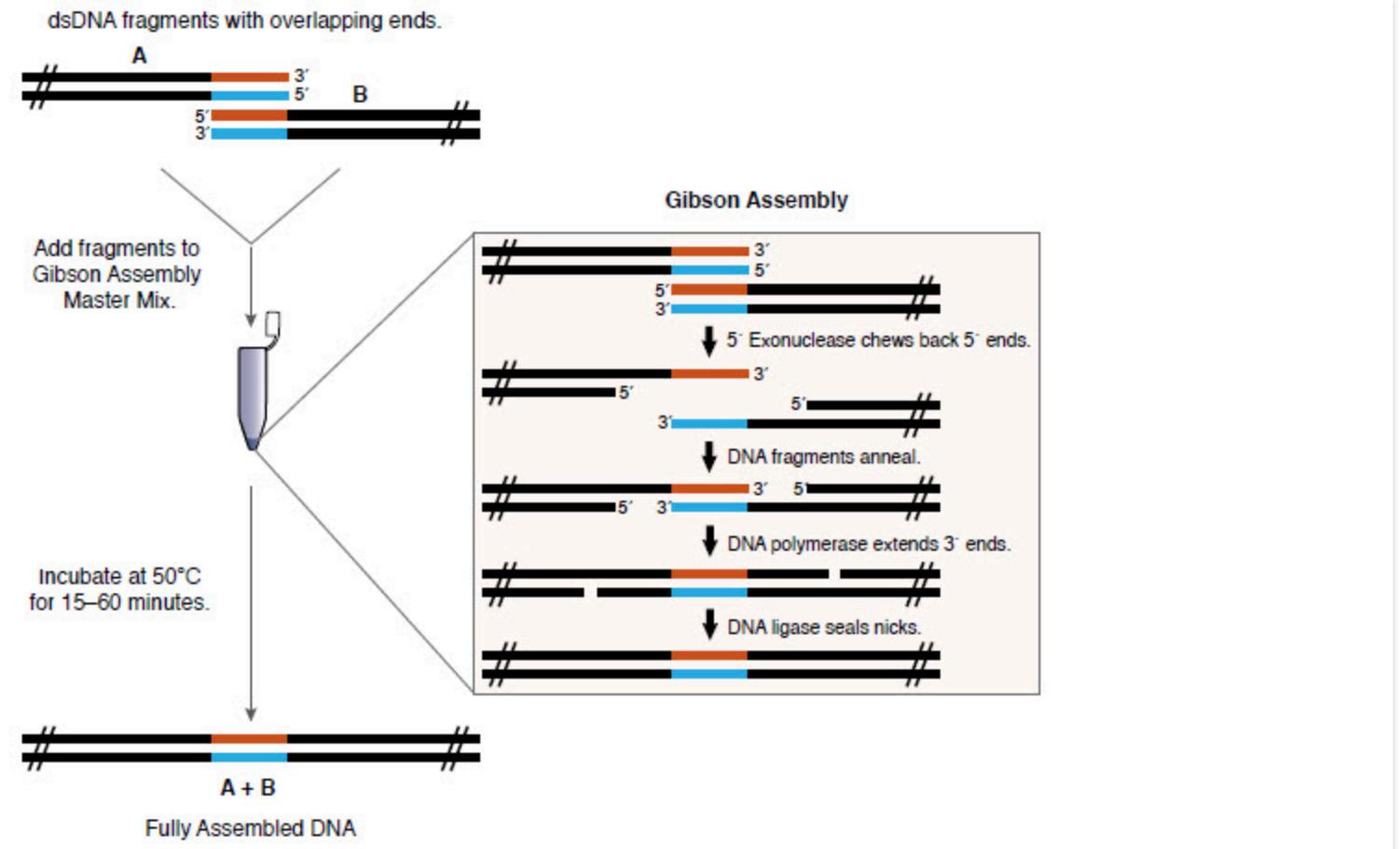

Overview of Gibson Assembly Cloning Kit Protocol:

- Design primers to amplify fragments (and/or vector) with appropriate overlaps
- PCR amplify fragments using a high-fidelity DNA polymerase.
- Prepare linearized vector by PCR amplification using a high-fidelity DNA polymerase or by restriction digestion.
- Confirm and determine concentration of fragments and linearized vector using agarose gel electrophoresis, a NanoDrop™ instrument or other method.
- Add fragments and linearized vector to Gibson Assembly Master Mix and incubate at 50°C for 15 minutes to 1 hour, depending on number of fragments being assembled.
- Transform into NEB 5-alpha Competent E. coli (provided) or use directly in other applications.

NOTES:

1. We highly recommend using our web tool, [NEBuilder®](#) to design PCR primers with overlapping sequences between the adjacent DNA fragments and for their assembly into a cloning vector.

2. Storage Note: The kit is shipped on dry ice. Upon arrival, store kit at -80°C. After first use, store the kit components at indicated temperatures.

3. Usage notes:

To ensure the successful assembly and subsequent transformation of assembled DNAs, NEB recommends the following:

- *DNA*: PCR product purification is not necessary if the total volume of all PCR products in the Gibson Assembly reaction is 20% or less of the Gibson Assembly reaction volume. Higher volumes of PCR products may reduce the efficiency of Gibson Assembly and transformation due to the elevated carryover amounts of PCR reaction buffer and unused primers present in the PCR product. Column purification of PCR products may increase the efficiency of both Gibson Assembly and transformation by 2–10 fold and is highly recommended when performing assemblies of three or more PCR fragments or assembling longer than 5 kb fragments. Purified DNA for assembly can be dissolved in ddH₂O (Milli-Q® water or equivalent is preferable), TE or other dilution buffers.
- *Insert*: When directly assembling fragments into a cloning vector, the concentration of assembly fragments should be at least 2–3 times higher than the concentration of vector. **For assembly of multiple fragments into a vector, we recommend using equimolar ratio of fragments.**
- *Transformation*: NEB 5-alpha Competent E. coli (High Efficiency, NEB #C2987) provided with the kit are recommended for use for assembled products of less than 20 kb in size. It is also possible to use other NEB competent E. coli strains, with the exception of BL21, BL21(DE3), Lemo21(DE3) and Nico21(DE3). For example, Shuffle T7 Express Competent E. coli can be used for the expression of a difficult to express protein. When using competent E. coli from a vendor other than NEB, we have seen decreased robustness of transformation with the Gibson Assembly reaction.
- *Electroporation*: Electroporation can increase transformation efficiency by several logs. When using the Gibson Assembly Master Mix product for electroporation, it is necessary to dilute the reaction 3-fold and use 1 µl for transformation. Should you require the use of Electrocompetent cells, please use the Electrocompetent Cells Transformation Protocol.

REFERENCES:

1. Gibson, D.G. et.al. (2009). Nature Methods. 343-345.
2. Gibson, D.G. et al. (2010). Nature Methods. 901-903.
3. Barnes, W.M. (1994). Proc. Natl. Acad. Sci.. 91, 2216-220.

Materials

MATERIALS

 Gibson Assembly Cloning Kit - 10 rxns **New England Biolabs Catalog #E5510S**

Troubleshooting

1 Set up the following reaction on ice:

	2-3 Fragment Assembly	4-6 Fragment Assembly	Positive Control
Total Amount of Fragments	.02-0.5 pmol X ul	0.2-1 pmol X ul	10 ul
Gibson Assembly Master Mix (2X)	10 µl	10 µl	10 µl
Deionized H ₂ O	10-X µl	10-X µl	0
Total Volume	20 µl	20 µl	20 µl

Note

Control reagents are provided for 5 experiments.

Note

If greater numbers of fragments are assembled, additional Gibson Assembly Master Mix may be required.

Note

Optimized cloning efficiency is 50–100 ng of vectors with 2–3 fold of excess inserts. Use 5 times more of inserts if size is less than 200 bps. Total volume of unpurified PCR fragments in Gibson Assembly reaction should not exceed 20%.

Protocol



NAME

Mixture for E5510 Gibson

CREATED BY

New England Biolabs

Preview

1.1 DNA fragments

1.2 Gibson Assembly Master Mix (2X)

1.3 Deionized H₂O

- 2 Incubate samples in a thermocycler at 50°C for 15 minutes when 2 or 3 fragments are being assembled or 60 minutes when 4-6 fragments are being assembled.

Note

Extended incubation up to 60 minutes may help to improve assembly efficiency in some cases (for further details see [FAQ section](#)).

- 3 Store samples on ice or at -20°C for subsequent transformation.
- 4 Transform NEB 5-alpha Competent E. coli cells (provided with the kit) with 2 µl of the assembly reaction, following the transformation protocol.

 2 µL