



Sep 27, 2019

Version 2

## Agarose gel electrophoresis V.2

DOI

[dx.doi.org/10.17504/protocols.io.7q7hmzn](https://dx.doi.org/10.17504/protocols.io.7q7hmzn)

Cristian Riccio<sup>1</sup>

<sup>1</sup>University of Cambridge



Cristian Riccio

University of Cambridge

### Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.7q7hmzn>

**Protocol Citation:** Cristian Riccio 2019. Agarose gel electrophoresis. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.7q7hmzn>

**License:** This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**We use this protocol and it's working**

**Created:** September 27, 2019

**Last Modified:** September 27, 2019

**Protocol Integer ID:** 28159

**Keywords:** agarose, agarose gel, agarose gel electrophoresis, electrophoresis, agarose gel electrophoresis agarose gel for dna electrophoresis, agarose gel electrophoresis agarose gel, dna electrophoresis, dna fragment, gdna extraction, smaller than for the dna, hmw gdna extraction, gdna, dna, primer, separating pcr amplicon, fragments of different size, extraction, pcr amplicon, gel, pcr

## Abstract


Agarose gel for DNA electrophoresis

This protocol is useful to check the approximate size of a mixture of DNA fragments.

For example, after a HMW gDNA extraction, it is good to check that the DNA we extracted is actually high-molecular weight. Another use is for separating PCR amplicons when checking for a mutation. Primers can be designed in such a way that PCR will yield two fragments of different sizes depending on whether a sequence is present or not. If a deletion is present, then the amplified fragment will be smaller than for the DNA where there is no deletion.

## Materials

### MATERIALS

-  Gel Loading Dye Orange (6X) - 4.0 ml **New England Biolabs Catalog #B7022S**
-  TAE Buffer (Tris-acetate-EDTA) **Catalog #B49**
-  Agarose **Merck MilliporeSigma (Sigma-Aldrich)**
-  Image Lab Software **Bio-Rad Laboratories**
-  Autoclave Indicator Tapes **Bio Basic Inc. Catalog #AT011.SIZE.1ROLLS**
-  Ethidium bromide 10 mg/ml **Merck MilliporeSigma (Sigma-Aldrich) Catalog #E1510**
-  PowerPac™ HC Power Supply **Bio-Rad Laboratories Catalog #1645052**
-  agarose gel container
-  Bio-Rad electrophoresis machine
-  Bio-Rad Gel Doc XR with Image Lab Software
-  Microwave oven
-  HyperLadder 100 bp **Bioline Catalog #BIO-33030**
-  HyperLadder 1 kb **Bioline Catalog #BIO-33026**

## Troubleshooting




## Safety warnings

! Ethidium bromide is a carcinogen and you should be trained on how to handle it before use. EtBr also needs to be disposed of in a special way.



1 Get a 250 ml flask and a 250 ml measuring cylinder.

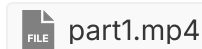

2 Add 1 gram agarose to the flask. Add 130 ml TAE in measuring cylinder. Taper container. Put comb in container.

3 Microwave flask with agarose + TAE for 1 minute. After one minute, check that the solution is homogeneous. If not, microwave in steps of 15 seconds until it is homogeneous.

4 Move taped and combed container to the cold room along with the quasi-boiling agarose solution. Add 1 microlitre ethidium bromide to the agarose solution. Pour the agarose solution into the container. Make sure there is no spillage. If there is spillage, stop pouring, clean the mess and start again.

5 Add a small amount of water to the flask, rinse it and dispose of the contaminated water in the liquid EtBr waste through the white material.

6 Wash flask with water and leave it in the dirty dishes trolley.

7 Leave the gel to solidify in the cold room. If you want to use the gel the following day, wrap it in a plastic sheet to protect it from desiccation.

2h

8 Remove tape, put gel in electrophoresis tank. Add 5 microlitres ladder.

9 Add 2 microlitres loading buffer to 5 microlitres of DNA sample (~ 100 nanograms) and load into gel.

10 Attach the cables from the tank to the Bio-Rad power machine. Set the programme to 150 V, 400 mA and 30 minutes. Start run.

30m

11 Check gel every 10 minutes. Put gel in Imager. Open Image Lab. New Protocol, Position gel, zoom in/out. Start protocol. Image is ready on computer. Invert contrast, change contrast, print image. Don't touch computer without gloves. It is contaminated.





## Acknowledgments

---

- 12 Jérémie Le Pen and other colleagues in the Miska lab  
Ulrika Yuan for recording