

Aug 01, 2025

Analysis of Islet Function in Dynamic Cell Perifusion System

DOI

<https://dx.doi.org/10.17504/protocols.io.eq2ly4nxrlx9/v1>



IIDP-HIPP Brissova¹

¹IIDP-HIPP

IIDP-HIPP Brissova: Integrated Islet Distribution Program and Human Islet Phenotyping Program

Integrated Islet Distribution Program and Human Islet Phenotyping Program
Tech. support email: anastasia.coldren@vumc.org



Anastasia Coldren

Vanderbilt University Medical Center

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.eq2ly4nxrlx9/v1>

Protocol Citation: IIDP-HIPP Brissova 2025. Analysis of Islet Function in Dynamic Cell Perifusion System. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.eq2ly4nxrlx9/v1>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: July 30, 2025

Last Modified: August 01, 2025

Protocol Integer ID: 223723

Keywords: IIDP-HIPP, Islet Function, Dynamic Cell Perifusion System, insulin, glucose, human islet phenotyping program, analysis of islet function, islet functional analysis, potency of the human islet preparation, human islet preparation, vanderbilt human islet phenotyping program, research in the integrated islet distribution program, islet function, integrated islet distribution program, purified human pancreatic islet product, insulin, hormone secretion measurement, method for quantitative determination, assay method, hipp procedure, dynamic cell perifusion system this standard operating procedure, diabetes, hipp procedure for dynamic perifusion, glucose stimulation, dynamic cell perifusion system, quantitative determination, goal of this sop

Abstract

This Standard Operating Procedure (SOP) is based on the Vanderbilt Human Islet Phenotyping Program (HIPP) Islet Functional Analysis. This SOP provides HIPP procedure for dynamic perifusion and hormone secretion measurement to assess islet function.

This SOP defines the assay method used by the Human Islet Phenotyping Program (HIPP) for qualitative determination of the Purified Human Pancreatic Islet product, post-shipment, manufactured for use in the National Institute of Diabetes and Digestive and Kidney Diseases (NIDDK)-sponsored research in the Integrated Islet Distribution Program (IIDP).

The goal of this SOP is to define the method for quantitative determination of insulin released after glucose stimulation for proving the potency of the human islet preparation shipped by the IIDP.

Guidelines

- **Integrated Islet Distribution Program (IIDP) (RRID:SCR_014387):** The IIDP is a grant funded program commissioned and funded by the NIDDK to provide quality human islets to the diabetes research community to advance scientific discoveries and translational medicine. The IIDP consists of the NIDDK Project Scientist and Program Official, the External Evaluation Committee and the CC at City of Hope (COH). The IIDP CC integrates an interactive group of academic laboratories including the subcontracted IIDP centers.
- **IIDP Coordinating Center (CC):** Joyce Niland, Ph.D. and Carmella Evans-Molina, M.D., Ph.D. serve as Co-Principal Investigators (Co-PIs) for the IIDP Program located within the Department of Diabetes and Cancer Discovery Science at COH to coordinate the activities of the IIDP and Human Islet Phenotyping Program (HIPP). Dr. Niland, contact PI, oversees the daily activity of the IIDP staff, provides informatics/ biostatistical input, and subcontracts with the Islet Isolation Centers (IICs) to ensure the delivery of the highest quality human islets to IIDP-approved investigators. Dr. Evans-Molina serves as the liaison to the HIPP, interacting closely to ensure that extensive, high quality phenotypic data are collected on islets distributed by the IICs. She also facilitates the delivery of this information to both the IICs and the IIDP-approved investigators, while responding to questions, issues, or suggestions for further HIPP enhancements.
- **Human Islet Phenotyping Program (HIPP):** The HIPP is a subcontracted entity of the IIDP through the COH and Vanderbilt University. The HIPP is directed by Marcela Brissova, Ph.D. and is responsible for performing specific standardized phenotyping assays agreed upon by both the IIDP and the HIPP, in order to provide enhanced, quality data on the human islets post-shipment, to the IIDP. The results of these assays will be approved by the CC and posted on the IIDP website for both the centers and the approved investigators.
- **Actual Islets (AI):** The actual number of islets counted.
- **Islet Equivalent (IEQ):** An islet with a diameter of 150 μm determined mathematically by compensating for islet shape.
- **Islet Perfusion Assay:** A functional assay that acquires dynamic hormone secretory profiles simultaneously from islet cell types such as β and α cells in response to their respective secretagogues. Insulin and glucagon are detected in perfusion fractions by hormone assay. The islet hormone secretory profile is generated by graphing hormone concentration over time with respect to islet volume and/or hormone content.
- **Secretagogue:** A substance that elicits another substance to be secreted.

References:



Citation

Dai C, Brissova M, Hang Y, Thompson C, Poffenberger G, Shostak A, et al
. Islet-enriched gene expression and glucose-induced insulin secretion in human and mouse islets.
Diabetologia. 2012 Mar;55(3):707–18. PMID: PMC3268985.

<https://pubmed.ncbi.nlm.nih.gov/22167125/>

LINK

Citation

Kayton NS, Poffenberger G, Henske J, Dai C, Thompson C, Aramandla R, et al
. Human islet preparations distributed for research exhibit a variety of insulin-secretory profiles.
Am J Physiol Endocrinol Metab. 2015 Apr 1;308(7):E592–602. PMID: PMC4385877.

<https://pubmed.ncbi.nlm.nih.gov/25648831/>

LINK

Materials

1. The following equipment is necessary to assess human islet function by perfusion.

- 1.1 Biological Safety Cabinet (BSC)
- 1.2

 37 °C

 CO² Incubator
- 1.3 Analytical Balance
- 1.4 Micropipettes (10-100 µL, 20-200 µL, and 100-1000 µL ranges)
- 1.5 Pipette Aid automatic pipettor
- 1.6 Light Microscope
- 1.7 Eyepiece with Calibrated reticle, 1 mm
- 1.8 Computer with Excel Counting Worksheet or equivalent
- 1.9 Manual or Electronic Cell Counter
- 1.10 Olympus SZX12 stereomicroscope equipped with an Olympus DP-80 high-resolution digital camera
- 1.11 Olympus cellSens™ image acquisition and analysis software
- 1.12 Perifusion System

Equipment	
A2 Biosafety Cabinet	NAME
4' REDISHIP Purifier Logic+	TYPE
Labconco	BRAND
302419000	SKU

Equipment

37° C CO2 Incubator	NAME
Thermo Fisher	BRAND
N/A	SKU
https://www.thermofisher.com/us/en/home/life-science/lab-equipment/lab-automation/lab-automation-incubators/37c-incubation.html	LINK



Equipment

Analytical Balance	NAME
Mettler Toledo	BRAND
https://www.fishersci.com/us/en/browse/90165061/analytical-balances	LINK



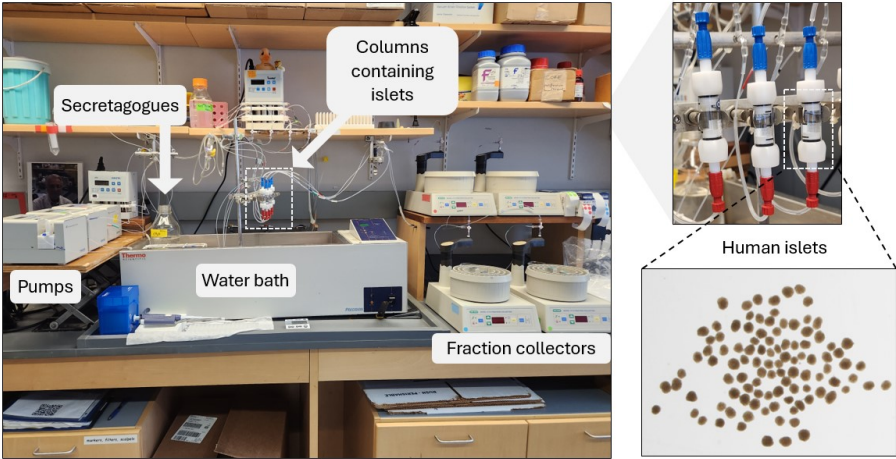
Equipment	
Micropipettes	NAME
10-100 µL, 20-200 µL, and 100-1000 µL	TYPE
Eppendorf	BRAND
N/A	SKU
https://www.pipette.com/Eppendorfpipettes?sbr=785&sbn=Eppendorf%20Research%20Plus%20Pipette	LINK
	

Equipment	
Drummond Pipette Aid automatic pipettor	NAME
Drummond	BRAND
N/A	SKU
https://www.drummondsci.com/pipette-controllers/portable-pipet-aid-xp/	LINK
	

Equipment	
Light Microscope	NAME
Leica Microsystems	BRAND
N/A	SKU
https://www.leica-microsystems.com/products/light-microscopes/ ^{LINK}	
	

Equipment	
Olympus DP-80 high resolution digital camera	NAME
DP80	TYPE
Olympus	BRAND
N/A	SKU
https://www.olympus-lifescience.com/en/camera/color/dp80/ ^{LINK}	
	

Perifusion System



Equipment

Circulating Water Bath

NAME

Thermo Fisher

BRAND

Model TSCIR89

SKU

<https://www.thermofisher.com/order/catalog/product/TSCIR89?SID=srch-srp-TSCIR89>^{LINK}



Equipment

Fraction Collector

NAME

Bio-Rad

BRAND

Model 2110

SKU

<https://www.bio-rad.com/en-us/product/model-2110-fraction-collector?ID=d67ac4c5-1f12-4a11-a1ea-ea6cebef40f8>

LIN
K



Equipment

Peristaltic Pump

NAME

Ismatec

BRAND

Model C.P. 78001-00

SKU

<https://www.avantorsciences.com/in/en/product/38670064/masterflex-ismatec-ipc-and-ipc-n-digital-multichannel-peristaltic-pumps-avantor>

LIN
K



2. The following supplies and materials are necessary for the perfusion system.

- 2.1  5-Luer Caps **Thomas Scientific Catalog #Western Analytical BC-125**
- 2.2  25 μ M Frits Polyethylene **Cole-Parmer Catalog #11945-04** Mfr # 006FR-10-25-PE
- 2.3  13 X 100 mm Polyethylene Tubes **Fisher Scientific Catalog #149567A**



- 2.4  Caps for 13 X 100 mm Tubes **Fisher Scientific Catalog #02681204**
- 2.5  Glucometer **Chemglass Life Sciences Catalog #CLS-1322-02**
- 2.6  Glucose Strips **Chemglass Life Sciences Catalog #CLS-1324-01**
- 2.7  Pressurized Gas, 95% O₂, 5% CO₂
- 2.8  Omnifit Glass Columns with 2 fixed-end pieces **Cole-Parmer Catalog #006CC-10-05-FF**
- 2.9  Tubing, 10 m, i.d. 1.0 mm, o.d. 1.8 mm, PE **GE Healthcare Catalog #19004001**
- 2.10  Tubing, Blue/Blue **Fisher Scientific Catalog #14-190-511**
- 2.11  Tubing, Orange/White **Fisher Scientific Catalog #14-190-503**
- 2.12  Teflon® Tubing, 1/16 OD X 0.03 ID X 20 **Fisher Scientific Catalog #50-453-683**

3. The following supplies and materials are necessary to culture human islets.

- 3.1  CMRL 1066 **Corning Catalog #15-110-CV**
- 3.2  L-Glutamine **Thermo Fisher Scientific Catalog #Invitrogen 25030-081**
- 3.3  Penicillin/Streptomycin **Thermo Fisher Scientific Catalog #Invitrogen 15140-122**
- 3.4  Fetal Bovine Serum **Merck MilliporeSigma (Sigma-Aldrich) Catalog #TMS-013-B**
- 3.5  P-1000 Pipet Tips Sterile **Thermo Fisher Scientific Catalog #2079-HR**
- 3.6  10 cm Petri Dishes SureGrip **Sycamore Life Sciences Catalog #Sarstedt 83.3902.500**
- 3.7  6 cm Petri Dishes SureGrip **Sycamore Life Sciences Catalog #Sarstedt 83.3901.500**

4. The following supplies and materials are necessary for the functional assessment of human islets.

- 4.1  Phosphate buffered saline (PBS) without Ca/Mg **Thermo Fisher Scientific Catalog #14190144**
- 4.2  1.5 mL Centrifuge Tube **Fisher Scientific Catalog #05-408-129**
- 4.3  NaHCO₃ **Merck MilliporeSigma (Sigma-Aldrich) Catalog #S6014-500G**
- 4.4  L-Glutamine **Merck MilliporeSigma (Sigma-Aldrich) Catalog #G8540-100G**
- 4.5  Sodium Pyruvate **Merck MilliporeSigma (Sigma-Aldrich) Catalog #P2256-25G**
- 4.6  HEPES **Merck MilliporeSigma (Sigma-Aldrich) Catalog #H7523-50G**
- 4.7  DMEM **Merck MilliporeSigma (Sigma-Aldrich) Catalog #D5030-10X1L**
- 4.8  RIA-grade BSA **Merck MilliporeSigma (Sigma-Aldrich) Catalog #A7888**
- 4.9  Glucose **Fisher Scientific Catalog #D1610**
- 4.10  Isobutylmethylxanthine (IBMX) **Merck MilliporeSigma (Sigma-Aldrich) Catalog #I5879**

- 4.11  KCl Merck MilliporeSigma (Sigma-Aldrich) Catalog #BP366-500
- 4.12  Epinephrine HCl Merck MilliporeSigma (Sigma-Aldrich) Catalog #E4642
- 4.13  Ascorbate (L-Ascorbic acid) Merck MilliporeSigma (Sigma-Aldrich) Catalog #A5960
- 4.14  200 µL (**P-200 ART**) and  1000 µL pipette tips (**P-1250 ART**)
- 4.15  0.22 µM Vacuum Filtration System Merck Millipore (EMD Millipore) Catalog #SCGPU05RE/S2GPU05RE

Safety warnings

Epinephrine HCl (Sigma E4642)

Safety information

 Epinephrine HCl (Sigma E4642) .pdf

Procedures

1 General Perifusion Startup

- 1.1 Fill water bath with deionized water to about 1 inch from the top, and set the temperature to  37 °C .
- 1.2 Label perifusion tubes with date, islet type, and any other identifying information, position the fraction collector trays for perifusion, and load the tubes.
- 1.3 Rinse the tubing with deionized water at max pump speed for 15 minutes, then place 2 new frits into the islet column.

2 Preparation of Base Perifusion Medium and Secretagogues

- 2.1 Prepare Base Perifusion Medium by combining compounds below in a 1-liter Erlenmeyer flask. Add 1L of deionized water and mix for at least 15 minutes until dissolved.

2.1.1  3.2 g NaHCO_3

2.1.2  0.58 L L-Glutamine

2.1.3  0.11 g Sodium Pyruvate

2.1.4  1.11 g HEPES

2.1.5 1 bottle DMEM for  1 L of media

2.1.6  1 g RIA-Grade BSA

- 2.2 Add  70 mg Ascorbate, Sigma A5960



2.3 Use a vacuum-filtration system to filter medium, transfer to a side-arm flask, and de-gas at 37 °C for at least 30 minutes.

2.4 Prepare 1.7 millimolar (mM) Low Glucose Medium (for use with Epinephrine)

2.5.1 Add 0.0613 g glucose to 200 mL Base perfusion medium to 250 mL bottle and mix until dissolved.

2.5.2 Wait 30 min and check glucose levels using a glucose meter.

2.5 Prepare 5.6 millimolar (mM) Glucose Medium

2.6.1 Add 0.5549 g glucose to 550 mL Base perfusion medium to 500 mL bottle and mix until dissolved.

2.6.2 Wait 30 min and check glucose levels using a glucose meter.

2.6.3 Reserve some 5.6mM medium in a 50 mL conical tube for islet loading and unloading.

2.6 Prepare 16.7 millimolar (mM) High Glucose Medium

2.7.1 Add 0.7522 g glucose to 250 mL Base perfusion medium to 250 mL bottle and mix until dissolved.

2.7.2 Wait 30 min and check glucose levels using a glucose meter.

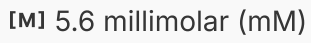
2.7 Prepare 1 micromolar (μM) Epinephrine/HCl in 1.7 mM Glucose (Store at -20 °C)

2.8.1 Prepare 200mM Epi Stock by adding 0.043934 g Epinephrine in 1 mL 1.7 mM Glucose medium.

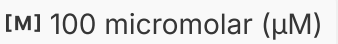
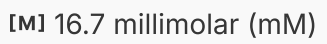
2.8.2 Prepare 0.4 mM Intermediate Epi Dilution (1:500): Dissolve 20 μL 200 mM Epi stock into 10 mL 1.7 mM Glucose medium.



2.8.3 Prepare 1 μ M Final Epi Concentration (1:400). Dissolve  250 μ L Intermediate Epi Dilution into  100 mL 1.7mM Glucose medium.

2.8 Prepare  20 millimolar (mM) KCl in  5.6 millimolar (mM) Glucose

2.9.1 Add  0.149 g KCl to  100 mL 5.6 mM Glucose medium.

2.9 Prepare  100 micromolar (μ M) IBMX in  16.7 millimolar (mM) Glucose

2.10.1 Weigh out around  10 mg IBMX

2.10.2 Make 100 mM Stock IBM by dividing the IBMX weight in mg by 22.22 and add "X" mL of Dimethyl Sulfoxide (DMSO). **Caution:** Wear gloves and avoid spilling any DMSO on skin.

2.10.3 Make 100 μ M Final Concentration (1:1000) by adding  100 μ L of Stock into  100 mL of 16.7 mM Glucose medium.

3 Setup of Secretagogues in Perifusion Water Bath

3.1 Place the bottles of media in the water bath to warm up for at least 10 minutes before beginning the perifusion.

3.2 Replace Pyrex orange caps with 4-Luer (+1) caps on every bottle of media to be used, and tape over the holes so that gas cannot escape.

3.3 Turn on gas (95% O₂, 5% CO₂), and place one gas catheter into each of the bottles with the new caps. Put tape over the other holes on the caps. Make sure that gas catheter is suspended above the media, and not inside the liquid.

3.4 Place intake catheters into baseline media bottle, making sure that they reach to the bottom of the bottle. Run the media through the chambers for about 10 minutes while islets are being aliquoted for perifusion, *disposing of media*.

4 Preparing Islets for Perifusion

- 4.1 Label a  1.5 mL clear Eppendorf tube for islets.
- 4.2 Place petri dish with cultured islets on the stage of inverted microscope and view the islets under 4x objective. Use the micrometer in the 10x eyepiece to size the islets.
- 4.3 Under microscope guidance handpick 267–300 IEQs using the chart below and record islet size and count in the *Islet Cell Calculation Excel Worksheet*.

Number of Tick marks	Diameter of Islet (μm)	Islet Size Category
>4 and ≤6	>100 and ≤150	Small
>6 and ≤8	>150 and ≤200	Medium
>8 and ≤10	>200 and ≤250	Large
>10 and ≤12	>250 and ≤300	Extra Large

- 4.4 Place the dish containing islets on the stage of a stereomicroscope equipped with a high-resolution camera and swirl until all islets are in the camera field of view at 10x magnification. Capture brightfield images at approximately 12-ms exposure and darkfield images at approximately 1.2-s exposure, each at 10x magnification. Ensure all islets are present in image. Save all image files.
- 4.5 Transfer all islets from the center of the dish to a labeled  1.5 mL clear Eppendorf tube for loading into perfusion chamber.
- 4.6 Open the darkfield image in the cellSens software. Using the manual HSV threshold function, segment the islet tissue channel.
- 4.7 Use the custom Count and Measure algorithm to determine islet count and mean islet diameter. Split adjacent but discrete islets using the Manually Split Objects tool to get an accurate islet count and mean diameter measurements.

Note

Note: Perfusion IEQ analysis by cellSens began 06/19/2018.

4.8 Use the mean diameter measurements to assign islets to a diameter group using the chart provided in 4.4 above.

5 Islet Loading into the Perifusion Chamber

5.1 Turn off the pump.

5.2 Turn the stopcock on the air bubble uptake lines so that the waste pathway is open, and close the stopcock on the outlet line at the fraction collector.

5.3 Remove the chamber from its mounting, turn it upside down, and remove the red end piece (inlet).

5.4 Remove and discard two thirds of the media from the chamber.

5.5 Using a  1 mL pipette set to  700 μ L transfer the slurry of islets from the Eppendorf tube to the perifusion chamber. Rinse the tube at least 3 times with  700 μ L of baseline medium, and transfer to the chamber.

5.6 Place the chamber back onto the mounting, and fill it up with baseline medium until there is a convex meniscus. Tap the sides of the column lightly to dislodge air bubbles from the walls, and collect and discard any bubbles from the top of the meniscus.

5.7 When all the bubbles have been removed, carefully replace the inlet plunger.
NOTE- make sure that no bubbles are introduced into the chamber during this process; if bubbles get in, remove the plunger, and repeat step 5.6.

5.8 Turn the chamber right side up, and put it back on the mounting rack.

5.9 Open the outlet line at the fraction collector and close the waste line on the air bubble uptake line.

5.10 Turn on the pump, set the fraction collector to 3 minutes and flip the collector arm so that it is over the first collection tube. Press "Start" on the collector to start the timer.

5.11 Tighten column end pieces, and end fittings, and make sure there are no leaks.

5.12 Lower the column mounting rack into the water bath, and tighten the clamp to prevent wobbling.

6 Fraction Collector Startup and Islet Wash Period

6.1 Collect 10 preliminary fractions to synchronize pump speed to deliver  3 mL 3-minute fractions, and to rinse the islets.

6.2 Record each pump speed on the perfusion worksheet in the perfusion logbook.

7 Collection of Perifusate Fractions

7.1 Begin to collect fractions.

7.2 Change secretagogues at predetermined fractions.

7.3 As soon as the fraction collector moves, switch off the pump.

7.4 Move the needles from one secretagogue to the other, making sure to not tangle the tubing, and inserting the needle all the way to the bottom of the bottle.

7.5 When the needles have been moved, restart the pump.

8 Recovery of Islets from Chamber and Perfusion System Cleanup

8.1 At the end of the perfusion raise the mounting rack from the water.

- 8.2 Stop pumps and close all outlets. Wait about 2 minutes, so that the islets can drift down to the bottom of the chamber.
- 8.3 Carefully remove the blue end piece (outlet) from the first chamber.
- 8.4 Pour the medium from the chamber into a 60 mm untreated dish, and rinse out the chamber and the blue column end piece with  1 mL of baseline medium 5 times each into the dish.
- 8.5 Rinse the chamber with deionized water, and remove frit with a frit removal stage and tool.
- 8.6 Put the column back together, and run 10% bleach through them at maximum pump speed for about 15 minutes. Make sure that all the tubing gets bleached, including all the waste lines. After the bleach, run deionized water through the entire system at maximum speed for 1.5 hours.
- 8.7 After the system cleanup, turn off the pump using the **master power switch**, turn off and drain the water bath, and log the perfusion in the notebook.
- 8.8 Perform any scheduled maintenance.

9 Islet Hormone Extraction

- 9.1 After retrieving islets from the perfusion chamber, size and count retrieved islets (step 4.4) to determine IEQ and *record islet size and count in Islet Cell Calculation Excel Worksheet*. Transfer the islets to a  1.5 mL Eppendorf tube.
Note: The average expected retrieval rate is approximately 93%. Lower islet recovery has been observed for islets showing hallmarks of disintegration.
- 9.2 Centrifuge the tubes for 3 minutes at 200 rcf and aspirate the supernatant using a pipette, being careful not to disturb the islet pellet.
- 9.3 Remove as much supernatant as possible from the tube, using a  1 mL tip, followed by a  200 μ L tip and place tube on ice.



- 9.4 Prepare fresh acid alcohol for hormone extraction by adding $50\ \mu\text{L}$ of concentrated HCl to $5.5\ \text{mL}$ of 95% ethanol.
- 9.5 Add $200\ \mu\text{L}$ of acid alcohol to tube containing islets.
- 9.6 Incubate sample at $4\ ^\circ\text{C}$ for 24 hours.
- 9.7 Spin samples down for 5 minutes at 3000 rcf, and transfer three $50\ \mu\text{L}$ supernatant aliquots into pre-labeled $2\ \text{mL}$ screwcap tubes and store at $-80\ ^\circ\text{C}$.

Data Storage and Reporting

- 10 To facilitate data management and ensure data security, the Vanderbilt HIPP uses an institutional server-based platform for data storage and analysis.

Deviations and Resolutions

- 11 Document any deviations that occurred during this protocol that affect the final results and report with the analysis of the assay.