


Dec 20, 2019

Version 2

🌐 Chloral Hydrate Seed Clearing V.2

DOI

dx.doi.org/10.17504/protocols.io.bakaicse

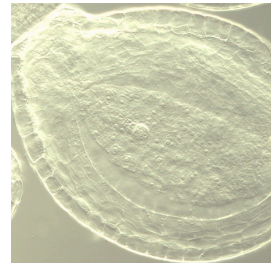

Gabrielle Sandstedt¹, Andrea Sweigart¹

¹University of Georgia

Mimulus



Gabrielle Sandstedt

Utah State University

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.bakaicse>

Protocol Citation: Gabrielle Sandstedt, Andrea Sweigart 2019. Chloral Hydrate Seed Clearing. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bakaicse>



License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: December 16, 2019

Last Modified: December 20, 2019

Protocol Integer ID: 31074

Keywords: Chloral Hydrate, Seed Clearing, Embryo, Endosperm, early seed development in mimulus, seeds with chloral hydrate, chloral hydrate seed clearing, early seed development, images of embryo, pollination, embryo, chloral hydrate, days after pollination, seed, differential interference contrast microscope, mimulus, endosperm development

Abstract

To characterize early seed development in Mimulus (1-5 days after pollination), we clear seeds with chloral hydrate and quickly obtain images of embryo and endosperm development using a Differential Interference Contrast Microscope.

Materials

Hoyer's Solution: 19g gum Arabic, 12g glycerol, 250g Chloral Hydrate, 75ml Water.

Diluted Hoyer's solution: 3 parts Hoyer's solution: 1 part 10% gum Arabic.

Forceps, glass microscope slides, coverslip, access to Differential Interference Contrast Microscope.

Troubleshooting

Safety warnings

- ❗ Chloral Hydrate is DEA regulated and you either need a permit or access to a lab that has a permit with this substance.

This substance is also hazardous and acutely toxic. SDS:





<https://www.caymanchem.com/msdss/21843m.pdf>




Before start

Prepare Hoyer's solution: 19g gum Arabic, 12g glycerol, 250g Chloral Hydrate, 75ml Water.



Prepare a diluted Hoyer's solution: 3 parts Hoyer's solution: 1 part 10% gum Arabic.

Prepare Hoyer's Solution

- 1 Hoyer's Solution:  19 g gum Arabic,  12 g glycerol,  250 g Chloral Hydrate,  75 mL Water.

 10 mL Diluted Hoyer's solution:  7.5 mL Hoyer's solution:  2.5 mL 10% gum Arabic

Dissecting fruit and clearing

- 2 Emasculate a bud from some maternal plant. 2 to 3 days later pollinate by selfing/outcrossing or use an unfertilized fruit.
- 3 Remove the developing fruit 1 to 5 days after pollination or 2 to 3 days after emasculation. In *Mimulus*, this protocol is useful for capturing early seed development (0- to 5 days). After 5 days, the seed tissue thickens and becomes difficult for viewing.
- 4 Pipette  10 μ L of diluted Hoyer's solution onto a glass slide and dissect developing ovules from the fruit directly onto the glass slide using sharp forceps
- 5 After dissection, pipette about  30 μ L of the diluted Hoyer's over the developing ovules and place a coverslip on top. Then, set the slide flat, upright in a 4°C fridge.
- 6 Depending on how developed the ovules are, clear for at least 1 to 12 hours before viewing with a Differential Interference Contrast Microscope