

Sep 08, 2020

COVIDscanDX mScanner Protocol (Confidential)

DOI

dx.doi.org/10.17504/protocols.io.bkyckxsw

don.cooper¹

¹Neuroganics

Mobile Assay

COVIDScanDx



don.cooper

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.bkyckxsw>

Protocol Citation: don.cooper 2020. COVIDscanDX mScanner Protocol (Confidential). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bkyckxsw>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: In development

We are still developing and optimizing this protocol



Created: September 07, 2020

Last Modified: September 08, 2020

Protocol Integer ID: 41700

Keywords: covidscandx mscanner protocol, protocol, confidential

Disclaimer

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to **protocols.io** is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with **protocols.io**, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.


Troubleshooting

COVIDscanDX mScanner Instructions for Use

1m

- 1 For iOS iPad and iPhone devices newer than 2018 use the following procedure. To scan the result of the developed test open the Safari browser on the iPad or iPhone device and navigate to www.register.covidscandx.com.
For Android Samsung Galaxy Tab A 10.1 inch (2019) use the following procedure. To scan the result of the developed test open the default browser on the tablet and navigate to www.register.covidscandx.com.
- 2 Enter the unique single-use product ID code included in the kit on the COVIDscanDX card. *Note the product code will expire after a valid image is automatically captured and a result is provided. Keep the COVIDscanDX card and product code to access the detailed results in the future to share with your electronic record of your result with your health care provider, travel entry authority, employer etc.

20s

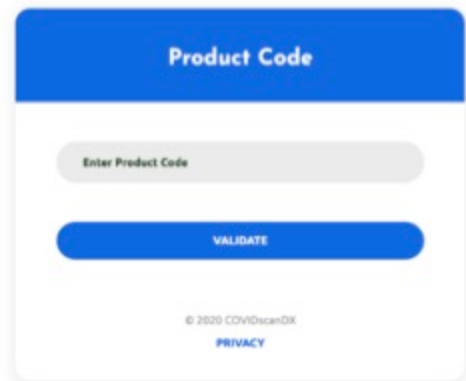


- 3 Following valid product code entry a new page is displayed on the device that gathers information about testing location and conditions and/or test subject symptoms. The user enters the information in the fields and provide a visual score of the test region signal ranging from 0-5 according to the COVIDscanDX card examples. A score of 0 indicates no test signal is visible while a score of 5 indicates the strongest test signal.
- 4 The application will ask for permission to open the device camera application. The user gives permission and positions the displayed blue box landmarks over the aruco targets and QR code until the blue boxes turn green and test image capture occurs within 1 second. If the lighting conditions are not optimal for image capture an error message will be shown that instructs the user to adjust the lighting or position of the test cartridge and image again.

10s

20s

Instructions

COVIDscanDX

1. Peel mounting tab and stick strip over template to the right. Place the test Control Line centered between the **Control Line** template brackets []. Do not cover black box targets or dashed lines at the TOP of template
2. Enter Product ID in device browser at register.covids candx.com
3. Score the Test line using the Test Line Score scale below.

4. Hold the camera over the test until the **blue** squares on the device screen turn **green** with the black boxes inside.

Product ID:x x x x x x

TOP
-----

Test Line



Control Line



Test Line Score

0 1 2 3 4 5

T**C**

Enter Product ID code at:
<https://register.covids candx.com>

LOT:xxxxxxx



- 5 Once the images have been captured the user is notified and server-side algorithms check the images for lighting conditions and position. A Positive, Negative or Invalid result will appear on the device display within 5 seconds (assuming web connectivity). Depending on the result the user is given instructions on further support to establish the diagnosis of COVID-19 via telemedicine portal.
- 6 Post test results are accessible for review and integration via JSON or CSV format to partner database by site approved users through the partner dashboard.

10s