

Dec 22, 2022

CZI Pediatric Nasopharyngeal Swab Processing for 10X scRNA-seq (Illustrated Protocol)

 Forked from [Human Nasopharyngeal Swab Processing for Viable Single-Cell Suspension](#)

DOI

<https://dx.doi.org/10.17504/protocols.io.6qpvr49jogmk/v1>

Jaclyn Long¹, erica.langan¹, Ying Tang¹, Carly G.K. Ziegler², Vincent N. Miao², Andrew W. Navia², Joshua D. Bromley², Kenneth J. Wilson³, Yilianys Pride³, Mohammad Hasan³, Taylor Christian³, Hannah Laird³, Anna Owings³, Meredith Sloan³, Haley B. Williams³, Tanya O. Robinson³, George E. Abraham III³, Michal Senitko³, Sarah C. Glover³, Bruce Horwitz¹, Alex K. Shalek², Jose Ordovas-Montanes¹

¹Boston Children's Hospital; ²Massachusetts Institute of Technology; ³University of Mississippi Medical Center

Human Cell Atlas Metho...

Coronavirus Method De...



Erica M Langan

Boston Children's Hospital

Create & collaborate more with a free account


Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.6qpvr49jogmk/v1>

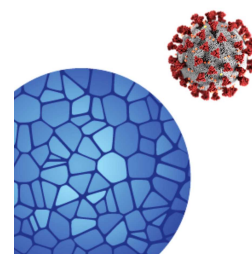



Protocol Citation: Jaclyn Long, erica.langan , Ying Tang, Carly G.K. Ziegler, Vincent N. Miao, Andrew W. Navia, Joshua D. Bromley, Kenneth J. Wilson, Yilianys Pride, Mohammad Hasan, Taylor Christian, Hannah Laird, Anna Owings, Meredith Sloan, Haley B. Williams, Tanya O. Robinson, George E. Abraham III, Michal Senitko, Sarah C. Glover, Bruce Horwitz, Alex K. Shalek, Jose Ordovas-Montanes 2022. CZI Pediatric Nasopharyngeal Swab Processing for 10X scRNA-seq (Illustrated Protocol). **protocols.io** <https://dx.doi.org/10.17504/protocols.io.6qpvr49jogmk/v1>

License: This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: December 19, 2022

Last Modified: December 22, 2022

Protocol Integer ID: 74234

Keywords: Nasal Swab, COVID-19, Nasopharyngeal Swab, single-cell RNA-seq, nasal epithelia, nasopharynx, sars-cov-2, cryopreserved human nasopharyngeal swab, czi pediatric nasopharyngeal swab, recovering viable single cell suspension, individual bulk rna, cell rna, 10x scrna, rna, swab, seq, cell


Disclaimer

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to **protocols.io** is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with **protocols.io**, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

Abstract

A protocol for recovering viable single cell suspensions from cryopreserved human nasopharyngeal swabs for pooled single-cell RNA-seq and individual bulk RNA-seq. The illustrated schematic below details the process.

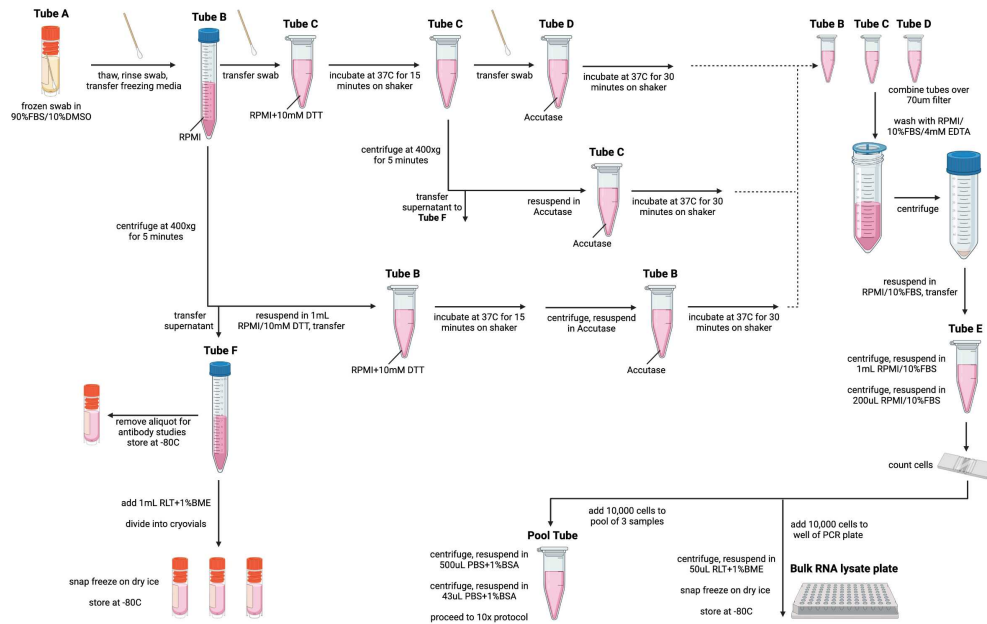

Attachments



Cell Isolation from ...

2MB

Image Attribution

Created with BioRender

Guidelines

Samples should be collected by a trained medical professional using a nasal swab (FLOQSwabs, Copan flocked swabs) in accordance with the manufacturer's instructions. Briefly, the process was performed as follows. First, the patient's head was tilted back slightly, and the swab was inserted along the nasal septum, above the floor of the nasal passage to the nasopharynx until slight resistance was felt. The swab was then left in place for several seconds to absorb secretions and slowly removed while rotating. The swab was placed in a 1.7 mL cryovial containing 90% fetal bovine serum (FBS)/10% dimethyl sulfoxide (DMSO) and frozen using a slow-cooling device (Thermo Scientific Mr. Frosty Freezing Container) at -80 °C, and stored in liquid nitrogen.

Materials

Materials & Reagents for 1 swab (with overages):

- 2.5 mL of RPMI + 10 mM Dithiothreitol (DTT) (made fresh)
- 3.5 mL of Accutase
- 6 mL RPMI
- 8 mL quenching buffer containing RPMI + 10% fetal bovine serum (FBS) + 4 mM EDTA
- 2 mL RPMI + 10% FBS
- 15 mL conical labeled **Tube B** containing 2 mL RPMI (previously 5 mL)
- 1.5 mL tube labeled **Tube B**, empty
- 1.5 mL tube labeled **Tube C**, with 1 mL RPMI + 10 mM DTT
- 1.5 mL tube labeled **Tube D**, with 1 mL Accutase
- 50 mL conical
- 70 µm cell strainer that fits 50 mL conical
- 1.5 mL tube labeled **Tube E**
- 1.5 mL tube for counting
- 1.5 mL tube for pooling cells (if doing 3 samples, need 1 tube)
- 10 µL trypan blue
- 1 mL PBS.+ 1% BSA
- 15 mL conical for collecting viral supernatant labeled **Tube F**
- 3 cryovials for viral supernatant lysate
- 3 96-well PCR plates or 3 cryovials for bulk RNAseq lysate
- RLT + 1% 2-mercaptoethanol (BME)

Equipment

- Forceps and scissors
- Thermomixer set to 37°C, agitating at 300 rpm
- Hemocytometer slides

Safety warnings

- ❗ For hazard information and safety warnings regarding nasopharyngeal swabs or any listed materials, please refer to the SDS (Safety Data Sheet).

For samples obtained from individuals diagnosed with, or at any risk of, an infection (e.g. SARS-CoV-2), additional precautions should be taken in accordance with your institute's regulations on biosafety. These include elimination of aerosol generating steps where possible (e.g., no vacuum aspiration), all steps prior to cell lysis should be carried out in a biosafety cabinet, including thermomixing and centrifugation where possible. When necessary, samples should only be removed from the biosafety cabinet in decontaminated and sealed secondary containment. Personal protective equipment including a gown, two pairs of non-sterile gloves, a protective surgical or N95 mask, and a face shield should be worn during sample processing.


Before You Start

1 Record sample characteristics in Table 1


	Site ID	Study ID	Case/Control	Virus	Date of Enrollment	Sex	Age	Sample ID	Batch	Pool

2 Prepare 50 mL conicals with necessary reagents:

- 1X RPMI + 10 mM DTT
- 1X Accutase
- RPMI + 10% FBS + 4 mM EDTA (quenching buffer)
- RPMI + 10% FBS

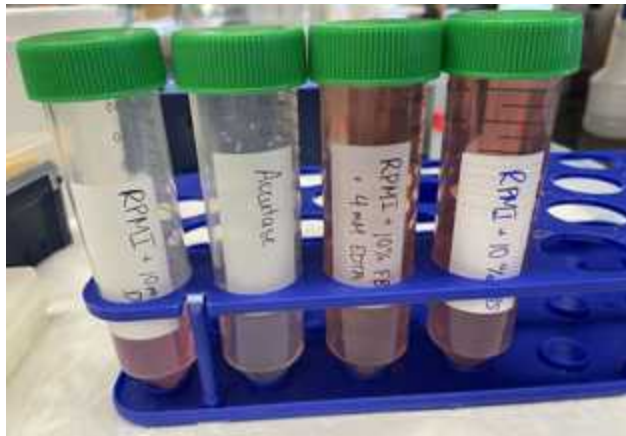

3 Label 2 sets of 15 mL conicals (**Tube B & Tube F**). Prepare **Tube B** with 2 mL RPMI .

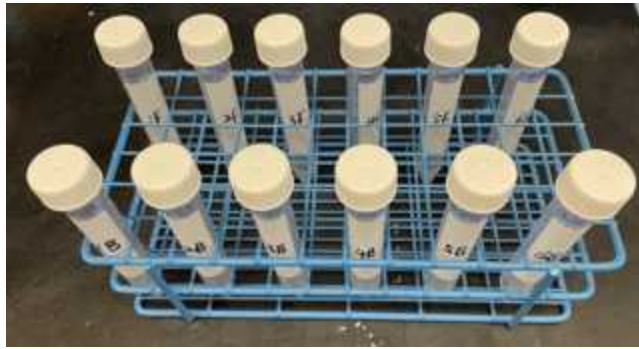


Tube B and Tube F setup for 6 swabs

- 4 Label 4 sets of 1.5 mL tubes (**Tube B, Tube C, Tube D, Tube E**). Prepare **Tube C** with  1 mL RPMI + 10 mM DTT and **Tube D** with  1 mL Accutase



- 5 Label 1 set of 50 mL tubes for filtering at later step (can set aside for beginning of protocol)

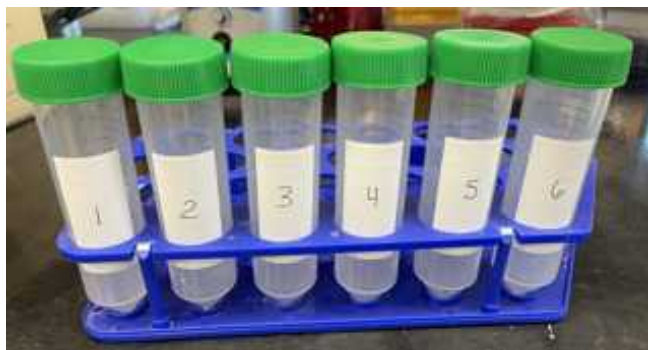

Final Set-Up

6



Timing

7 Fill out timing chart at indicated steps of protocol

Step	Time
Start time:	
Time after preparation of single cell suspension:	
Time after counting:	
Time 10X controller run started:	
Time 2nd 10X controller run started (if applicable):	
Time lysis buffer added for bulk lysates:	



Tube A (the cryovial)

- 8 Rapidly (within 1-2 minutes) thaw cryovial (**Tube A**) in thermal block set to  37 °C .

Safety information

Carefully perform in accordance with your institute's safety guidelines. If handling potentially infectious material, inspect for cracks or leaks during warming

- 9 Remove swab from **Tube A** using clean forceps.



Step 9-11 - first swab transfer.mp4

- 10 Dip swab in **Tube B** (15 mL conical) once. Swirl briefly (approximately 3 seconds) to rinse swab.

- 11 Move swab from **Tube B** (15 mL conical) to **Tube C**. Trim swab handle using scissors if necessary so that it can fit inside a 1.5mL tube (see video). Proceed directly to step 27 for **Tube C**. In parallel, continue processing **Tube B**.

- 12 Transfer liquid in **Tube A** to **Tube B** (15 mL conical).



- 13 Using  1 mL RPMI from **Tube B** (15 mL conical), wash **Tube A**.



- 13.1 Collect washing from **Tube A** in **Tube B** (15 mL conical). Continue below with **Tube B**

- 14 Discard **Tube A**.

Tube B

- 15 Centrifuge **Tube B** (15 mL):  400 x g, 4°C, 00:05:00 .



- 16 Remove supernatant with serological pipette or *P1000*. Do **not** discard supernatant.



- 17 Transfer supernatant to **Tube F**.
- 18 Resuspend pellet in  1 mL RPMI + 10mM DTT .
- 19 Transfer suspended cells from **Tube B** (15 mL) to **Tube B** (1.5 mL). Discard empty 15 mL conical. 
- 20 Place **Tube B** (1.5 mL) on thermomixer (37°C, 300 rpm).
- 21 Incubate for  00:15:00 . In parallel, continue processing **Tube C**. 
- 22 Centrifuge **Tube B** (1.5 mL):  400 x g, 4°C, 00:05:00 . 
- 23 Remove supernatant with P1000 pipette. Discard supernatant. 
- 24 Resuspend pellet in  1 mL Accutase .
- 25 Place **Tube B** (1.5 mL) on thermomixer (37°C, 300 rpm).
- 26 Incubate for  00:30:00 . 

Tube C

- 27 Place **Tube C** on thermomixer (37°C, 300 rpm).
- 28 Incubate for  00:15:00 . 



29 Place swab in **Tube D**. Proceed directly to step 36 for **Tube D**.



Step 29 - second swab transfer.mp4

30 Centrifuge remaining liquid at  400 x g, 4°C, 00:05:00 .



31 Remove supernatant with P1000 pipette. Do not discard supernatant.



32 Place supernatant in **Tube F**. Keep **Tube F** on ice until you proceed to Step 71.

33 Resuspend pellet in  1 mL Accutase .

34 Place **Tube C** on thermomixer (37°C, 300 rpm).

35 Incubate for  00:30:00 .



Tube D

36 Place **Tube D** on thermomixer (37°C, 300 rpm).

37 Incubate for  00:30:00 .



Note

Take 10X reagents out at this step so that they have time to equilibrate to the appropriate temperature.

After **Tube B**, **C**, and **D**'s 30 minute Incubations

38 After **Tube B**, **Tube C**, and **Tube D** have each finished their 30 minute incubations:

**Note**

In practice, we wait until all tubes have finished their 30 minute incubation in Accutase to synchronize. We leave tubes on incubation for a maximum 50 minutes.

39 Place 70 μ m filter in 50 mL conical.

40 Wet filter with  3 mL quenching buffer (RPMI + 10% FBS + 4 mM EDTA) .

41 Pipette contents of **Tube B**, **Tube C**, and **Tube D** onto filter with P1000 pipette

 Step 41 - adding Tubes B,C,D conte...

**Note**

Do not discard original tubes.

42 Use  1 mL quenching buffer to wash each **Tube B**, **Tube C**, and **Tube D**.



42.1 Manually agitate the swab in **Tube D** in the quenching media with forceps to ensure full rinse.

43 Add quenching buffer from washes to filter. Discard **Tubes B**, **C**, and **D**.



44 Wash filter with additional  2 mL quenching buffer .



45 Discard filter, cap 50 mL conical.

46 Centrifuge 50 mL conical  400 x g, 4°C, 00:10:00 .





Note

It will likely be very challenging to see a pellet in the 50 mL conical at this point

- 47 Remove supernatant with serological pipette. Leave ~500 μ L in the bottom of the tube.



Approximate residual volume to aim for

Note

Remove supernatant carefully at this step! You don't want to pipette up your cells but you will need to add the volume you leave in this 50 mL conical and 1 mL RPMI + 10% FBS to a 1.5 mL tube in the next steps, so try not to exceed 500 μ L residual volume.

48 Add 500 μ L RPMI + 10% FBS to the tube to resuspend cells.



49 Transfer resuspended cells (~ 1 mL) from 50 mL conical to **Tube E** (1.5 mL tube).

50 Wash 50 mL conical with  500 μ L RPMI + 10% FBS .



50.1 Transfer washing from 50 mL conical to **Tube E**.



51 Centrifuge **Tube E**  400 x g, 4°C, 00:05:00 .



52 Remove supernatant with P1000 pipette.



53 Resuspend pellet in  1 mL RPMI + 10% FBS .



Example pellet in Tube E after centrifuging

Note

At this point, you should be able to see a reasonably-sized pellet!

54 Centrifuge **Tube E**  400 x g, 4°C, 00:05:00

5m

55 Resuspend pellet in 200 μ L RPMI + 10% FBS

Count cells from **Tube E**

56 In a new empty 1.5 mL tube, add  10 μ L trypan blue .



57 Add  10 μ L cells from Tube E to 1.5 mL tube containing trypan blue.



58 Pipette to mix cells in trypan blue, transfer  10 μ L to hemocytometer port.



59 Count viable cells across 8 quadrants. One quadrant is the corner of a hemocytometer grid, consisting of 16 smaller squares.




60 Record total cell number and calculate cell concentration using the "Cell Counts & 10X Volumes" sheet of the Calculation Tables template.



Calculation Tables template.xlsx

61 Take photo of cells at 4X. This can be done using a camera attached to the microscope or with a phone placed adjacent to the microscope lens.

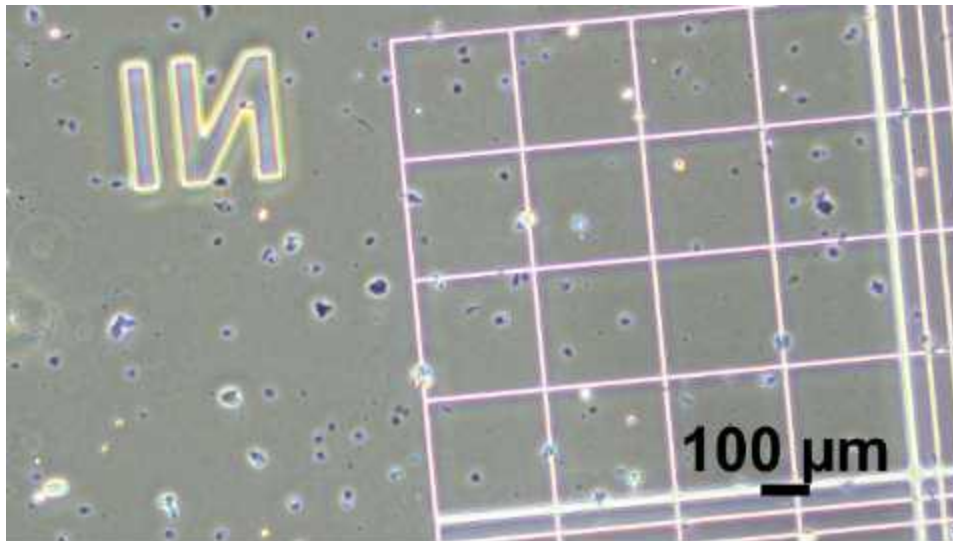

Representative hemocytometer image at 4X

Prepare cells for 10X

62 Prepare 1.5 mL tube for each pool

63 For each sample, add "Volume Necessary for Desired Number of Cells" from "Cell Counts & 10X Volumes" sheet of the Calculation Tables template to correct pool.

64 Pipette cells to mix

65 If volume of pool is $> 100 \mu\text{L}$, centrifuge cells  $400 \times g$, 4°C , 00:05:00 , then resuspend in  $500 \mu\text{L}$ PBS + 1% BSA and proceed to next step.

5m

If volume of pool is $< 100 \mu\text{L}$, add  $500 \mu\text{L}$ PBS + 1% BSA directly to pool and proceed to next step.

66 Centrifuge cells  $400 \times g$, 4°C , 00:05:00

5m

67 Resuspend cells in  $43.3 \mu\text{L}$ PBS + 1% BSA (Total volume of cell suspension + water on page 27 of 10X Protocol).

68 Add  $43.3 \mu\text{L}$ cell suspension directly to  $31.9 \mu\text{L}$ master mix (Step 1.2b of 10X Protocol)

69 Proceed with instructions in 10X Protocol to load the chip and run the controller (through Step 1.3)

70 At Step 1.4f, take a picture of the GEMs in the pipette tips.



Example photo of successful run

Tube F Processing (Viral Lysates)

71 Remove a 200uL aliquot from **Tube F** and transfer to a new, empty 1.5mL tube. Store at -80°C for antibody studies.

72

**Note**

In **Tube F**, you should have 2 mL from original RPMI + 0.5-1 mL from swab cryopreservative + 1 mL from **Tube C** supernatant = ~4 mL total

Add  1 mL RLT + 1% BME to **Tube F**.

73 Distribute contents of **Tube F** into 3 cryovials per sample

Note

You can save time here by labeling cryovials/plates for viral and bulk lysates the day before!

74 Snap freeze on dry ice

75 Store at  -80 °C

Lysates for Bulk RNA-seq

40m

76 For each lysate, add "Volume Necessary for Each Lysate" from "Bulk RNA Lysate Volumes" sheet of the Calculation Tables template to one well of a 96-well PCR plate according to the Lysate Storage Plate map.

Note

Typically, we aim to store 3 lysates per sample. To do so, we label 3 plates and fill 1 well per sample in each plate, if cell numbers allow. Make sure to record which well has which sample.

77 Seal plates and centrifuge  400 x g, 4°C, 00:05:00

5m

78 Aspirate media

**Note**

It may be difficult to see a pellet here, so it's okay to leave some residual volume before resuspending in lysis buffer in the next step

- 79 Resuspend cells in  50 μ L RLT + 1% BME . Mix with pipette. Bubbles are okay here.
- 80 Seal plate with foil seal and spin down briefly.
- 81 Place plate on dry ice for  00:15:00 -  00:20:00 to snap-freeze lysate.
- 82 Store at  -80 °C until ready to use.

35m