


Nov 25, 2024

Version 3

Illumina library preparation and dual hybridization protocol of ERC GLOBAL V.3

DOI

<https://dx.doi.org/10.17504/protocols.io.n92ldp787l5b/v3>

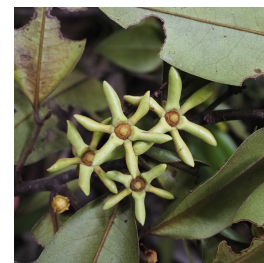

Vincent R C Soulé¹, Theo Delaigue¹, Cedric Mariac¹, Couvreur Thomas¹

¹DIADÉ, Univ Montpellier, CIRAD, IRD, Montpellier, France



Couvreur LP Thomas

Institut de Recherche pour le Développement (IRD)

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.n92ldp787l5b/v3>

Protocol Citation: Vincent R C Soulé, Theo Delaigue, Cedric Mariac, Couvreur Thomas 2024. Illumina library preparation and dual hybridization protocol of ERC GLOBAL. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.n92ldp787l5b/v3> Version created by **Couvreur LP Thomas**

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: November 25, 2024

Last Modified: November 25, 2024

Protocol Integer ID: 112712

Keywords: library preperation, herbarium specimens, Annonaceae, hybridization, Angiosperm353, Small fragment DNA, illumina, novaseq, annonaceae dna extraction protocol, annonaceae dna extraction protocol from silicagel, annonaceae species within the erc global project, daicel arbor biosciences mybaits angiosperm353 kit, dna extraction of plant, annonaceae species, daicel arbor biosciences mybait, silicagel specimens of annonacea, dna extraction, illumina library preparation, custom annonacea, particular from herbarium, dual hybridization protocol of erc global step, hybridization capture, silicagel specimen, diade research unit ird montpellier france, hybridisation protocol, annonacea, herbarium, dual hybridization protocol, step protocol for library preparation, hybridization capture for targeted ngs user manual version, steps for the dual hybridization, extraction, purification, baiting kit, plant, silicagel, illumina, extraction protocol

Funders Acknowledgements:

ERC Consolidator

Grant ID: 865787

Disclaimer

This project has received funding from the European Research Council (ERC) under the European Union's Horizon 2020 research and innovation program (grant agreement No. 865787).

Bill Baker and Donovan Bailey are thanked for drawing our attention to this option and sharing preliminary data and results.

Abstract

Step by step protocol for library preparation of Annonaceae species within the ERC GLOBAL project. This protocol was followed for all libraries prepared within the project. Each step involves a purification phase (AMPure XP protocol), and is detailed at the start of this protocol. The protocol is done per series of 48 samples at a time with indexing.

The protocol also describes the steps for the dual hybridization of Daicel Arbor Biosciences Mybaits Angiosperm353 kit + Custom Annonacea baiting kit. The protocol is based on hybridisation protocols used at DIADE research unit IRD Montpellier France and on the Hybridization Capture for Targeted NGS User Manual version 5 protocol "standard" of Daicel Arbor Biosciences Mybaits.

To proceed DNA extraction of plant and in particular from herbarium and silicagel specimens of annonacea read the Annonaceae DNA extraction protocol from silicagel dried and herbarium preserved leaves

- Vincent Soulé¹,
- Théo Delaigue¹,
- **Thomas LP Couvreur**¹,
- **Cedric Mariac**¹
- ¹DIADE, Univ Montpellier, CIRAD, IRD, Montpellier, France
- [dx.doi.org/10.17504/protocols.io.5qpvorqx9v4o/v1](https://doi.org/10.17504/protocols.io.5qpvorqx9v4o/v1)

Attachments



[TAG ADAPT INDEX](#)

[sequ...](#)

5KB

Image Attribution

Photo by Thomas L.P. Couvreur

Guidelines

This library preparation & dual hybridisation protocol follows and is inspired by Henriks et al (2021) using both the Annonaceae specific baits and the Angiosperms353 baits (see references). In the end we sequence both the Annonaceae specific baits and the Angiosperms353 universal ones.

Both the Annonaceae (Couvreur et al. 2019) and Angiosperms353 (Johnson et al., 2019) bait kits are available as Arbor Biosciences "myBaits" kits (Daicel Arbor Biosciences, Ann Arbor, Michigan, USA; <https://arborbiosci.com/genomics/targeted-sequencing/mybaits/>).



Materials

- Pipettes and cones with "low bind" filters (1000-200-20-10µL)
- Eppendorf tubes 1.5mL-2mL-0.2mL, aluminum films, refrigerated holder, magnetic holder, centrifuge and vortex
- PLATES RT-PCR 96 (Roche)
- VWR® , Centrifugeuse de microplaques, PCR Plate Spinner
- Lightcycler 480 (Roche)
- QIAxcel "Advanced system" analyze + QIAxcel DNA screening kit (2400) or other bioanalyser
- Personal protective equipment (PPE): gown, gloves, goggles
- DNA Dilution Buffer
- Agencourt AMPure XP solution
- 70% ethanol (fresh)
- DNase/RNase Free Water
- Thermo Scientific™ ADN ligase T4 (5 U/µl) EL0013 Code produit.10548730
- Bst DNA Polymerase, Large Fragment (NEB#: M0275)
- DNTP Mix , REF U151B 10mM U1515 or 1511 (200 or 1000µL)
- KAPA HiFi HS Real-Time Master Mix (2X) + Real-Time Standards 1 to 4
- Kit MYbaits® "In-Solution Sequence Capture for Targeted High-Throughput Sequencing".

Protocol materials

⊗ Agencourt AMPure XP **Beckman Coulter Catalog #A63880**

⊗ 70% EtOH

⊗ PCR Grade Water **Catalog #AM9935**

Safety warnings

! Wear protections

Before start

Clean bench and pipettes with water and then ethanol. Follow the instructions for use of the Arbor Biosciences V5 kit.

check if you have all the material needed for all the step or equivalent tools



Introduction and AMPure XP protocol

30m

- 1 Library preparation protocol is carried out after the DNA extraction detailed in :

Protocol



NAME

Annonaceae DNA extraction protocol from silicagel dried and herbarium preserved leaves

CREATED BY

Couvreur LP Thomas

[Preview](#)

To start we need purified DNA diluted in water or TE buffer 50µL or more at 2 to 900ng/µL

below 60ng/µL samples are considered low quantity if only diluted in 50µL

samples with size less than 200pb is considered low quality/small fragments dna (use 1.8X ampure)

2 **AMPURE XP DNA Purification protocol**

On each steps of the protocol a purification step is need to change solution, concentrate or purify dna.

Each purification step is done with



Agencourt AMPure XP **Beckman Coulter Catalog #A63880**

The only changes between the different purification steps is the start volumes as such: volume ratio 1X to 1.8X AMPure. For example for 1.8X: for 100 µl to be purified added 180 µl Ampure

Below we provide a protocol for 1X



Agencourt AMPure XP **Beckman Coulter Catalog #A63880**

for start volume of  50 µL of DNA

resuspension (elution) volume =  20 µL .

The protocol is described for 96 well pcr plates, but can also be done in 1.5/ 2mL tubes.



- 2.1 Use  Agencourt AMPure XP **Beckman Coulter Catalog #A63880** at  Room temperature , take out 30 minutes in advance.
Prepare a 70% ethanol solution.
- Add  50 μ L AMPure XP in each well.
- 2.2 Mix by pipetting or shake and wait  00:10:00 min. quickly spin down 10sec using pcr plate spinner if pcr plate or mini centrifuge for tubes
DNA will bind to the beads 
- 2.3 Place the plate on the magnetic stand and wait  00:05:00 min.
Remove liquid phase by pipetting , **avoid pipetting the beads.**  
- 2.4 Keep on magnetic stand add  150 μ L  70% EtOH . 
Remove liquid phase by pipetting, **avoid pipetting the beads.** 
- 2.5 Repeat wash step on magnetic stand and add  150 μ L  70% EtOH 
Remove liquid phase by pipetting, **avoid pipetting the beads.** 
Allow to dry for 5 minutes at  Room temperature
- 2.6 Take PCR plate out of the magnetic stand and add the needed "elution volume"; here  20 μ L of ultra pure water  PCR Grade Water **Catalog #AM9935** 
Mix by pipetting.   
- Wait  00:10:00 minutes at  Room temperature for beads to resuspend in water.
- 2.7 The AMPure purification step is over.
Make sure to place the plate on the magnetic stand to avoid pipetting beads for next usage or quantification.
Mind that purification step reduces DNA quantity by 20 to 50% and removes small DNA fragments.
To keep small fragments use 1.8X AMPure ratio instead of 1X.



Part 1: library preparation , DNA preparation and shearing

2h

- 3 Depending on concentration of each sample, dilute the DNA in  0.2 mL microtubes Diagenode suitable for shearing with ultrasound. Target amount of DNA should, if possible, be  3-4 µg in  50 µL (to be supplemented with water). The number of shearing cycles and their duration can be modified depending on the initial degradation state of the DNA and/or the desired fragment size. Shearing is done with the below equipment:

Equipment

Bioruptor Pico sonication device	NAME
Sonicator	TYPE
Diagenode	BRAND
B01060010	SKU

- 3.1 Turn ON the Bioruptor Pico  00:30:00 before use to cool down water at  4 °C 30m
- 3.2 Pulse centrifuge the tubes. 1m
- 3.3 Deposit  50 µL of DNA and water if needed to dilute to match concentration to the tubes in the Pico Bioruptor (Diagenode).
- 3.4 Start N cycles  00:00:15 ON /  00:00:30 OFF at  4 °C . 45s
Depending on the DNA extraction control gel $0 < N < 10$.



Expected result

For low quality small fragments herbarium specimen 0-3 cycles
For high quality herbarium specimen 3-8 cycles
For Silicagel specimens 5-10 cycles

- 3.5 Check size using  3 μL DNA + 10 μL DNA dilution buffer AL420 protocol with QIAxCEL, if necessary undertake additional cycles.

Pattern must be between 150 and 600 bp, but depending on the initial degradation, it may be centred on 100 bp.

QIAxcel QIAGEN (DNA Screening Kit)

- 3.6 Restart N+B cycles  00:00:15 ON /  00:00:30 OFF at  4 $^{\circ}\text{C}$.

45s

*

- 3.7 Measure DNA size again using QIAGEN QIAxCEL Bioanalyser
 3 μL DNA + 10 μL DNA dilution buffer AL420 protocol analysis.

*

or other bioanalyser , you can also make a DNA gel

Expected result

We aim for 4 μg DNA post shearing and optimal DNA size of 400pb ranging between 100-600pb.

volume after shearing step 40-50 μL (40 needed for next step)



Fast End Repair & phosphorylation

1h

- 4 Prepare a buffer and enzyme mix according to the number of samples to be processed

A	B
For 1 sample:	Volume
Sheared DNA	40 μ L
End Repair Reaction Buffer 10X	5 μ L
End Repair Enzyme Mix	5 μ L
Final volume	50 μ L

Buffer mix needed per sample

NEBNext End Repair Module E6050L or E6050S

- 4.1 Distribute  10 μ L of enzyme and buffer mix in each well. 5m
- 4.2 Add  40 μ L of DNA from previous step, vortex and Spin down the sample (using pcr plate spinner). 5m
- 4.3 Leave  00:30:00 at  20 $^{\circ}$ C on the thermocycler. 30m
- 4.4 Stop the reaction using a purification step AMPure XP 1X (50 μ L) or 1.8X (90 μ L) if you want to retain small DNA fragments. See introduction for details of the purification step using  Agencourt AMPure XP **Beckman Coulter Catalog #A63880**

4.5 Elute the samples in  30 µL water and quantify using Nanoquant.



Expected result

optional if low quantity DNA elute on 16µl (2µL for nanoquant and 14 for next step)

Ligation

3h

5 Depending on the previous dosage, calculate the volume to be taken per sample to ideally obtain  200 ng of DNA in  13.8 µL .

If necessary, dilute samples that are too concentrated .

Careful not to vortex the ligase. Mix by flicking the tube.

5.1 Take the samples, place them on a 96-well plate on a refrigerated stand, and adjust all the volumes to  13 µL with water.

5.2 Then add to each sample  2 µL of TAG-P5 adapt (4 µM).

5.3 Prepare a volume of Mix according to the number of samples:

A	B
1 sample (DNA = 200 ng)	
DNA + water	13.8 µL
adaptateur P5 with tag 4µM	2 µL
adaptateur P7 40µM	0.2 µL
T4 DNA ligase 5U/µL	2 µL
Buffer T4 10X	2 µL
Final Volume	20 µL

Mix volume for 1 sample



Update the thawing/freezing number of P7 P5 adapters on their storage box.

NB; Tag refers as P5 6pb barcoded adapter used as index but not read by illumina ; demultiplexed after sequencing.

adapter p7 & p5 sequences on attached file

Thermo Scientific™ ADN ligase T4 (5 U/μl) EL0013

5.4 Deposit  4.2 μL of Mix by well, vortex briefly and pulse centrifuge.

5.5 Place the plate on the thermocycler and launch the "LIG" program:

	A	B	C
	Step	Temperature	Duration
	1	+22°C	2H30
	2	+65°C	10 min
	3 optional	+4°C	20 min

LIG program on thermocycler

5.6 Purify 20μL post ligation using

 Agencourt AMPure XP Beckman Coulter Catalog #A63880 with XP1.8X ratio

( 36 μL of ampure);

Then elute in  25 μL ultra pure water.



**Expected result**

25µL of post ligation purified DNA solution 23.5 is needed for next step

optional you can conserve for 24h sample after purification step at -20°C

Nick Fill-in (Bst DNA Polymerase)

45m

- 6 Position the 96-well plate on a refrigerated stand.
Prepare a volume of Mix according to the number of samples:

5m

A	B
For 1 sample:	
DNA	23.5 µL
DNTPs 5mM	1.5 µL
Bst DNA POLYMERASE 8U/ µL	2 µL
Thermo Pol Buffer 10X	3 µL
Final Volume	30 µL

Mix for 1 sample

Bst DNA Polymerase, Large Fragment
(NEB#: M0275)

DNTP Mix , REF U151B 10mM U1515 or 1511 (200 or 1000µL)

- 6.1 Deposit  6.5 µL of Mix per well in a AB1400 pcr plate and pulse centrifuge.

1m

- 6.2 Deposit  23.5 µL of DNA per well, mix with the pipette and pulse centrifuge.

- 6.3 Place the plate on the thermocycler and launch the "BST" program:

35m

A	B	C
Step	Temperature	Duration

	A	B	C
	1	+37°C	15 min
	2	+10°C	20 min

BST program on thermocycler

6.4 Purify using  Agencourt AMPure XP **Beckman Coulter Catalog #A63880** XP 1X (

 30 µL); see section 1.

Elute in  15 µL ultra pure water.

II

Expected result

15µL of post BST purified solution, 12µL is needed for PCR step
optional you can use the remaining 3µL for quantify or check DNA/library size and quality pre PCR

optional you can conserve for 24h sample after purification step at -20°C and done PCR step the next day

PCR1 amplification + Purification AMPure

1h

7 There is the DNA construction we want to obtain at the end of this PCR step :



7.1 In a LightCycler plate® 480 Multiwell plate 96 white arranged on refrigerated stand: deposit  30 μL of each fluorescent standards ref KK2701 (1 to 4) in duplicates in 8 separate wells. Keep the plate away from light.

7.2 Prepare a volume of Mix according to the number of samples:

A	B
For 1 sample	volume
DNA	12 μL
adapter PE 34 solution 10 μM	1.5 μL
INDEX 10 μM	1.5 μL
KAPA HiFi Master Mix 2X	15 μL
Final Volume =	30 μL

Mix for 1 sample

Roche KAPA Hifi PCR KIT 250rx

Index and adapter sequence on attached file

Index used: AD001; AD002; AD003; AD004; AD005; AD006; AD007; AD008; AD009; AD0010; AD0011; AD0012.

7.3 Deposit  18 μL of Mix per well and pulse centrifuge.

7.4 Deposit  12 µL of DNA per well, mix with the pipette and pulse centrifuge.

7.5 Place the plate in the lightcycler and launch the following the "RTLib1" program.




	A	B	C	D	E	F
	Step	temperatur e	acquisitio n	Duratio n	spleen	Cycles
	Pre-incubation	98°C	None	45sec	4.4°C/sec	1 cycle
	Amplification	98°C	None	15sec	4.4°C/sec	20 cycles
		62°C	None	30sec	2.2°C/sec	
		72°C	Single	30sec	4.4°C/sec	

RTLib1 program on lightcycler

PCR should be stopped when signal is between standards 1 and 3 (optimal library amplification range), generally between cycles 10-14. more for low quality or low quantity DNA

Important : Do not over amplify, it will reduce you sequencing quality, stop before reachin plateau phase of amplification

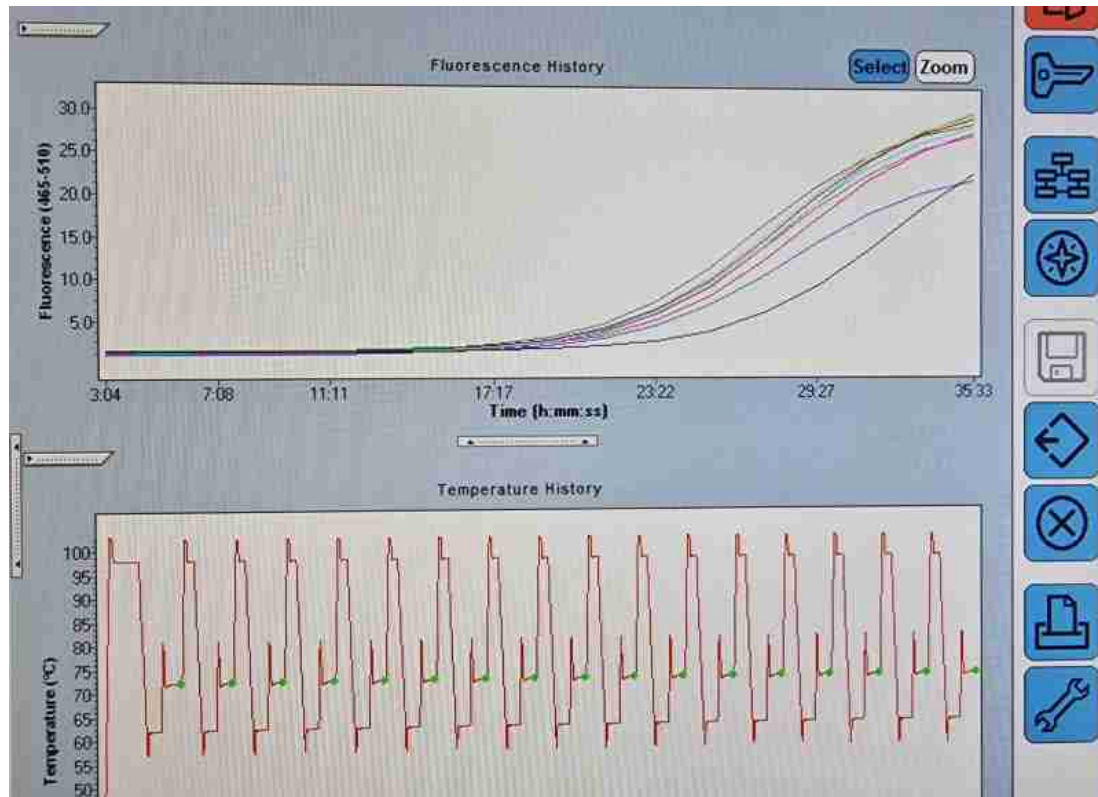

Asset URL: stop before plateau phase

- 7.6 Stop the run just after the last acquisition (green dot on the graph) when the double-stranded DNA is full length. If necessary, transfer the "passed" samples to another AB 1400 plate by piercing the film with a scalpel blade at each well concerned.
- 7.7 Repeat if necessary for samples requiring more subsequent cycles, using the same program but without pre-incubation.
- 7.8 Purify using  Agencourt AMPure XP **Beckman Coulter Catalog #A63880** XP 1X ( 30 μL); see section 2.
Eluate the samples in  30 μL of ultra pure water.
- 7.9 Perform an analysis on Qiaxcel ( 3 μL) and dose on Nanoquant ( 2 μL).

7.10

II

Expected result

after 10-16 cycles of amplification + purification
Mass concentration from 2 to 110ng/ μ L adn target 30ng/ μ L
DNA size from 250 to 550pb target 400pb peak

optional you can conserve for 24h sample after purification step at -20°C and continue the next day

Preparation of Bulk

1h 30m


- 8 Prepare an equimolar mix of all 48 sample by adding $\text{1 } \mu\text{L}$ to $\text{25 } \mu\text{L}$ of each purified sample from PCR1 step. Volume will depend of concentration and size of DNA.

- 8.1 Match 0.5 pmol per samples ex : 133 ng for 400pb DNA.

1h

Each volume is added in a 1.5 mL DNA low bind tube, total volume between

$\text{100 } \mu\text{L}$ to $\text{1000 } \mu\text{L}$

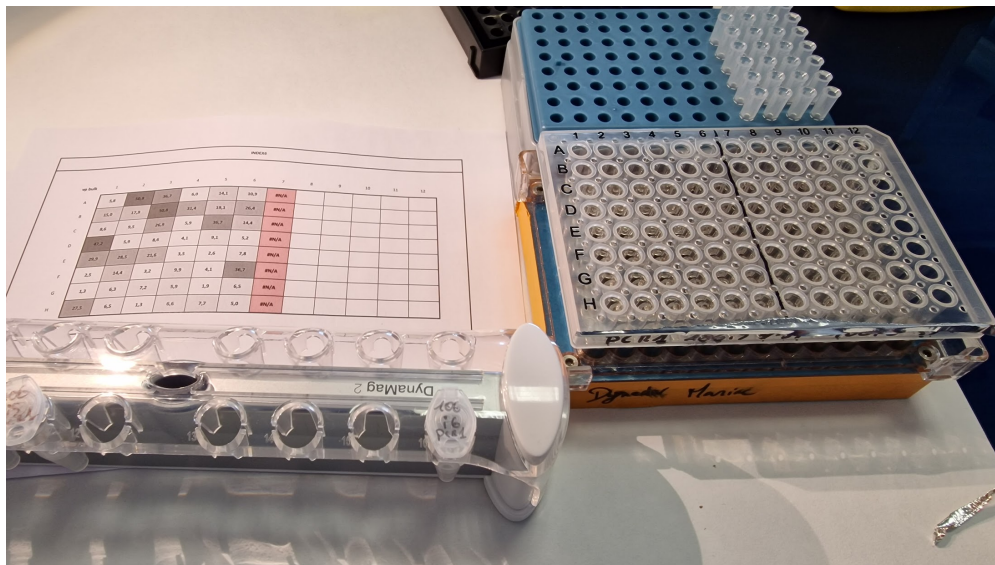



Bulk tubes before elution step.

8.3

Expected result

30 μ L of 50-200ng/ μ L DNA bulk, size 400pb peak
1-7 μ L is needed for Hybridization step

you can conserve DNA bulks for weeks at -20°C

you can repeat library preparations protocols for multiple bulk and when you have enough N bulks of 48/specimens so you can proceed multi hybridization step in same time and make N*48 specimens per sequencing.

Part 2: Hybridization step for targeted DNA sequencing

1d

- 9 This part of the protocol describes dual hybridization process for targeted DNA sequencing using Daicel Arbor Biosciences Mybait KIT V5 Angiosperm353 and custom Annonacea kit.
- In this protocol we use a bulk of 48 samples for each reaction (Arbor recommend max 8 sample per reactions). This protocol describes the quantity for one reaction or 8 reactions (48 or 384 samples)

II

A total of 5.5 μL of baits is required per reaction. We used a ratio of a **3/4** Angiosperms353 baits and **1/4** Annonaceae baits per reaction following Henriks et al. (2021):

4.125 μL of **Angiosperms353 baits**

1.375 μL of **Annonaceae baits**

Bait mix preparation

15m

- 10
- Program the Thermo Mixer so that it reaches 60 °C and a second Thermo mixer to 95 °C . Place the Arbor Biosciences Baits on ice.

Heat the Hyb N and Hyb S bait reactives at 60 °C and vortex to dissolve the aggregates.

Prepare a volume of Mix "**HYB**" according to the number of reactions planned, vortex and centrifuge:

	A	B	C
	Component	μL for 1 reaction	μL for 8 reactions
	Hyb N (+4°C)	9.25	74
	Hyb S (+4°C)	0.5	4
	Hyb D (-20°C)	3.5	28
	Hyb R (-20°C)	1.25	10
	Baits (-80°C)	4.125 Angio353 + 1.375 Annonaceae	33 Angio353 + 11 Annonaceae
	TOTAL	20 μL	160 μL



Mix "HYB" preparation

11 (-20°C)

15m

Incubate the **HYB** mix  00:10:00 at  60 °C and vortex to recover the condensate

For each capture reaction, aliquot  18.5 µL in  1.5 µL of DNA in a low bind tube.

We refer to these reaction aliquots of hybridization Mix as "HYBs"

Let stand  00:05:00 at  Room temperature

Prepare a volume of Mix "**Block**" according to the number of reactions of capture, vortex and centrifuge:

	A	B	C
	BLOCK	per capture	for 8 captures
	Block O (-20°C)	5 µl	40 µl
	Block X (-20°C)	0.5 µl	4 µl
	Total volume	5.5 µl	44 µl

Mix "**Block**" preparation

For each capture reaction, aliquot  5.0 µL of "Block" in a 0.5mL tube.

Hybridization

1d

12 Add  7 µL of BULK DNA library 50-100 ng/µL to each  5 µL aliquot of "Block" and homogenize by pipetting. We then obtain a Mix of  12 µL

10m



These libraries with Blocker Mix aliquots are now called "LIB" (block + libraries)

Place the "LIB" on the second Thermo Mixer at  95 °C for  00:05:00 .

Transfer the "LIB" and the "HYBs" mix into the Thermo Mixer at 65 °C

00:05:00

800 rpm

12.1 Transfer 18.5 μL of each "HYBs" to each corresponding "LIB" and homogenize by pipetting 5 times.

5m

Incubate at hybridization temperature of 65 °C Overnight (16h-24h).

800 rpm



eppendorf Thermomixer with ThermoTop , 65°C 800RPM

Capture and washing

1h

13 Bring Hyb S (the reactive, not the mix) and wash buffer to Room temperature in order to prepare the washing solutions about **1h30** before stopping the hybridization reaction.

Prepare the Wash Buffer X in a falcon tube of 15 mL

A	B
HYB S (+4°C)	40 µL
Water	3.96 mL
Wash Buffer	1 mL
Final volume	5 mL

Wash Buffer X preparation for 3 reactions

Wash Buffer X can be stored at 4 °C for up to 1 month

- 13.1 Vortex and place the Wash Buffer X at hybridization temperature (+ 65 °C in a water bath) at least 00:45:00 before use.

45m

Vortex to homogenise capture beads (kit , box at 4°C).

For each capture reaction, aliquot 30 µL of beads (Box at 4 °C)

TIP: Beads can be prepared in 8 (or fewer) reaction batches (240 µL) in a 1.5 µL tube. Multiply all volumes by the number of reactions in the batch; i.e., for 8 reactions-worth, wash with 1.6 mL and resuspend in 560 µL Binding Buffer, then aliquot 70 µL suspension to individual tubes.

- 13.2 Arrange tubes on a magnetic stand for 00:01:00 to 00:02:00 until the aggregate is formed and remove the supernatant.

3m 3s

Repeat 3 consecutive times:

Add 200 µL of Binding Buffer in each tube, vortex 00:00:03 , centrifuge briefly and finally eliminate the supernatant on magnetic stand.

Resuspend beads in 70 µL of Binding Buffer.

- 13.3 Place the beads at hybridization temperature on Thermo Mixer at least 00:02:00 before transfer.

7m

Quickly transfer each capture reaction (30 µL) into a tube containing the beads and homogenize by pipetting.

Incubate the library and beads mixture at hybridization temperature of **65 °C** on Thermo Mixer for **00:05:00** at 1200 rpm.

Eliminate the supernatant on magnetic stand.

13.4 Repeat 3 consecutive times:

Add **375 µL** of Wash Buffer X (heated) in each tube, vortex, centrifuge briefly, incubate at hybridization temperature on Thermo Mixer for **00:05:00** minutes at 1200 rpm and finally eliminate the supernatant on magnetic stand.

5m

Elution

15m

- 14 Add **25 µL** of Buffer E box -20° (bring back to **Room temperature** before hand) Incubate at **98 °C** for **00:05:00** on ThermoMixer or Heating Block (use lid). Place the tubes on a magnetic stand and quickly recover the supernatant containing the enriched libraries and transfer it into a new PCR tube / plate.

5m

RT-PCR 2 (final)

30m

- 15 In a LightCycler plate® 480 Multiwell plate with 96 arranged on refrigerated plate.



Deposit **50 µL** of each fluorescent marker 1 to 4 repeated twice in 8 separate wells on the plate. Keep the plate away from the light.
Prepare a Mix for each INDEX used. There are 12 Illumina INDEX that can be used.

A	B
By well:	
Bulk DNA	20µL
adapter PE 34 solution 10µM	2.5µL
INDEX 10µM	2.5µL
KAPA HiFi Master Mix 2X	25µL
Final volume	50µL

Mix for each INDEX per well

index and adapter sequence on attached file

Index used: AD001; AD002; AD003; AD004 AD005 AD006 AD007 AD008 AD009
AD0010 AD0011 AD0012

Deposit  30 µL of Mix (index+PE34+kapa pcr mix) per well and pulse centrifuge.

Deposit  20 µL of DNA per well, mix with the pipette and pulse centrifuge.

Place the plate at the lightcycler and launch the following program:

	A	B	C	D	E	F
	Step	Temperature	acquisition	Duration	speed	Cycles
	Pre-incubation	98°C	None	45sec	4.4°C/sec	1 cycle
	Amplification	98°C	None	15sec	4.4°C/sec	55 cycles
		62°C	None	30sec	2.2°C/sec	
		72°C	Single	30sec	4.4°C/sec	

Lightcycler program

Stop the run before the plateau phase. Generally 10 to 16 cycles. If necessary, add some extra cycles for some samples. Important : Do not over amplify, it will reduce you sequencing quality.

Carry out a purification step using

 Agencourt AMPure XP **Beckman Coulter Catalog #A63880** XP 1X on half the volume **25µL** , keep the leftover volume (approx.  25 µL) at  -20 °C as backup.

15.1 Proceed with the purification with  1.5 mL DNA low bind tube as follows:

5m

1. Prepare the sample: Mix the 25µL DNA sample with 25µL AMPure.
2. Bind: Incubate the mixture for 10 minutes at room temperature to allow the DNA to bind to the beads and put on magnetic stand;
3. Wash: after 5 minutes on magnetic stand remove supernatant;
4. Wash the beads twice with 70% ethanol to remove impurities;
5. Let dry for 5 minutes to remove residual ethanol and dry the purified DNA;
6. Remove from magnetic stand and resuspend the dried DNA by adding  20 µL water.



Let DNA elute at least  00:05:00 min.

Perform a quantification with Nanodrop and quality analysis on Qiaxcel.

Preparation of bulk for illumina sequencing

30m

- 16 Prepare an equimolar mix of all samples.
We end up with a single tube with all indexes mixed, and ready to sequence.

30m

Refer to your sequencing platform for final volume and concentration.

In our case we needed a minimum of 50 µL and DNA fragments should be at 400pb.

All samples were sequenced with Novogen using a NovaSeq PE150.

Protocol references

Couvreur TLP; Helmstetter AJ, et al. (2019) Phylogenomics of the Major Tropical Plant Family Annonaceae Using Targeted Enrichment of Nuclear Genes. *Frontiers in Plant Science* 9, <https://doi.org/10.3389/fpls.2018.01941>

Hendriks KP et al. The best of both worlds: Combining lineage-specific and universal bait sets in target-enrichment hybridization reactions (2021). *Applications in Plant Sciences* 9. <https://doi.org/10.1002/aps3.11438>.

Annonaceae DNA extraction protocol from silicagel dried and herbarium preserved leaves

DOI

dx.doi.org/10.17504/protocols.io.5qpvorqx9v4o/v1

- Théo Delaigue¹,
- Vincent Soulé¹,
- [Thomas LP Couvreur](#)¹,
- [Cedric Mariac](#)¹
- ¹DIADE, Univ Montpellier, CIRAD, IRD, Montpellier, France

Acknowledgements

The authors thank Kevin Bethune who worked on a previous version of this protocol (published in Couvreur et al. 2019)