

May 05, 2018 Version 1

Light/dark preference test for adult zebrafish (*Danio rerio*) V.1

DOI

dx.doi.org/10.17504/protocols.io.puydnxw

Caio Maximino¹


¹Universidade Federal do Sul e Sudeste do Pará

Fish behavior and physi...



Caio Maximino

Universidade Federal do Sul e Sudeste do Pará

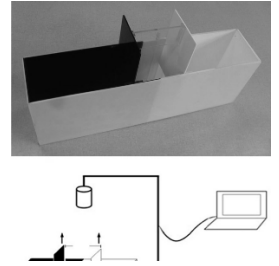

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.puydnxw

Protocol Citation: Caio Maximino 2018. Light/dark preference test for adult zebrafish (*Danio rerio*). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.puydnxw>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: April 30, 2018

Last Modified: May 05, 2018

Protocol Integer ID: 11896

Keywords: Zebrafish, Anxiety, Scototaxis, Light/dark preference test, LDT


Abstract

The light/dark test is increasingly being adopted as one of the main behavioral assays in zebrafish neuroscience research. It is based on the innate preference of adult zebrafish for a black vs. a white compartment. Modifications and extensions of the protocol increased its breadth, but also its heterogeneity. The protocol presented here involves confining the animal in a central compartment for a initial 3-min acclimation period, followed by free exploration of an acrylic tank with a black and a white compartment. In addition to observing the spatial distribution of the animal, "ethological" variables (risk assessment, erratic swimming, freezing, thigmotaxis) are also manually scored in the white compartment; these endpoints are selectively sensitive to drug treatments. The protocol is intended to measure anxiety-like behavior, with anxiogenic effects decreasing time on white and increasing risk assessment and thigmotaxis; some treatments also increase erratic swimming and freezing.

Guidelines

Standard guidelines for the ethical care of fish in research apply. Before beginning experiments, researchers should obtain all legally binding licenses, including IACUC approval.

Safety warnings



 The behavioral protocol is usually very innocuous for researchers. Care must be taken while preparing drugs for injection, following standard biosecurity practices. Moreover, care must be taken while capturing fish to not expose open wounds to tank water, reducing the risk of acquiring waterborne infections (cf. Murray et al., 2016, doi: 10.1089/zeb.2015.1206, for a primer on zebrafish biosecurity).

Before start

1. Room illumination is a critical variable, as it can alter fish behavior and the aversiveness of the white compartment. Before each trial, light levels should be measured; variations above or below 2 SDs from the average should be considered as too discrepant for an experiment. Ideally, lighting above the tank should be on the range of 500-600 lux.
2. Use Gaussian white noise throughout the experiment to cover unwarranted noise in the experimental room.
3. Bring the animal in a temporary transport recipient into the behavioral test room.

Apparatus

- 1
 - An acrylic tank (15 cm × 10 cm × 45 cm height × width × length) is used that is divided equally into one-half black and one-half white.
 - Walls and bottom are either black or white, so as to warrant uniform substrata for each compartment.
 - The water column is kept at 10 cm, yielding a final volume of 4.5 liters.
 - Water should be de-chlorinated.
 - The colored material chosen should not be reflective, to avoid the tendency of those animals that present shoaling tendencies to behave in relation to their own reflection.
 - The tank contains central sliding doors, colored with the same color of the aquarium side, thereby defining an uncolored central compartment measuring 15 cm × 10 cm × 10 cm.

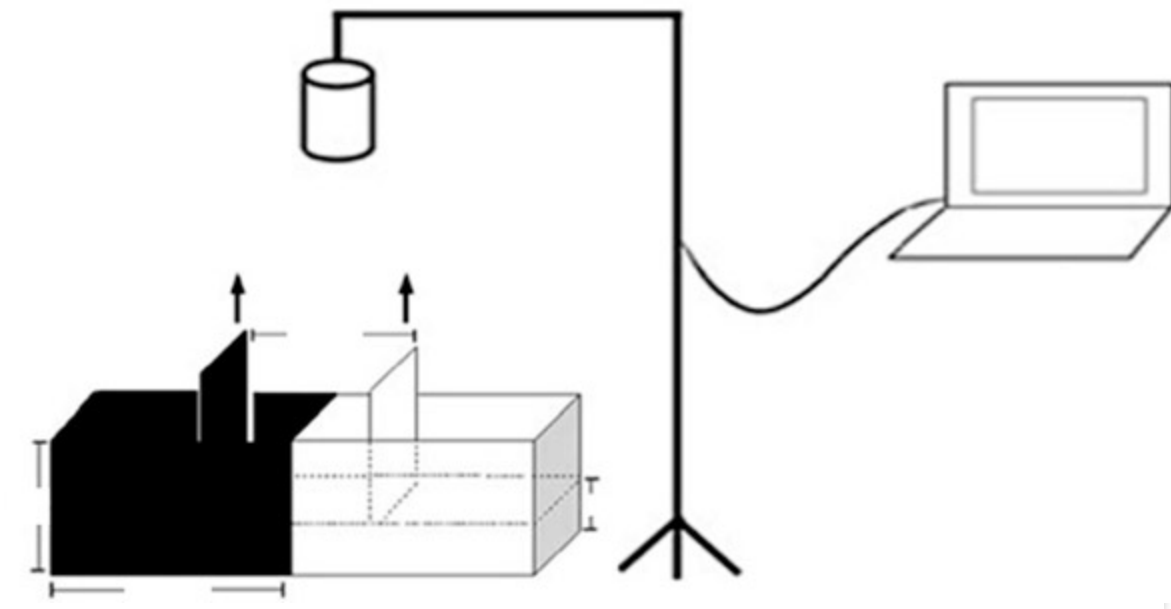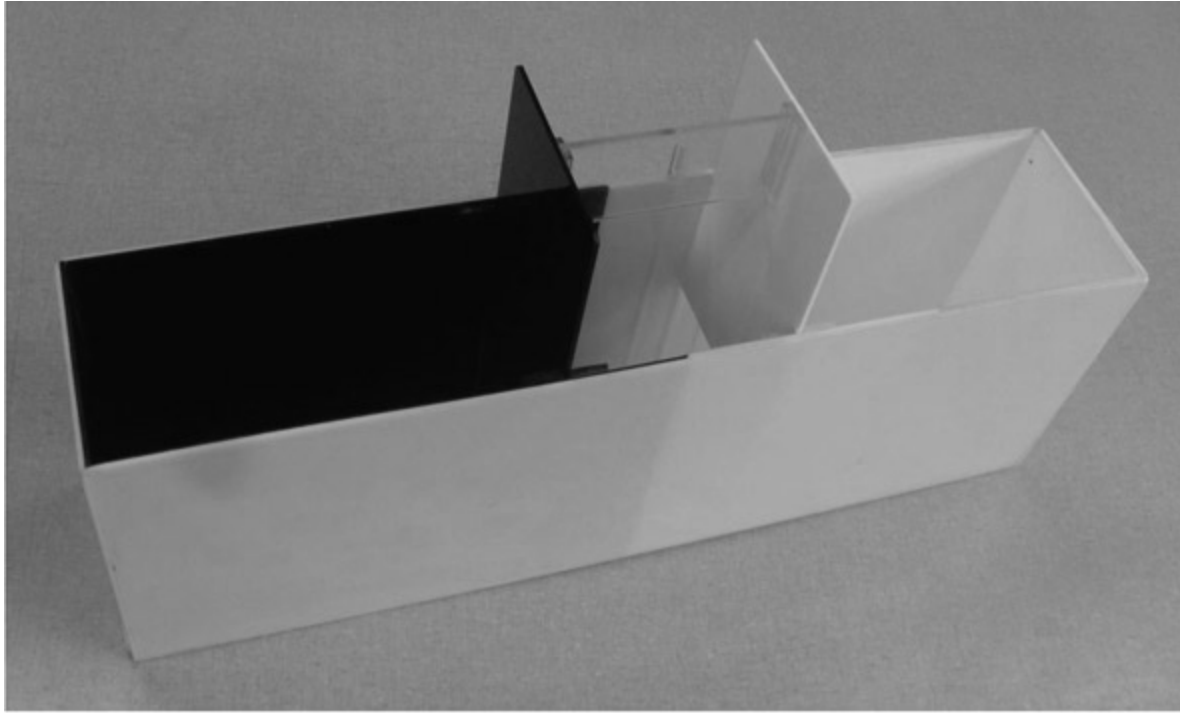



Preparation steps

- 2 Room illumination is a critical variable, as it can alter fish behavior and the aversiveness of the white compartment. 2. Before each trial, light levels should be measured; variations above or below 2 SDs from the average should be considered as too discrepant for an experiment. Ideally, lighting above the tank should be on the range of 500-600 lux.



- 3 Use Gaussian white noise throughout the experiment to cover unwarranted noise in the experimental room.
- 4 Bring the animal in a temporary transport recipient into the behavioral test room.

Behavioral assay

- 5 Transfer the animal to the central compartment of the tank with a net. Allow the animal to acclimate to the tank for 3 min
 00:03:00 Duration of acclimation
- 6 Remove the sliding doors, allowing the fish to swim freely in the apparatus for 15 min. Video recording is recommended.
 00:15:00 Duration of observation
- 7 Transfer the animal to a discard tank with a net.

Endpoints

- 8 Videos should be analyzed blindly (i.e., the experimenter should not know the treatment/genotype of the animal) by two experienced researchers. An event/spatial distribution logger, such as X-Plo-Rat (<https://github.com/lanec-unifesspa/x-plo-rat>).
- 9 Variables to be recorded (and ZBC references) are:-Time spent in the white compartment (ZBC 1.137)-Latency to enter the white compartment-Number of entries in the white compartment-Duration of entries in the white compartment-"Risk assessment" events (a fast (<1 s) entry in the white compartment followed by re-entry in the black compartment, or a partial entry in the white compartment (i.e., the pectoral fin does not cross the midline))-Duration and number of erratic swimming (ZBC 1.51): a zig-zag, fast, unpredictable course of swimming of short duration.-Duration and number of freezing events (ZBC 1.68): complete cessation of movements except eye and opercular movements.-Time on thigmotaxis (ZBC 1.173): swimming in a distance of 2 cm or less from the nearest wall.
- 10 If X-Plo-Rat is used, risk assessment, erratic swimming, freezing, and thigmotaxis are assigned to keyboard keys (A, S, D, and F are a good choice). The I-maze pre-configuration can be used, with the "walled" squares equivalent to the black compartment.



Software

X-Plo-Rat

NAME

Windows

OS

Fernando Cardenas

DEVELOPER

<https://github.com/lanec-unifesspa/x-plo-rat>

SOURCE LINK

Expected result

Under these conditions, adult zebrafish present preference for the dark compartment, spending about 80% of a 15 min session on that portion of the tank. Drugs or treatments with an "anxiogenic"-like effect increase the preference, and can also increase risk assessment, thigmotaxis, erratic swimming, and/or freezing in the white compartment. Conversely, drugs with an "anxiolytic"-like effect reduce some or all of these endpoints.