

Dec 21, 2023

Protocol for Neuronal Live-imaging of primary cultures

 [Science Advances](#)

 In 2 collections

DOI

<https://dx.doi.org/10.17504/protocols.io.q26g7pn4qgwz/v1>

Camila Pulido^{1,2}, Timothy A. Ryan^{1,2}

¹Department of Biochemistry, Weill Cornell Medicine, New York, NY 10065, USA;

²Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, Maryland 20815, USA



Camila Pulido

Weill Cornell Medicine

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.q26g7pn4qgwz/v1>

External link: <https://doi.org/10.1126/sciadv.adn6016>

Protocol Citation: Camila Pulido, Timothy A. Ryan 2023. Protocol for Neuronal Live-imaging of primary cultures. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.q26g7pn4qgwz/v1>

Manuscript citation:

Kokotos AC, Antoniazzi AM, Unda SR, Ko MS, Park D, Eliezer D, Kaplitt MG, Camilli PD, Ryan TA Phosphoglycerate kinase is a central leverage point in Parkinson's disease-driven neuronal metabolic deficits. Science Advances 10(34). doi: [10.1126/sciadv.adn6016](https://doi.org/10.1126/sciadv.adn6016)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: December 13, 2023

Last Modified: May 31, 2024

Protocol Integer ID: 92279

Keywords: ASAPCRN, Imaging, Fluorescence, Microscopy, Lasers, Neurons, Primary Culture, Biosensors, imaging of primary cultures neuron, imaging techniques on primary neuron, imaging neuronal activity, primary cultures neuron, protocol for neuronal life, neuronal life, primary neuron, neuron, neuronal activity, neuronal function, engineered biosensor, comprehension of neuronal function, cell, imaging, imaging technique, intricate cell

Funders Acknowledgements:

ASAP

Grant ID: 000580

Abstract

Neurons stand out as remarkably intricate cells. The development of live-imaging techniques on primary neurons expressing genetically engineered biosensor has proven pivotal in advancing our comprehension of neuronal function at both cellular and molecular level. This protocol introduces a concise method for live-imaging neuronal activity.

Attachments



[Protocol for neurona...](#)

284KB

Guidelines

The animals should be used in accordance with protocols approved by national and institutional regulatory organizations.



Materials

Media Supplies

- 6-cyano-7-nitroquinoxaline-2,3-dione (CNQX; Alomone Labs, Cat. No. C-141)
- Calcium Chloride (Millipore-Sigma, Cat. No. C3881)
- DL-2-Amino-5-phosphonovaleric acid (APV; Millipore-Sigma A5282)
- Glucose (Millipore-Sigma, Cat. No. G7021)
- HEPES (Millipore-Sigma, Cat. No. H3375)
- Magnesium Chloride (Millipore-Sigma, Cat. No. M0250)
- Potassium Chloride (Millipore-Sigma, Cat. No. P9333)
- Sodium Chloride (Millipore-Sigma, Cat. No. S5761)

Other Supplies

- Ethanol (200 proof; Koptec, Cat. No. 6175)
- Grease (Millipore-Sigma, Cat. No. Z273554-1EA)
- Wipes (Kimtech Science)
- P35 dishes (Falcon, Cat. No. 351008)

Stock solutions

APV (25mM Stock)

1. Add 20L ddH₂O to whole bottle of APV
2. Sonicate at ~55C for 30 minutes.
3. Make 1mL aliquots and store at -20C.

CNQX (10mM Stock)

1. Add 10 mL of ultrapure H₂O to the whole bottle of CNQX.
2. Vortex until dissolved.
3. Make 1 mL aliquots and store at -20 °C.

Media solution

Imaging Solution (119mM NaCl, 2.5mM KCl, 5mM D-Glucose, 50mM HEPES, 2mM CaCl₂, 2mM MgCl₂, 50uM APV, 10uM CNQX)

1. Add 450mL of ddH₂O
2. Add 3.47g of NaCl
3. Add 0.45g of Glucose
4. Add 5.96g of HEPES
5. Add 1.25mL of KCl (1mM stock prepared in ddH₂O)
6. Add 500uL of CaCl₂ (2mM stock prepare in ddH₂O)
7. Add 500uL of MgCl₂ (2mM stock prepare in ddH₂O)

8. Add 1mL of APV (25mM stock)
9. Add 500uL of CNQX (10mM stock)
10. Adjust pH to ~7.4 with NaOH
11. Add water to 500mL
12. Confirm pH is correct
13. Filter to a 500mL bottle

Equipment

- 40x, 1.3 N.A Objective (like Zeiss-Fluar)
- 488nm Laser (like Coherent-OBIS 488nm LX Laser)
- CO₂ Incubator (like VWR Symphony)
- EMCCD Camera (like Andor iXon Ultra 897)
- Epifluorescence Microscope (like Zeiss ObserverA1)
- Filter cube set for GFP excitation and emission (like Chroma filters)
- Laser speckle reducer (like Optotune LSR 3005)
- Laser/stage/stimulation-Master board controller (like a custom-designed **Arduino board**)
- Light transmission microscope (like Olympus CK40)
- Mercury arc lamp (like Zeiss HBO 100)
- Perfusion system (like Automate Scientific Valvelink 8.2)
- Microscope XYZ Stage (like Applied Scientific Instrumentation stage)
- Stimulus isolator (like World precision Instruments A385)
- Temperature controller for Objective (like Minco)

Safety warnings

- ! ▪ Laser light, because of its special properties, can cause safety hazards not related to light from common sources. The safe use of lasers requires that all laser users and all persons near the laser system understand the possible danger. The safe use of the laser depends on the user understanding the instrument and the properties of coherent, strong beams of light.
- Direct eye contact with the output beam from the laser will cause damage and possible injury to the eyes.

Before start

- Turn on live-imaging setup: Lasers, microscope, mercury arc lamp, stimulus isolator, XYZ stage, imaging computer.
- Open acquisition and controller programs (i.e. camera program, laser's controller, ImageJ).
- Make sure the EMCCD camera has reach the working temperature (~ -80°C) before start imaging.

Mounting the coverslip to recording chamber

3m

1

Add  2 mL of 'Imaging Solution' into a P35 dish.

2

Cut a piece of a new coverslip and fix it with some grease into the 'Stimulator insert' (Figure 1).

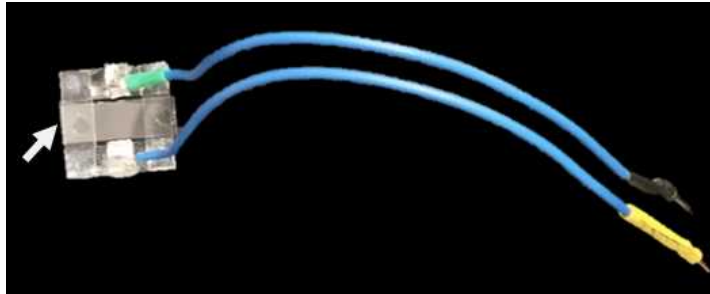


Figure 1. Stimulator Insert with a coverslip piece (white arrow)

3

Add a drop of grease to every corner (bottom face) of the 'Stimulator insert' (Figure 2).

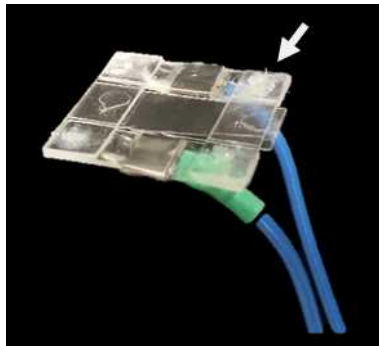


Figure 2. Stimulator Insert: Bottom face corners with a drop of grease (white arrow).

4

Add grease around the circle (facing up) of the 'platform chamber' (Figure 3).

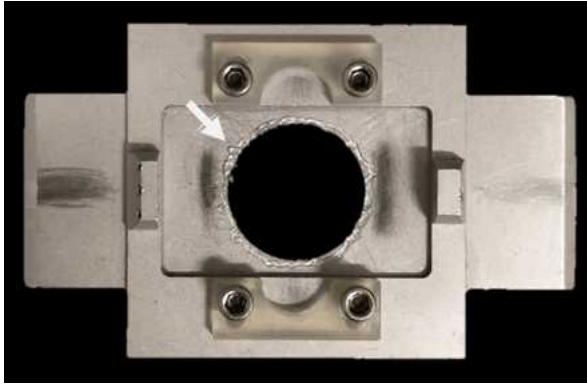


Figure 3. Platform chamber with grease (white arrow).

- 5 Take a DIV14 – DIV18 neuronal cultured dish from the incubator.
- 6 Assess for cell density, health and overall conditions using a basic light transmission microscope to ensure the suitability of the dish before starting the experiment.
- 7 Place the coverslip containing plated cells into the petri dish prepared at [go to step #1](#), and leave it for [00:00:10](#) to [00:00:20](#) to equilibrate with the 'Imaging Solution'.
- 8 Attach the coverslip to the 'Stimulator insert'; with the help of some tweezers make sure every corner of the 'Stimulator insert' is well attached to the coverslip by the grease.
- 9 Seal the 'Stimulator insert' to the 'platform chamber', clean any excess of grease that might be around and make sure there is no leak.

30s

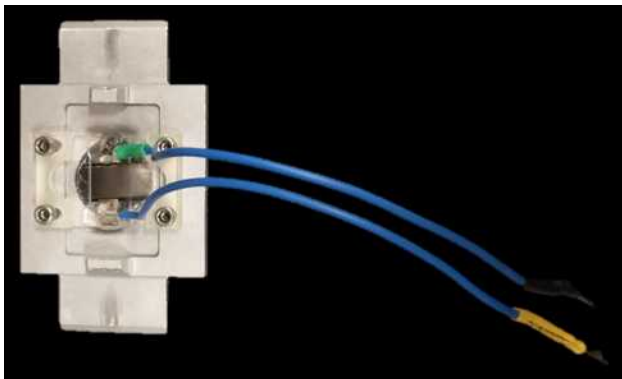


Figure 4. Recording Chamber assembled.



- 10 Clean the bottom of the coverslip with a wipe with  Sample , making it ready for imaging.

Checking for transfected neurons under the microscope

5m

- 11 Place the 'recording chamber' into the microscope's stage.
- 12 If electrically stimulating neurons, connect the 'recording chamber' to the stimulus isolator.
- 13 Start perfusing 'Imaging Solution' through the coverslip at a speed of ~  100 $\mu\text{L min}^{-1}$. Make sure all perfusion lines with additional test solutions are calibrated to run at the same speed.
- 14 Using transmission light find the focal plane where neurons are located.
- 15 Change microscope illumination to the mercury arc lamp, making sure the desired dichroic filter is in place. (i.e. GFP-filter cube).
- 16 Verify the presence of transfected neurons based on their baseline fluorescence emission. If the biosensor's baseline fluorescence is insufficiently bright due to its inherent nature (i.e. vGLUT-pH), add a suitable solution (i.e. NH_4Cl -Tyrode's) into the chamber that may enhance the biosensor signal, facilitating the identification of transfected neurons.
- 17 Save all localizations where transfected neurons are found.
- 18 Proceed with the imaging protocol.

Neuronal Imaging Protocol

- 19 Flip the microscope mirror from the mercury arc lamp towards the laser path.



- 20 Make sure the following programs are open at the 'Imaging-setup computer': EMCCD camera control software (like Andor-Solis program), laser control software (like Coherent connection), master control program (like a custom-designed Arduino controller program) and image analysis program (like ImageJ with Time-Analyzer plugin).
- 21 Load the desired acquisition settings protocol for imaging in the EMCCD camera control software.
- 22 Select one of the previously saved stage localizations, re-focus and move the neuron/boutons to a desire field of view, ready to be imaged.
- 23 If electrically stimulating, make sure that the correct settings are selected (i.e. number of pulses, stimulation frequency, time to first pulse) and ready to Go when triggered by camera action.
- 24 When ready for imaging, proceed to record.

Additional notes

- 25 This live-imaging protocol can be adapted to be used with ratiometric biosensors (i.e. iATPSnFR2-HALO), by using the appropriate dual filter cube corresponding to the different fluorescent emissions and, by switching lasers between frames when image acquisition is performed.

Protocol references

Ryan, T.A. (1999). Inhibitors of myosin light chain kinase block synaptic vesicle pool mobilization during action potential firing. *J. Neurosci.* 19, 1317–1323.