

Feb 22, 2026

Sequencing bacterial and fungal metabarcode with PCR primer based barcoding and Nanopore v2026

Forked from [Sequencing bacterial and fungal metabarcode with PCR primer based barcoding and Nanopore](#)

DOI

<https://dx.doi.org/10.17504/protocols.io.kqdg3nw2pv25/v1>

Abigail Graetz¹, Ashley Jones², Benjamin Schwessinger²


¹ANU; ²Australian National University

BIOL3106/6106



Benjamin Schwessinger

Australian National University

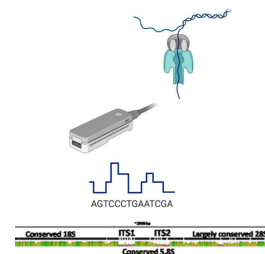

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.kqdg3nw2pv25/v1>



Protocol Citation: Abigail Graetz, Ashley Jones, Benjamin Schwessinger 2026. Sequencing bacterial and fungal metabarcode with PCR primer based barcoding and Nanopore v2026. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.kqdg3nw2pv25/v1>

Manuscript citation:

License: This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: February 22, 2026

Last Modified: February 22, 2026

Protocol Integer ID: 243820

Keywords: fungal metabarcode with pcr primer, different wheat pathogen, plant pathogen diagnostic, wheat pathogen, fungal metabarcode, major pathogen, uninfected wheat plant, pathogen, purification of the pcr amplicon, dna metabarcoding, step guide for pcr barcoding, pcr primer, biorxiv, repeated pcr, repeat pcr reaction, dna barcode, repeat pcr reactions on the sample, genomic, pcr barcoding, sequencing adapter, fungal taxa with high phylogenetic resolution, sequencing, pcr amplicon, nanopore this protocol, unique barcode sequence, fungal taxa, pcr reaction, single nanopore run, nanopore, sequencing adaptor, old wheat plant, infective agent for each treatment group, week old wheat plant, pcr, using primer, fidelity polymerase, sequencing run, dna sample, purification, molecular biology, high phylogenetic resolution, oxford nanopore dna sequencing, sequencing method, nanopore v2026, nanopore v2026 this protocol, minion sequencer

Abstract

This protocols is part of the ANU Biosecurity mini-research project #1 "Plant Pathogen Diagnostics: Visuals, subcultures, and genomics".

You will be provided four pots of 3-4 week old wheat plants that have been infected with different wheat pathogens. Each pot has been infected with one major pathogen. You will not know which pot has been infected with which pathogen. However, you will be provided a compendium of 10-15 wheat pathogens that will guide you to identify the infective agent for each treatment group. The fifth treatment group will be uninfected wheat plants which will be clearly identified. You can use treatment group #5 as negative control for your experiments.

In total, each group will obtain five pots each:

	A	B
	Treatment group 1	Unknown infective agent
	Treatment group 2	Unknown infective agent
	Treatment group 3	Unknown infective agent
	Treatment group 4	Unknown infective agent
	Treatment group 5	Uninfected control


This specific protocol describes the molecular biology and a step-by-step guide for PCR barcoding of the 16S and ITS regions using primers that contain barcodes, followed by ligation based amplicon sequencing with Nanopore. The initial steps of repeated PCR with high-fidelity polymerase and purification of the PCR amplicons has been performed by your demonstrator due to time limitations. We selected DNA samples of three research groups in week 5 so that we have three replicates per treatment group. This gives us 2x (16S and ITS) 15 PCR reactions plus 1 negative control. We skip the extraction control for this practical mini-research project. This is step 1 of the protocol.

This protocol is for reference only in 2026 as we opted to re-use previous years' data for this larger class. In previous year only 2-3 (out of 25) research groups could participate in the actual hands on sequencing, hence the actual learning effect was very limited across. We will bring along some MinION sequencers, flowcells, and an actual life sequencing run for you to have a closer look.

If you are interested to learn more about the molecular biology of DNA sequencing and Oxford Nanopore DNA sequencing/analysis, please reach out to Benjamin and/or Mareike as their might will be an opportunity for a research project as part of your degree, Honours or Master. We are always happy to chat about science.

We also teach into <https://programsandcourses.anu.edu.au/course/biol3210> where you can experience more DNA sequencing methods and analysis.

Conceptual overview:

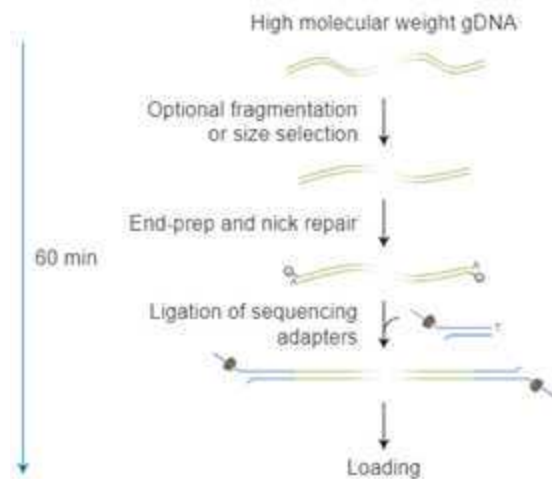

Overview of amplicon library prep we will not perform fragmentation and size selection.

- Repeat PCR reactions on the samples of three research groups using primers containing barcodes in them (see below).
- Clean and pool amplicons at equimolar ratios.
- End-prep the pooled amplicons to make them ready to ligate sequencing adapters.
- Ligate sequencing adaptars with DNA ligase.
- Make amplicon library ready for loading.
- Prepare for loading onto the flowcell.
- Prime flowcell to make it ready for loading.
- Load library.
- Start of sequencing run.
- Basecall during or after the sequencing run.

You can cite this protocol in the methods section of your report as for all other protocols. No need to write it all up again :).

Primers used:

We used ITS and 16S primers that directly include DNA barcodes as follows, where N₁₅₋₁₆ represents a unique barcode sequence



ITSFor (a.k.a. ITS1-F_KYO2_BARCODE): 5'-N₁₅₋₁₆-TAGAGGAASTAAAAGTCGTAA-3'

ITSRev (a.k.a. LR6_BARCODE): 5'-N₁₅₋₁₆-CGCCAGTTCTGCTTACC-3'

16SFor (a.k.a. 27F_BARCODE): 5'-N₁₅₋₁₆-AGAGTTTGATCCTGGCTCAG-3'

16SRev (a.k.a. 1492R_BARCODE): 5'-N₁₅₋₁₆-ACCTTGTTACGACTT-3'

References:

- Ohta, Nishi, Hirota and Matsuo, 2023, 'DNA metabarcoding workflow utilising nanopore long-read sequencing and consensus generation for rapid identification of fungal taxa with high phylogenetic resolution, BioRxiv, doi: <https://doi.org/10.1101/2023.04.14.536971>
- Youssef, Sheik, Krumholz, Najjar, Roe and Elshahed, 2009, 'Comparison of species richness estimates obtained using nearly complete fragments and simulated pyrosequencing-generated fragments in 16S rRNA gene-based environmental surveys', App. Env. Microbiol., vol. 75 (16)
- Hebert, P.D.N., Floyd, R., Jafarpour, S., Prosser, S.W.J., 2023. Barcode 100K Specimens: In a Single Nanopore Run, BioRxiv, doi: <https://doi.org/10.1101/2023.11.29.569282>

Attachments



[Ligation sequencing ...](#)

111KB



[ligation-sequencing-...](#)

3.2MB

Image Attribution

The logo image was generated with biorender.



Guidelines

This protocol describes how to carry out native barcoding of PCR amplicons that already contain barcodes using the Ligation Sequencing Kit (SQK-RBKLSK114). The PCR barcoding allows for 10–100,000 barcoded samples to be sequenced in one reaction. We will barcode 32 PCR amplicons (16 for ITS and 16 for 16S).

This approach is recommended for users who:

- wish to multiplex 10 to 10,000 samples to reduce sequencing cost per sample
- want to avoid costly (time and money) native barcoding reactions and two step PCR barcoding reactions
- want to conserve amplicon length

As always you need to bring a lab notebook, a printed version of this protocol, and a pen to record your adventures in the lab.

You must have read, understood, and follow the health and safety instructions provided in the "Overview Mini-Research Project #1 BIOL3106/6106" provided on Canvas (ANU learning portal).

You must have signed and returned one copy of the "Student Safety Declaration Form For Practical Class Work" before starting any laboratory work.

You must have read and understood the Hazard Sheets (Risk assessment) of all chemicals listed bellow in the "Safety Warnings" section. These Hazard Sheets are provided on Canvas as part of the "Overview Mini-Research Project #1 BIOL3106/6106" document.

Materials

Materials

- Pure DNA

Consumables

- R10.4.1 MinION flowcell **Flow Cell (R10.4.1) (nanoporetech.com)**
- SQK-LSK114 **Ligation Sequencing Kit V14 (nanoporetech.com)**
- Primers as described.
- Q5 High-Fidelity Polymerase
- NEBNext Ultra II End repair/dA-tailing Module (NEB, E7546)
- NEBNext Quick Ligation Module (NEB, E6056)
- Bovine Serum Albumin (BSA) (50 mg/ml) (e.g. Invitrogen™ UltraPure™ BSA 50 mg/ml, AM2616)
- magnetic beads (AMPure or equivalent) <https://www.beckman.com/reagents/genomic/cleanup-and-size-selection/pcr>
- 1.5 ml Eppendorf DNA LoBind tubes
- 0.2 ml thin-walled PCR tubes
- Nuclease-free water (e.g. ThermoFisher)
- Freshly prepared 80% ethanol in nuclease-free water
- Qubit Assay Tubes (ThermoFisher)
- Qubit dsDNA BR Assay Kit (ThermoFisher)

Equipment

- Hula mixer (gentle rotator mixer)
- Microfuge Magnetic rack
- Vortex mixer
- Thermal cycler/PCR machine
- Multichannel pipette
- P1000 pipette and tips
- P200 pipette and tips
- P100 pipette and tips
- P20 pipette and tips
- P10 pipette and tips
- P2 pipette and tips
- Ice bucket with ice Timer
- Qubit fluorometer (or equivalent for QC check)
- Computer with the latest version of MinKNOW installed.
- MinION Nanopore sequencer



Safety warnings

! This protocol does not require any hazardous substances.

You need to wear safety equipment at all times including lab coats, gloves, and safety goggles when handling chemicals and biological agents. While the major biological agents used for the wheat infection are pathogens commonly found in Australia, you must treat them as they were infective agents of general concern. Treat them with care. Do not remove them from the laboratory. Do not spread them via clothing. Use a dedicated notebook and pen to make notes during the mini-research project. Do not put anything into your mouth while in the laboratory. Wash your hands each time you leave the laboratory.

Preparation of input DNA

- 1 This first step has been performed by a demonstrator using DNA from the groups we selected in week 5. We selected three research groups so that we have three replicates per treatment group. This gives us 15 PCR reactions plus 1 negative control each for 16S and ITS PCRs. We skip the extraction control for this practical mini-research project.

We repeated the PCR like here **PCR reaction of marker regions (a.k.a metabarcodes) for bacteria and fungi** for these samples with the two following adjustments:

1. We used proof-reading high fidelity polymerase like Q5 for the PCR reaction.
2. We used ITS/16S primers that directly include DNA barcodes as follows, where N₁₅₋₁₆ represents a unique barcode sequence

ITSFor (a.k.a. ITS1-F_KYO2_BARCODE): 5'-N₁₅₋₁₆-TAGAGGAASTAAAAGTCGTAA-3'

ITSRev (a.k.a.LR6_BARCODE): 5'-N₁₅₋₁₆-CGCCAGTTCTGCTTACC-3'

16SFor (a.k.a. 27F_BARCODE): 5'-N₁₅₋₁₆-AGAGTTTGATCCTGGCTCAG-3'

16SRev (a.k.a. 1492R_BARCODE): 5'-N₁₅₋₁₆-ACCTTGTTACGACTT-3'

The PCR products were checked by agarose gel electrophoresis. These confirmed PCR products were then treated as follows.

- 1.1 Purify your PCR amplicons with magnetic beads (e.g. AMPure) by adding an equal volume of magnetic beads to the PCR reaction.

For example if your PCR reaction was 50 ul add 50 ul of beads.
- 1.2 Incubate while rotating the tube for 10 minutes at room temperature
- 1.3 Spin down the sample and pellet on a magnet. Keep the tube on the magnet, and pipette off the supernatant.
- 1.4 Keep the tube on the magnet and wash the beads with 200 ul of freshly prepared 80% ethanol without disturbing the pellet. Remove the ethanol using a pipette and discard.
- 1.5 Repeat the previous step 1.4.



- 1.6 Briefly spin down and place the tube back on the magnet. Pipette off any residual ethanol. Allow to dry for 1 - 2 mins, but do not dry the pellet to the point of cracking.
- 1.7 Remove the tube from the magnetic rack and resuspend the pellet in 20 ul of nuclease free water.
- 1.8 Incubate for 10 minutes.
- 1.9 Pellet the beads on a magnet until the eluate is clear and colourless, for at least 1 minute.
- 1.10 Remove and retain the full volume of eluate into a clean 200ul PCR tube.
- 1.11 Quantify 1ul of purified amplicon using Qubit BR as described in previous protocols. Note down the concentration.
- 1.12 Now we need to combine the 32 purified PCR amplicons [2 metabarcode PCR reactions {16S + ITS} * (3 reps * 5 treatment groups + 1 negative PCR control)] at equal molar ratios. We need about 400 ng DNA in total as input for the library prep as this equates to ~ 0.2 pmol when assuming an average PCR amplicon length between 1500 and 3000 bp.

This implies that we need to combine 12.5 ng per purified PCR amplicon and make it up to 50 ul total. These 50 ul contain our PCR amplicons in equal molar ratios and are the starting point for our library preparations.

To note is that we expect very low amounts for the negative PCR control and that it is best add as much as possible up to 12.5 ng to the DNA amplicon pool.

Native amplicon library preparation

- 2 These steps will be performed by the three groups we selected in week 5. Not much hands on but good fun to hang out.
- 2.1 Program the thermal cycler: 20°C for 5 minutes and 65°C for 5 minutes. Ramp it up and keep it at 20°C so we can start the program immediately.
- 2.2 In a 0.2 ml thin-walled PCR tube, mix the following:

Input Material	Volume [ul]
Pooled purified PCR amplicons (~400 ng)	50
Ultra II End-prep Reaction Buffer	7
Ultra II End-prep Enzyme Mix	3
Total	60

Between each addition, pipette mix 10-20 times.

- 2.3 Thoroughly mix the reaction by gently pipetting and briefly spinning down.
- 2.4 Incubate the tubes at 20°C for 5 minutes and 65°C for 5 minutes using the PCR machine we already set up. Once done, briefly put the tubes on ice to cool.
- 2.5 Resuspend the magnetic beads for DNA binding by vortexing.
- 2.6 Add an equal volume of 60 ul of magnetic beads to End-prep reaction.
- 2.7 Incubate while rotating the tube with your hand for 5 minutes at room temperature
- 2.8 Spin down the sample and pellet on a magnet. Keep the tube on the magnet, and pipette off the supernatant.

A demonstrator will help you with this step.
- 2.9 Keep the tube on the magnet and wash the beads with 200ul of freshly prepared 80% ethanol without disturbing the pellet. Remove the ethanol using a pipette and discard.
- 2.10 Repeat the previous step 2.9.



- 2.11 Briefly spin down and place the tube back on the magnet. Pipette off any residual ethanol. Allow to dry for 1 - 2 mins, but do not dry the pellet to the point of cracking
- 2.12 Remove the tube from the magnetic rack and resuspend the pellet in 61 μ l nuclease free water.
- 2.13 Incubate for 3 minutes while warming the tube with your fingers.
- 2.14 Pellet the beads on a magnet until the eluate is clear and colourless, for at least 1 minute.
- 2.15 Remove and retain the full volume of eluate into a clean 1.5ml Eppendorf PCR tube.
- 2.16 Now we will add the sequence adapters to the End-prepped PCR amplicons

Input Material	Volume [μ l]
End-prepped PCR amplicons	60
Ligation Buffer (LNB)	25
NEBNext Quick T4 DNA Ligase	10
Ligation Adapter (LA)	5
Total	100

Between each addition, pipette mix 10-20 times.

- 2.17 Thoroughly mix the reaction by gently pipetting and briefly spinning down.
- 2.18 Incubate the reaction for 10 minutes at room temperature.



- 2.19 Add 55 µl of resuspended magnetic beads to the reaction and mix by flicking the tube.
- 2.20 Incubate while rotating the tube with your hand for 5 minutes at room temperature
- 2.21 Spin down the sample and pellet on a magnet. Keep the tube on the magnet, and pipette off the supernatant.
- A demonstrator will help you with this step.
- 2.22 Wash the beads by adding 250 µl Short Fragment Buffer (SFB). Flick the beads to resuspend, spin down, then return the tube to the magnetic rack and allow the beads to pellet. Remove the supernatant using a pipette and discard.
- 2.23 Repeat the previous step 2.22.
- 2.24 Remove the tube from the magnetic rack and resuspend the pellet in 15 µl Elution Buffer (EB). Spin down and incubate for 5 minutes at room temperature.
- 2.25 Pellet the beads on a magnet until the eluate is clear and colourless, for at least 1 minute.
- 2.26 Remove and retain 15 µl of eluate containing the DNA library into a clean 1.5 ml Eppendorf DNA LoBind tube.
- 2.27 Quantify 1ul of purified amplicon using Qubit BR as described in previous protocols. Note down the concentration and calculate the volume you need to obtain 20 fmol assuming an average amplicon length of ~ 2500 bp (just an approximation between 1500 and 3000 bp on average).

Remember:

$n = c \cdot V$ with n being the amount in ng, c the concentration in ng/ul, and V the volume in ul.

and use the following online calculator.

[Biomath Calculators](#) | [DNA Calculator](#) | [Vector Insert Ratio \(promega.com.au\)](#)

- 2.28 Make up 20 fmol of final library in 12 ul of Elution Buffer.



- 2.29 Add 37.5 uL of Sequencing Buffer (SB) and 25.5uL of Library Beads (LIB) to the sample.

Note – the library beads settle on the bottom of the tube very quickly. Make sure to mix well before pipetting.

- 3 We now have our library ready to load and will wait for the rest of the class to come in.

Flowcell preparation for sequencing aka priming

- 4 These steps will be performed by one group member while the rest of the class is watching.

If you are not familiar with MinION flowcells check out this explanation video here [Community - Flow cell introduction \(nanoporetech.com\)](#).

This is a useful video for flowcell priming and loading [Community - Priming and loading your flow cell \(nanoporetech.com\)](#).


- 4.1 Perform a flowcell check as shown here [Community - MinKNOW: Flow cell check \(nanoporetech.com\)](#).

- 4.2 To prepare the flow cell priming mix with BSA, combine the following reagents in a fresh 1.5 ml Eppendorf DNA LoBind tube. Mix by inverting the tube and pipette mix at room temperature.


	Reagents	Volume [ul]
	Flow Cell Flush (FCF)	1170
	Bovine Serum Albumin (BSA) at 50 mg/ml	5
	Flow Cell Tether (FCT)	30

Reagents	Volume [ul]
Total volume	1205


- 4.3 Open the MinION device lid and slide the flow cell under the clip. Press down firmly on the flow cell to ensure correct thermal and electrical contact. This step might have been already performed if you did the flow cell check.

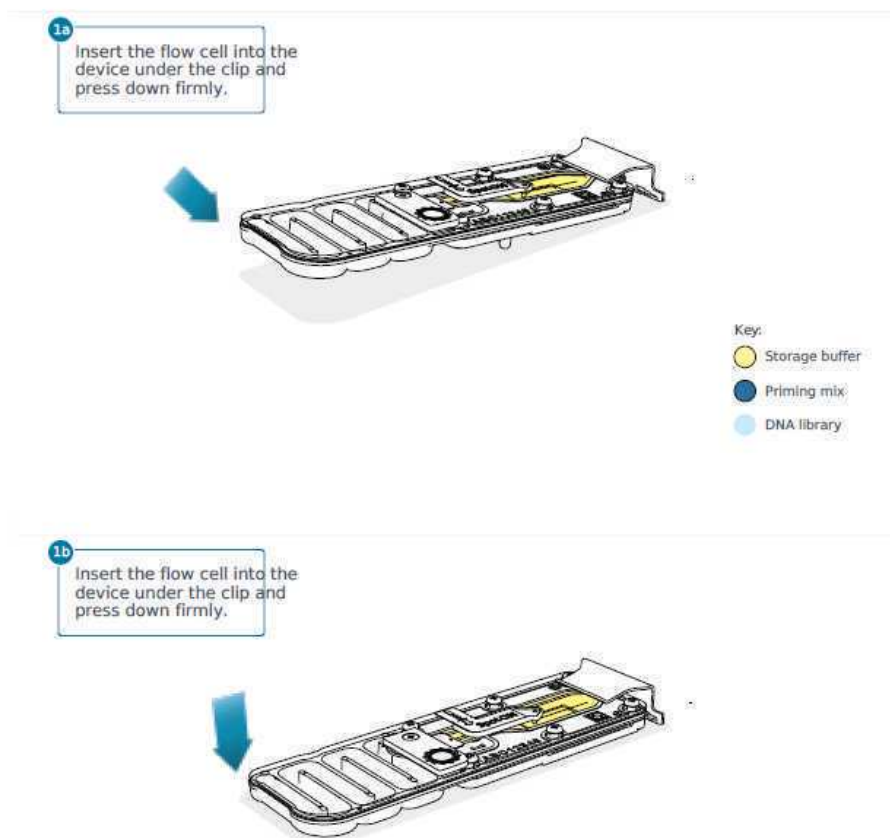

- 4.4 Slide the flow cell priming port cover clockwise to open the priming port. The priming port cover has the Nanopore logo.

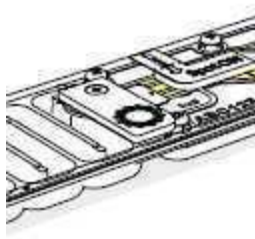

Closed priming port cover.

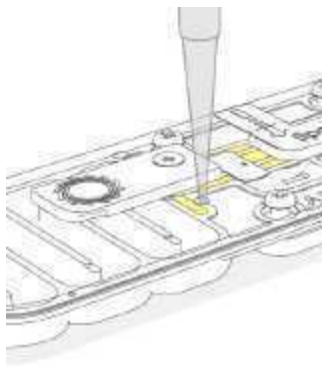

Open priming port cover with pipette entering the priming port.

4.5 IMPORTANT

Take care when drawing back buffer from the flow cell. Do not remove more than 20-30 μl , and make sure that the array of pores are covered by buffer at all times. Introducing air bubbles into the array can irreversibly damage pores.


4.6 After opening the priming port, check for a small air bubble under the cover. Draw back a small volume to remove any bubbles:

1. Set a P1000 pipette to 200 μl
2. Insert the tip into the priming port
3. Turn the wheel until the dial shows 220-230 μl , to draw back 20-30 μl , or until you can see a small volume of buffer entering the pipette tip.

Note: Visually check that there is continuous buffer from the priming port across the sensor array.




3 Insert a P1000 pipette with an empty tip into the **Priming port**. Turn the pipette wheel to draw back 20–30 μl or until you can see a small volume of buffer entering the pipette tip.

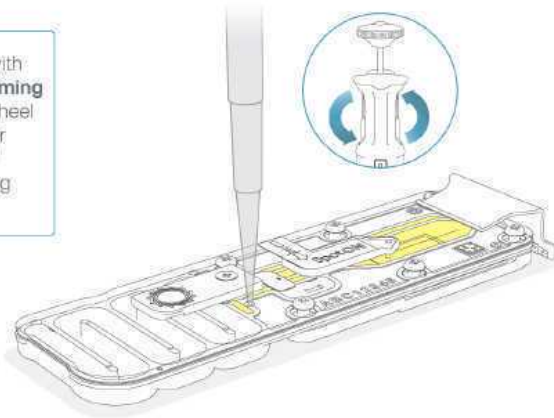

- 4.7 Load 800 μl of the priming mix into the flow cell via the priming port, avoiding the introduction of air bubbles. This is best done by turning the wheel slowly and ejecting the priming buffer slowly into the priming port. Close the priming port cover during incubation. Wait for five minutes.

4 Slowly load 800 μl of the priming mix into the **Priming port**. Ensure there are no air bubbles in the pipette tip.

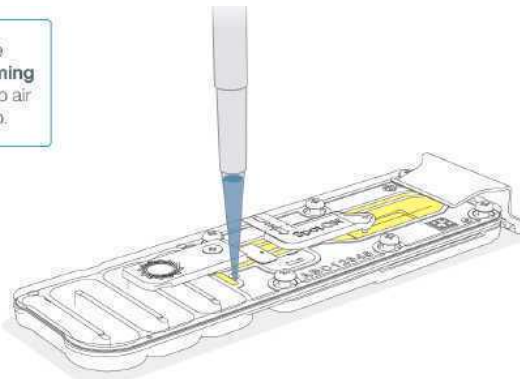

Wait 5 minutes before proceeding to the next step.

- 4.8 The flowcell is now ready for loading and sequencing.

Loading of DNA library onto MinION flowcell

- 5 Now we will load the DNA library from step 3.23 to start sequencing.


This is a useful video for flowcell priming and loading [Community - Priming and loading your flow cell \(nanoporetech.com\)](#).

5.1 IMPORTANT

The Library Beads (LIB) tube contains a suspension of beads. These beads settle very quickly. It is vital that they are mixed immediately before use.

5.2 Complete the flow cell priming:


1. priming port cover.
2. Gently lift the SpotON sample port cover to make the SpotON sample port accessible.
3. Load 200 μ l of the priming mix into the flow cell priming port (not the SpotON sample port), avoiding the introduction of air bubbles.

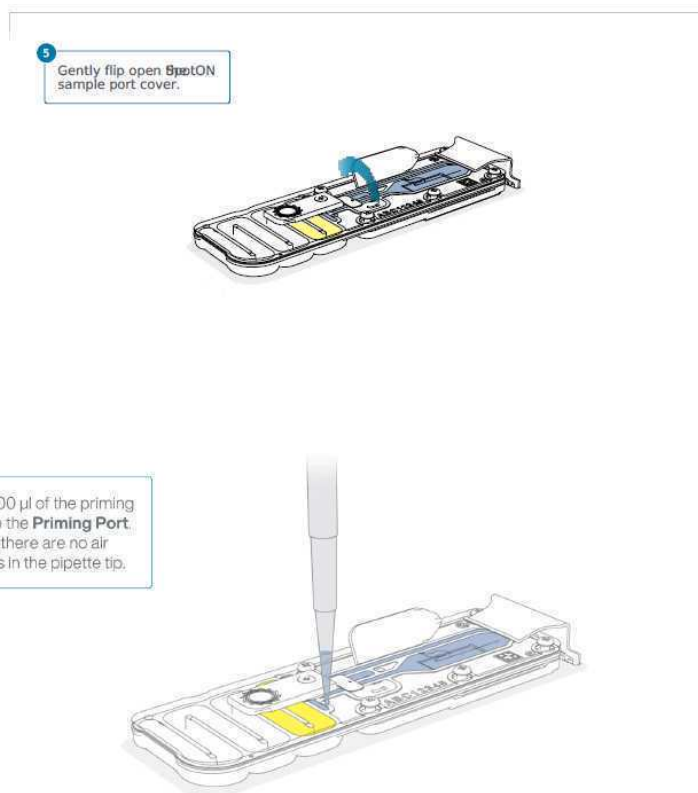

5.3 Mix the prepared library gently by pipetting up and down just prior to loading.


- 5.4 Add 75 μ l of the prepared library to the flow cell via the SpotON sample port in a dropwise fashion. Ensure each drop flows into the port before adding the next.

A demonstrator will assist you with this step.

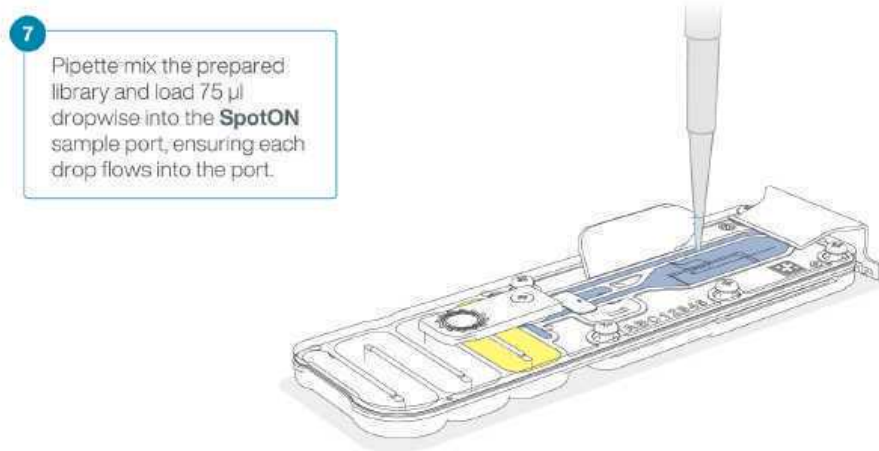


- 5.5 Gently replace the SpotON sample port cover, making sure the bung enters the SpotON port and close the priming port.

8 Gently replace the **SpotON** sample port cover.

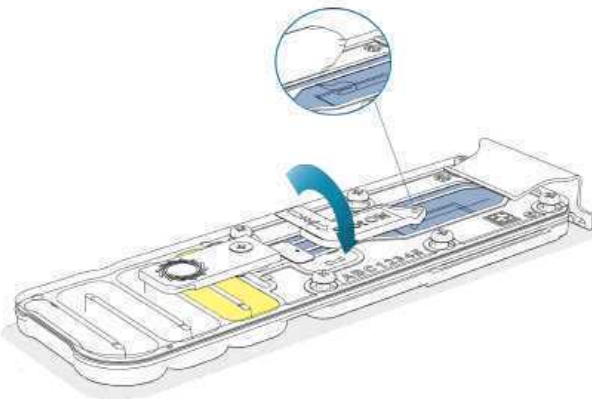

9 Gently close the **Priming port**.

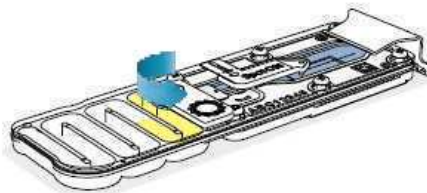

5.6 IMPORTANT

Install the light shield on your flow cell as soon as library has been loaded for optimal sequencing output.

We recommend leaving the light shield on the flow cell when library is loaded, including during any washing and reloading steps.

The shield can be removed when the library has been removed from the flow cell. Gently replace the SpotON sample port cover, making sure the bung enters the SpotON port and close the priming port.

5.7 Place the light shield onto the flow cell, as follows:


1. Carefully place the leading edge of the light shield against the clip.

Note: Do not force the light shield underneath the clip.



2. Gently lower the light shield onto the flow cell. The light shield should sit around the SpotON cover, covering the entire top section of the flow cell.

10
Carefully align the **light shield** against the clip and lower onto the flow cell.

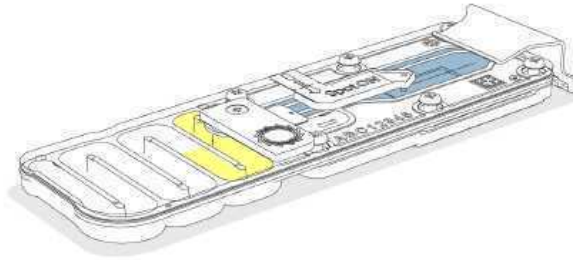

- 5.8 Now we are ready for DNA sequencing via the MinKNOW graphical user interface.

DNA sequencing and basecalling

- 6 The demonstrators will help setting up the sequencing runs. We will perform basecalling later on.

Here is a list of useful videos:

- Introduction to the MinKNOW graphic user interface [Community - MinKNOW: The graphical user interface \(GUI\) \(nanoporetech.com\)](#)
- How to set up a sequencing experiment [Community - MinKNOW: Setting up a new experiment \(nanoporetech.com\)](#)
- How to assess your run is going fine [Community - MinKNOW: Assessing your run \(nanoporetech.com\)](#)
- How MinKNOW selects the best channels [Community - MinKNOW: Mux scan and active channel selection \(nanoporetech.com\)](#)
- More about checking on if your sequencing run is going smoothly [Community - MinKNOW: Pore occupancy \(nanoporetech.com\)](#)

- 6.1 We will share sequencing results and initial analysis at the end of week 6. This way you have plenty of time to play around with the data.



Protocol references

This protocol is based on publicly available protocols from Oxford Nanopore Technology found here [Community - Ligation sequencing amplicons V14 \(SQK-LSK114\)Protocol - Ligation sequencing amplicons V14 \(SQK-LSK114\) \(nanoporetech.com\)](#) and the preprint here [Barcode 100K Specimens: In a Single Nanopore Run | bioRxiv](#).

Further reading on Nanopore ITS analysis can be found

<https://www.sciencedirect.com/science/article/pii/S1087184525000830?via%3Dihub>.