

Sep 01, 2022

Version 2

🌐 Setting a sequencing run with a nanopore MinION and the Rapid Sequencing gDNA kit (SQK-RAD004) V.2

DOI

dx.doi.org/10.17504/protocols.io.4r3l2oy53v1y/v2

Narjol Gonzalez-Escalona¹


¹Food and Drug Administration

Narjol GE - CFSAN



Narjol Gonzalez-Escalona

Food and Drug Administration

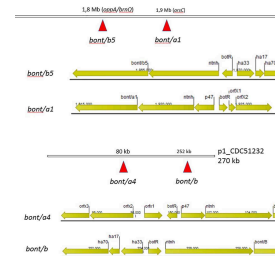

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.4r3l2oy53v1y/v2>

External link: <https://community.nanoporetech.com/>



Protocol Citation: Narjol Gonzalez-Escalona 2022. Setting a sequencing run with a nanopore MinION and the Rapid Sequencing gDNA kit (SQK-RAD004). **protocols.io** <https://dx.doi.org/10.17504/protocols.io.4r3l2oy53v1y/v2> Version created by **Narjol Gonzalez-Escalona**

Manuscript citation:

Closed genome sequences of two Clostridium botulinum strains obtained by nanopore sequencing N Gonzalez-Escalona, J Haendiges, JD Miller, SK Sharma Microbiology Resource Announcements 7 (9), e01075-18

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: September 01, 2022

Last Modified: September 01, 2022

Protocol Integer ID: 69415

Keywords: MinION, setting up a NANopore run, rapid sequencing kit, RAD004, SQK-RAD004 , rapid sequencing kit from nanopore, rapid sequencing gdna kit, nanopore minion, sequencing kit, nanopore, sequencing run, sequencing, minion, rapid preparation time, rad004, gdna, preparation time, kit

Abstract

This protocol is to help in setting up a MinION sequencing run using the rapid sequencing kit from Nanopore (SQK-RAD004). It contains all steps and material needed for a successful run. This kit is recommended for rapid preparation time and limited laboratory equipment. The output is lower than that generated by the ligation sequencing kit.

preparation time approximately 15 minutes and require a minimum of 400 ng per reaction.

Guidelines

The Rapid Sequencing Kit contains enough reagents to perform six sequencing runs:

Box 1: Contains all reagents for the sequencing step (store at -20C)

Box 2: Contains all the reagents for flowcell priming (store at -20C) - Flush Buffer (FB – blue caps) and Flush tether (FLT - purple cap)



Materials

Equipment:

- MinION device (MK1B)
- PC with 64 bit Windows 10 or 64 bit Linux (Intel i7, i9, Xeon or better, 16GB RAM, NVIDIA GPU at least 16GB GPU memory, 1T SSD)
- High speed internet high speed or ethernet connection with access to the internet
- Thermocycler
- Microcentrifuge

Supplies:

- SQK-RAD004 kit
- Flowcell (9.4.1 -lower accuracy or 10.4.1 - high accuracy)
- Pipette Tips, sterile, filtered (assorted volumes)
- Microcentrifuge for 0.2 ml and 0.5 ml tubes
- Tubes, 0.2 ml, sterile (ThermoFisher cat# AB0337 or equivalent)
- Nuclease-free water (e.g. ThermoFisher, cat # AM9937)
- Timer

Troubleshooting













Before start

Before using this protocol:

1. **Make sure that you have extracted the DNA of the bacteria or other organism** (if long reads needed, I suggest using a tested DNA extraction kit or in house DNA extraction that renders high quality less sheared DNA).
2. **Quantify the DNA** (I suggest using the Qubit Fluorometer for determining the DNA concentration).
3. **Check that the PC to be used for the run is compatible with the process.**
4. **Make sure that the MinKnow Software version** (software installed in the PC to run the MinION) **is up to date** (search the community and download the latest version if the version on the PC is not up to date).
5. **Make sure that you have all supplies.**




Preparation of the DNA Library:

- 1
Thaw all reagents in box 1 and 2 of the RAD004 kit.
- 2
Take a flow cell from  4-8 °C and leave it at RT.
- 3
Prepare the DNA to the suggested concentration of  400 ng total in  7.5 µL, or  54 µL in a  0.2 mL sterile thin-walled PCR tube. Adjust the volume using nuclease-free water. Mix by gently flicking the tube 5 times. DO NOT vortex. Spin down briefly in a micro-centrifuge.
- 4
Add  2.5 µL of FRA to the tube containing the prepared DNA from the previous step. Mix by gently flicking the tube as above and spin down. The FRA is a transposase that
- 5
Incubate in Thermocycler at  30 °C for  00:01:00 1 minute, then at  80 °C for  00:01:00. Let it cool at RT. 2m
- 6
Add  1 µL of RAP to the tube. Mix gently by flicking 5 times and spin down. Incubate for  00:05:00. After this step your library DNA is ready for mixing with the sequencing reagents and loading to the flow cell. While this step is running you can start setting up the priming of the flow cell. 5m

Flow cell priming:

- 7
Flow cell priming:


Mix FB and FLT (independently) by vortexing and spin down. Then add  30 µL of FLT to 1 tube of FB and mix by pipetting thoroughly.





8 Place the RT flow cell in the MinION device and proceed with the flow cell QC by clicking on flow cell check and wait till the program reports back how many live pores are in the flow cell. If live pores below 800, do not continue sequencing (if very important samples) . Contact the vendor if the flow cell is still under warranty for a replacement.



8.1

9 If flow cell passed QC then open the priming port by sliding the priming port cover clockwise.

9.1 By inserting a pipette tip on the priming port, remove approximately  30 μL avoiding introducing bubbles.

9.2 Add  800 μL of the FB+FLT solution prepared on step 7 to the priming port, avoiding introducing bubbles, by slowly dialing down the pipette. Wait  00:05:00 .




5m

9.3 After the  00:05:00 , gently open the flow cell SpotON cover and expose the loading port. Add  200 μL of FB+FLT to the priming port (NOT the SpotON port). Make sure that the loading port does not have any bubble.

5m


Preparation of the sequencing Library:

10 To the tube containing the DNA library add the following:

-  34 μL of the sequencing buffer (SQB) - red
-  25.5 μL of loading beads (LB) - pink (Thoroughly mixed before use)
-  4.5 μL of nuclease-free water

Mix thoroughly by pipetting and make sure that the solution looks homogenous. This is the sequencing reaction ready to be loaded to the flow cell.

Loading the flow cell

11 Add the  75 μL of the sequencing library to the SpotON port in a dropwise fashion. Make sure that every drop falls into the port and that it goes in, before adding the next




drop.

Replace the cover over the SpotON and close the priming lid port, and close the MinION device lid

Sequencing

2d

- 12 On the Minknow software select Start, then Start sequencing and follow the instructions for setting the sequencing run. Experiment name (I use the same as the flow cell serial number), then select the kit (RAD004), The run options I leave it mostly default except that I change the time to  48:00:00 , the basecalling to fast (or high accuracy, your choice), Leave the output default (you can point it to wherever you want the run to be save, then just click start at the end and your sequencing run will start.

2d