



Jul 25, 2025

Version 2

Bioloc3D (B3D) : Protocol for 3D automated analysis of colocalization in microscopy V.2

DOI

<https://dx.doi.org/10.17504/protocols.io.kxygx4p5zl8j/v2>



Thibault Dhellemmes¹, Rémi inet¹, Fabrice Cordelières², Marc Landry¹

¹Université de Bordeaux - Institut des Maladies Neurodégénératives (IMN) / UMR 5293;

²Bordeaux Imaging Center

Bioloc3D (B3D)



Thibault Dhellemmes

Institut des Maladies Neurodégénératives (IMN)

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.kxygx4p5zl8j/v2>

External link: <https://github.com/Bioloc3D>

Protocol Citation: Thibault Dhellemmes, Rémi inet, Fabrice Cordelières, Marc Landry 2025. Bioloc3D (B3D) : Protocol for 3D automated analysis of colocalization in microscopy. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.kxygx4p5zl8j/v2> Version created by **Thibault Dhellemmes**

License: This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: July 25, 2025

Last Modified: July 29, 2025

Protocol Integer ID: 223277

Keywords: Software, Microscopy, Automatization, Colocalization, Fluorescent, Quantification, Automated analysis, Statistics, bioloc3d, tools capable of colocalization quantification, microscopy this protocol, microscopy, fluorescent labeling across entire stack, fluorescent labeling, colocalization quantification, analysis of colocalization, colocalization, b3d, automate statistical analysis, automated analysis

Funders Acknowledgements:

ANR

Grant ID: ANR-19-CE37-0019 Relax

FRM

Grant ID: FDT202404018302

Abstract

This protocol introduces Bioloc3D, a fully open-source tool designed to quantify fluorescent labeling across entire stacks and automate statistical analysis. It offers a solution to the lack of tools capable of colocalization quantification between multiple identified elements in 3D. Expected results includes graphical representations of the results through images and estimation plots.

Image Attribution

Thibault Dhellemmes

Materials

Raw image stacks with fluorescent signals to analyze in any format adapted to ImageJ.

Linux/Mac/Windows computer with ImageJ, Excel and Bioloc3D-Plt installed. ImageJ installation file is available online at <https://imagej.net>.


Before start

To operate in the best condition, the minimal requirement for colocalization analysis is a stack of 2 channels with 2 images (i); in 8-bit color depth (ii); preferably in TIFF format (iii)


B3D workflow and protocol

- 1
 - Schematic representation of the analysis workflow using B3D
 - Diagram of the simplified procedure
 - User-interface (UI) of *B3D-Img* toolset in ImageJ and *B3D-Plt* standalone application

B3D workflow

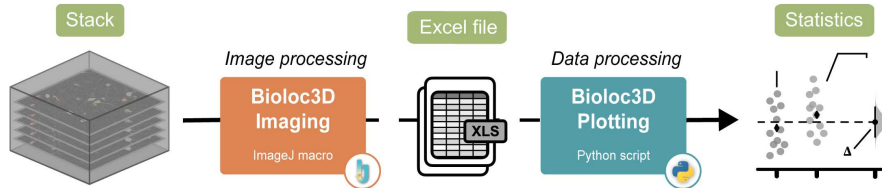

B3D simplified protocol & UI

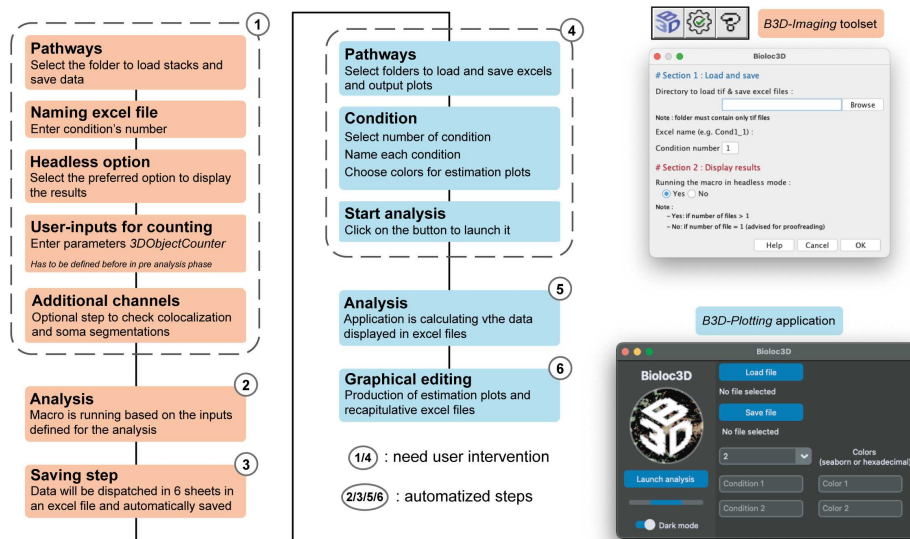

Illustrations of B3D workflow and protocol

- 2
 - **! CAUTION** Figures mentioned under brackets (e.g. (*Supplementary Fig. 1C*)) refer to the method publication in which the experimental design of B3D toolset is detailed (cf. "Manuscript citation" section)

Installation

3m

- 3 Download the *B3D-Imaging (B3D-Img)* ImageJ macro at <https://github.com/Bioloc3D/B3D-Img>
- 4 Install *B3D-Img* in ImageJ macro library, by drag/drop the .ijm file in the toolsets folder of the application (*Fiji/ImageJ* → *macros* → *toolsets*)
- 5 Install latest version of *3D Objects Counter* plugin available online on GitHub, by drag/drop the .jar file in the plugins folder of the application (*Fiji/ImageJ* → *plugins*)
- 6 Launch ImageJ software

Software

ImageJ (Fiji)	NAME
Windows 10	OS
National Institutes of Health (USA)	DEVELOPER
https://imagej.net/software/fiji/downloads	REPOSITORY
https://imagej.net/software/fiji/downloads	SOURCE LINK

- 7 Through ImageJ install *MorphoLibJ*, *3D ImageJ Suite* and *Read and Write Excel* additional plugins available in plugin list (*Help* → *Update...* → *Manage Update Sites* → *tick plugins* → *Apply and Close*)
- 8 Close ImageJ software
- 9 Download the *B3D-Plt* application depending of the user operating system at <https://github.com/Bioloc3D/B3D-Plt>

Software

Bioloc3D

NAME

Windows / MacOS

OS

Thibault Dhellemmes / Rémi Kinet

DEVELOPER

<https://github.com/Bioloc3D>

REPOSITORY

<https://github.com/Bioloc3D>

SOURCE LINK

Image analysis : *B3D-Imaging (B3D-Img)*

10 Launch B3D toolset in ImageJ software

30s

10.1 Start ImageJ software

10.2 Click on the "More Tools" menu (>>) and select *B3D-Img*

11 Defining the correct parameters

2m 30s



11.1 Click on B3D parameters icon (cf. Supplementary Fig. 1A)

11.2 Select the stack to test the parameters and the folder in which test data will be saved. Parameters must be tested in an iterative matter, select one channel from the stack (cf. Supplementary Fig. 1C, Step 1 of the diagram).



11.3 Click "OK"

11.4 Define the threshold adapted to the signal of the selected channel (cf. Supplementary Fig. 1C, Step 2 of the diagram).



11.5 Click "OK"

11.6 Use the threshold defined in sub-step 11.4 and define minimal and maximal size threshold to test (Supplementary Fig. 1C, Step 3 of the diagram). If needed, define minimal area for soma detection (cf. Supplementary Fig. 1C, Step 4 of the diagram)

11.7 Click "OK": Test analysis is launched

12 Analysis of one stack and raw data export

! CAUTION As discussed in the "advantages and limitations" section, time of analysis depends on several aspects, including computation power of the machine, size, potential number of objects to count per channel and signal-to-noise ratio of your images. To estimate the time of each step of the procedure, we used a stack presenting the following characteristics: size= 1024 × 1024; z= 25; objects counted per channel= 1000. In addition, sample is presenting a moderate signal-to-background (SBR) but a low signal-to-noise ratio (SNR).

12.1 Click on B3D analysis icon (cf. Fig. 2A)

12.2 Select the folder containing stack(s) to analyze and numbered your batch of analysis (Fig2C, section 1 on the diagram). Normalization of the name of exported data files is compulsory for the automatic computation. If more than one file is present in this folder, "headless" option must be defined on "Yes"

! CAUTION Select folder must contains only tiff file(s) before starting the analysis in ImageJ.

12.3 Click "OK"

12.4 Enter parameters defined in step 11, in the GUI (Fig. 2C, section 3 on the diagram) and select options for quantification based on expected information

12.5 Option A :

Simple quantification (3D-Simple) of primary channels, without quantification of additional features (control of colocalizations sites and soma)

(i) Select "Yes" for "Simple segmentation of C1 & C2:"



30s



(ii) Select "No" for "Check colocalization with an additional channel"

(iii) Select "No" for "Complete counting with an additional channel"

12.6 **Option B :**

Complete quantification (3D-OC) of primary channels, without quantification of additional features

(i) Select "No" for "Simple segmentation of C1 & C2:"

(ii) Select "No" for "Check colocalization with an additional channel"

(iii) Select "No" for "Complete counting with an additional channel"

3m

12.7 **Option C :**

Simple quantification (3D-Simple) of primary channels, including quantification of additional features

(i) Select "Yes" for "Simple segmentation of C1 & C2:"

(ii) Select "Yes" for "Check colocalization with an additional channel"

(iii) Enter parameters for colocalization control (threshold and size)

(iv) Select "Yes" for "Complete counting with an additional channel"

(v) Enter parameters for soma counting (minimum area and threshold)

40s

12.8 **Option D :**

Complete quantification (3D-OC) of primary channels, including quantification of additional features

(i) Select "No" for "Simple segmentation of C1 & C2:"

(ii) Select "Yes" for "Check colocalization with an additional channel"

(iii) Enter parameters for colocalization control (threshold and size)

(iv) Select "Yes" for "Complete counting with an additional channel"

(v) Enter parameters for soma counting (minimum area and threshold)

3m 30s

12.9 Click "OK"

12.10 (Optional) Click "Help" to be redirected to GitHub repository

13 **(Optional) Proofreading of 3D quantification on one stack**

5m

! CAUTION If only one stack is present in the folder, illustrations of the results can be displayed by selecting "No" for headless mode (Fig2C, section 2 on the diagram).

13.1 For a more precise proofreading apply option B or D (Sub-steps 12.6 or 12.8)



- 13.2 Overview the different output stacks with the excel of the data aside to check if it counts what you wanted

Computation, plotting and statistical analysis : *B3D-Plotting (B3D-Plt)*

- 14 Open the B3D-Plt application
- 15 Click on "Load file" and select the folder containing the Excels to analyze 
- 16 Click on "Save file" and select the folder in which graphs and analysis will be saved 
- 17 Click on the drop-down menu and select the number of conditions you have to analyze 
- 18 In the bars that appear, enter the name of the condition and the color you want to associate with it (seaborn or hexadecimal #FF0000 style, but can stay empty)
- Seaborn color link : https://seaborn.pydata.org/tutorial/color_palettes.html
- 19 Click on "Launch analysis" 15m
- ! CAUTION** Timing indicated for this step has been tested on a file containing 60 excel files 
- 20 The name of the document currently analyzed is print to monitor progress
- 21 Analysis end is notified through a pop-up window
- 22 Find the figures and Excel in the folder selected in step 16