

Jan 13, 2025

# Characterization of the Photon Conversion Factor, Noise, and Dynamic Range of Light Microscope Detection Systems



DOI

<https://dx.doi.org/10.17504/protocols.io.14egn61pyl5d/v1>

David McFadden<sup>1</sup>, Laszlo Barna<sup>2</sup>, Luis-Francisco Acevedo-Hueso<sup>3</sup>, Sergiy V Avilov<sup>4</sup>, Gert-Jan Bakker<sup>5</sup>, Sebastian Beer<sup>6</sup>, Ivan Belyaev<sup>7</sup>, Valeria Berno<sup>8</sup>, Mariana T Carvalho<sup>9</sup>, Yann Cesbron<sup>10</sup>, Catalin Chiritescu<sup>11</sup>, Orestis Faklaris<sup>12</sup>, Nathalie Gaudreault<sup>13</sup>, David Grunwald<sup>14</sup>, Thomas Guilbert<sup>15</sup>, Mathias Hammer<sup>16</sup>, Rainer Heintzmann<sup>17</sup>, Gerhard Holst<sup>18</sup>, Ayse Aslihan Koksoy<sup>19</sup>, Gabriel Martins<sup>20</sup>, Glyn Nelson<sup>21</sup>, Roland Nitschke<sup>22</sup>, Alex Payne-Dwyer<sup>23</sup>, Santosh Podder<sup>24</sup>, Sathya Srinivasan<sup>25</sup>, Kees van der Oord<sup>26</sup>, Andre Zeug<sup>27</sup>, Britta Schroth-Diez<sup>28</sup>

<sup>1</sup>Friedrich Schiller University Jena; <sup>2</sup>Institute of Experimental Medicine, Budapest;

<sup>3</sup>Neuro-AI Monitoring and Quantum Medical Diagnostics;

<sup>4</sup>Max Planck Institute of Immunobiology and Epigenetics;

<sup>5</sup>Radboud University Medical Center, Medical BioSciences department, Nijmegen, the Netherlands;

<sup>6</sup>Hamamatsu Photonics Deutschland GmbH; <sup>7</sup>The Microverse cluster, Friedrich Schiller University Jena;

<sup>8</sup>ALEMBIC Advanced Light and Electron Microscopy BiImaging Center, San Raffaele Scientific Institute;

<sup>9</sup>INL - International Iberian Nanotechnology Laboratory; <sup>10</sup>Institute of Science and Technology Austria;

<sup>11</sup>Phi Optics; <sup>12</sup>MRI, Biocampus, University of Montpellier, CNRS, INSERM, Montpellier, France;

<sup>13</sup>Allen Institute for Cell Science, Seattle, WA, USA; <sup>14</sup>UMass Chan Medical School;

<sup>15</sup>Institut Cochin, INSERM U1016, CNRS UMR 8104, Université Paris Cité, Paris, France;

<sup>16</sup>UMass Chan Medical School, RNA Therapeutics Institute, Worcester, MA, USA;

<sup>17</sup>Leibniz Institute of Photonic Technology; <sup>18</sup>Excelitas PCO GmbH; <sup>19</sup>MD Anderson Cancer Center;

<sup>20</sup>Gulbenkian Institute for Molecular Medicine; <sup>21</sup>Newcastle University, BiImaging Unit, UK;

<sup>22</sup>Life Imaging Center, University of Freiburg, Germany;

<sup>23</sup>School of Physics Engineering and Technology, University of York, UK;

<sup>24</sup>Indian Institute of Science Education and Research (IISER) Pune; <sup>25</sup>ONPRC; <sup>26</sup>Nikon Europe BV;

<sup>27</sup>Cellular Neurophysiology, Hannover Medical School, Hannover, Germany;

<sup>28</sup>Max Planck Institute of Molecular Cell Biology and Genetics, Dresden, Germany

QUAREP-LiMi

Tech. support email: [info@quarep.org](mailto:info@quarep.org)



**Britta Schroth-Diez**

Max Planck Institute of Molecular Cell Biology and Genetics,...

**Create & collaborate more with a free account**

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

[Create free account](#)

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.14egn61pyl5d/v1>

External link: <https://quarep.org/>

**Collection Citation:** David McFadden, Laszlo Barna, Luis-Francisco Acevedo-Hueso, Sergiy V Avilov, Gert-Jan Bakker, Sebastian Beer, Ivan Belyaev, Valeria Berno, Mariana T Carvalho, Yann Cesbron, Catalin Chiritescu, Orestis Faklaris, Nathalie Gaudreault, David Grunwald, Thomas Guilbert, Mathias Hammer, Rainer Heintzmann, Gerhard Holst, Ayse Aslihan Koksoy, Gabriel Martins, Glyn Nelson, Roland Nitschke, Alex Payne-Dwyer, Santosh Podder, Sathya Srinivasan, Kees van der Oord, Andre Zeug, Britta Schroth-Diez 2025. Characterization of the Photon Conversion Factor, Noise, and Dynamic Range of Light Microscope Detection Systems. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.14egn61pyl5d/v1>

**License:** This is an open access collection distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**We use this collection and it's working**

**Created:** September 08, 2024

**Last Modified:** September 01, 2025

**Collection Integer ID:** 107127

**Keywords:** protocol collection, detection system, characterization, sample preparation, in-homogeneous, photon conversion factor, dark noise, dynamic range, light microscope, photon transfer method, quality control, reproducibility, light microscopy detection system, light microscope detection system, monitoring of microscope performance, dynamic range of light microscope detection system, microscope performance, photon conversation factor, detector, area detector, point detector, detection setting, detection system, experiment qc, imaging scientist, homogeneous detector illumination scheme, characterization of the photon conversion factor, using detection setting, suitable sample for the acquisition, photon conversion factor, detection, pcf photo, uniform illumination, monitoring, full characterization of detection path performance, respective analysis protocol, instrument qc, detection path performance, protocol collection, image quality, protocols in this collection, aims this protocol, experiment, Quarep

## Disclaimer

This protocol collection was developed by members of **WG2 "Detection System Performance"** of the international consortium **QUAREP-LiMi**.

The Consortium for Quality Assessment and Reproducibility for Instruments and Images in Light Microscopy (QUAREP-LiMi), formed by the global community of practitioners, researchers, developers, service providers, funders, publishers, policy makers and industry related to the use of light microscopy, is committed to democratizing access to quantitative and reproducible light microscopy and the data generated by it.

This protocol collection has undergone the internal approval process of QUAREP-LiMi.

It is not the intention of this protocol collection to supplant the need for independent professional judgement, advice, diagnosis or treatment. Any action taken or withheld on the basis of the information presented here is undertaken at the user's own risk. The user agrees that neither the company nor any of the authors, contributors, administrators, or other individuals associated with protocols.io or QUAREP-LiMi can be held liable for any injuries, damages or losses incurred as a result of the user's use of the information contained in or linked to this protocol or any of our websites, applications, or services.

## Abstract

The protocols in this collection describe how to measure and analyze the photon conversation factor (PCF photo-electrons/count), readnoise, and dynamic range of a light microscopy detection system; which can either be a point detector or an area detector, using an in-homogeneous detector illumination scheme (as opposed to uniform illumination). The collection includes protocols on how to prepare a suitable sample for the acquisition, how to acquire data, as well as a respective analysis protocol.


There are three aims a detection system can be characterized for, briefly:

**Aim 1 - experiment QC:** Characterize the microscope performance using detection settings that match the experiment.

**Aim 2 - instrument QC.:** Monitoring of microscope performance over time for service purposes, to maintain image quality constant and at high level.

**Aim 3 - system characterization:** Full characterization of detection path performance under the range of settings applied by users.

For a more detailed description refer to protocol 1. Introduction - Background and Aims.

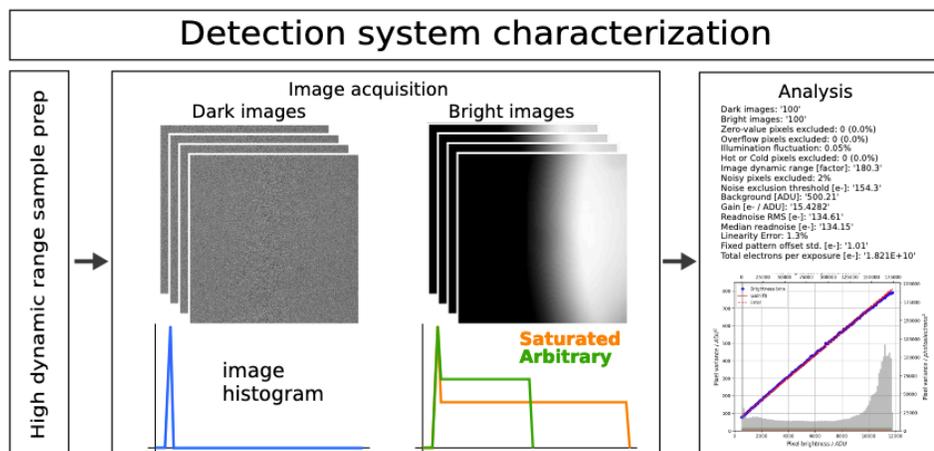

*Workflow of detection system characterization using the Light Microscopy Detection System Characterization toolset.*

This protocol collection represents the collective experience of over 100 imaging scientists and industry experts. Measurements made by our working group with these protocols will be available in a public database.

**Please note**, that this is an evolving document, to be versioned and updated, based on community feedback and new data.

## Materials

For materials please refer to the respective protocol.

## Troubleshooting

### Safety warnings

- ❗ Ensure you follow general lab safety guidelines for radiation sources and chemicals as outlined within your organisation.



#### Laser safety and regulations

- Please refer to the documentation provided by the manufacturer for additional warnings and preventive, protective equipment (PPE) requirements (e.g. laser safety goggles). Always consult your local Laser Safety Officer or Radiation Safety Officer and refer to your laboratory safety documentation for more information.
- You can also consult your Laser Safety Standards ANSI Z136 in North America, SUVA 66049.D in Europe, and BS EN 60825-1 in the UK. Additionally, laser safety standards and regulations are covered by IEC norm 60825-1, and LED eye safety standards and regulations are covered by IEC norm 62471 in Europe.
- Safety guideline: Hazardous, visible, or invisible radiation from lasers, lamps, and other light sources used for microscopy **can cause permanent damage to the retina, skin burns, and fire**. Always follow proper laser safety protocols for your equipment and situation.

## Before start

For a detailed introduction and description, please read [protocol 1](#).

For characterization of an area detector (e.g., sCMOS camera), please follow [protocol 2](#) , [protocol 3](#), and [protocol 5](#).

For characterization of a point detector (e.g., PMT or HyD detector), please follow [protocol 2](#), [protocol 4](#) and [protocol 5](#).

Make sure the microscope is turned on at least one hour prior to the measurements to allow the system to stabilize.

The test sample can be prepared well in advance.

## Disclaimer

This protocol collection was developed by members of WG2 "Detection System Performance" of the international consortium QUAREP-LiMi.

The Consortium for Quality Assessment and Reproducibility for Instruments and Images in Light Microscopy (QUAREP-LiMi), formed by the global community of practitioners, researchers, developers, service providers, funders, publishers, policy makers and industry related to the use of light microscopy, is committed to democratizing access to quantitative and reproducible light microscopy and the data generated by it.

This protocol collection has undergone the internal approval process of QUAREP-LiMi.

It is not the intention of this protocol collection to supplant the need for independent professional judgement, advice, diagnosis or treatment. Any action taken or withheld on the basis of the information presented here is undertaken at the user's own risk. The user agrees that neither the company nor any of the authors, contributors, administrators, or other individuals associated with protocols.io or QUAREP-LiMi can be held liable for any injuries, damages or losses incurred as a result of the user's use of the information contained in or linked to this protocol or any of our websites, applications, or services.

## Files

Q SEARCH

### Protocol



NAME

#### 1. Introduction - Background and Aims

VERSION 1

CREATED BY



**Britta Schroth-Diez**

Max Planck Institute of Molecular Cell Biology and Genetics, Dresden, Germany

OPEN →

### Protocol



NAME

#### 2. Sample Preparation - An Easy-to-Prepare Sample Slide

VERSION 1

CREATED BY



**Britta Schroth-Diez**

Max Planck Institute of Molecular Cell Biology and Genetics, Dresden, Germany

OPEN →

### Protocol



NAME

#### 3. Data Generation - Systems with an Area Detector

VERSION 1

CREATED BY



**Britta Schroth-Diez**

Max Planck Institute of Molecular Cell Biology and Genetics, Dresden, Germany

OPEN →

### Protocol



NAME

#### 4. Data Generation - Systems with a Point Detector

VERSION 1

CREATED BY



**Britta Schroth-Diez**  
Max Planck Institute of Molecular Cell Biology and Genetics, Dresden, Germany

OPEN →

## Protocol



NAME

**5. Analysis - Characterization of the Photon Conversion Factor, Noise, and Dynamic Range**

VERSION 1

CREATED BY



**Britta Schroth-Diez**  
Max Planck Institute of Molecular Cell Biology and Genetics, Dresden, Germany

OPEN →

## Protocol references

- Art J. 1995. Photon detectors for confocal microscopy. In: Pawley J, editor. Handbook of Biological Confocal Microscopy. New York, NY: Springer. pp. 183–195.
- Cho EH, Lockett SJ. 2006. Calibration and standardization of the emission light path of confocal microscopes. *Journal of Microscopy* **223**:15–25. doi:10.1111/j.1365-2818.2006.01598.x
- Deagle RC, Wee T-L (Erika), Brown CM. 2017. Reproducibility in light microscopy: Maintenance, standards and SOPs. *The International Journal of Biochemistry & Cell Biology* **89**:120–124. doi:10.1016/j.biocel.2017.06.008
- EMVA 1288. n.d. . *European Machine Vision Association*. <https://www.emva.org/standards-technology/emva-1288/>
- Heintzmann R, Relich PK, Nieuwenhuizen RPJ, Lidke KA, Rieger B. 2016. Calibrating photon counts from a single image. doi:10.48550/ARXIV.1611.05654
- Janesick JR. 2007. Photon Transfer. SPIE. doi:10.1117/3.725073
- McFadden D. 2022. gui calibration tool. *Releases · bionanoimaging, NanoImagingPack*. <https://github.com/bionanoimaging/NanoImagingPack/releases>
- McFadden D, Amos B, Heintzmann R. 2022. Quality control of image sensors using gaseous tritium light sources. *Phil Trans R Soc A* **380**:20210130. doi:10.1098/rsta.2021.0130
- Mullikin JC, Van Vliet LJ, Netten H, Boddeke FR, Van Der Feltz G, Young IT. 1994. Methods for CCD camera characterization In: Titus HC, Waks A, editors. Presented at the IS&T/SPIE 1994 International Symposium on Electronic Imaging: Science and Technology. San Jose, CA. pp. 73–84. doi:10.1117/12.175165
- Murray JM. 2013. Practical Aspects of Quantitative Confocal Microscopy Methods in Cell Biology. Elsevier. pp. 427–440. doi:10.1016/B978-0-12-407761-4.00018-X
- Olevsko I, Szederkenyi K, Corridon J, Au A, Delhomme B, Bastien T, Fernandes J, Yip C, Oheim M, Salomon A. 2021. A simple, inexpensive and multi-scale 3-D fluorescent test sample for optical sectioning microscopies. *Microscopy Res & Technique* **84**:2625–2635. doi:10.1002/jemt.23813
- Plakhotnik T, Chennu A, Zvyagin AV. 2006. Statistics of single-electron signals in electron-multiplying charge-coupled devices. *IEEE Trans Electron Devices* **53**:618–622. doi:10.1109/TED.2006.870572
- QUAREP-LiMi. n.d. Quality Assessment and Reproducibility for Instruments & Images in Light Microscopy. <https://quarep.org/>
- Ryan DP, Dunlap MK, Gelfand MP, Werner JH, Van Orden AK, Goodwin PM. 2021. A gain series method for

accurate EMCCD calibration. *Sci Rep* **11**:18348. doi:10.1038/s41598-021-97759-6

van Vliet L, Sudar D, Young I. 1998. Digital fluorescence imaging using cooled charge-coupled device array cameras In: Celis JE, editor. *Cell Biology: A Laboratory Handbook*. Vol. 3. New York: Acad. Press. pp. 109–120.

## Acknowledgements



CRediT\_author\_contributions\_Colle...



Funding acknowledgement.pdf 50KB