




Jun 01, 2023

CHROMagar Malassezia Medium (CHROM)

 In 1 collection

DOI

<https://dx.doi.org/10.17504/protocols.io.6qpvre38plmk/v1>

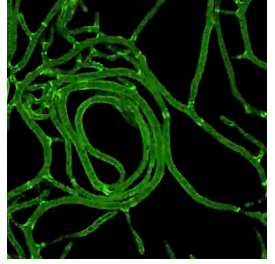

Amy Gladfelter¹

¹Duke University



Amy Gladfelter

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.6qpvre38plmk/v1>

Protocol Citation: Amy Gladfelter 2023. CHROMagar Malassezia Medium (CHROM). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.6qpvre38plmk/v1>

License: This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: November 03, 2018

Last Modified: June 01, 2023

Protocol Integer ID: 17458

Keywords: chromagar malassezia medium, growth of malassezia species, malassezia species, malassezia

Abstract

Media recipe for growth of Malassezia species.

Attachments



Malassezia media 1.p...

33KB



- 1 CHROMagar Malassezia basal medium  56.3 g

Note

CHROMagar, Paris, France

- 2 Tween 40  10 mL

- 3 dH₂O  1 L

- 4 Sterilize by autoclaving.