

Sep 21, 2024

DeepSlice Automated Alignment

DOI

<https://dx.doi.org/10.17504/protocols.io.3byl4946ogo5/v1>

Michael Henderson¹

¹Van Andel Institute



Michael Henderson

Van Andel Research Institute

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.3byl4946ogo5/v1>

Protocol Citation: Michael Henderson 2024. DeepSlice Automated Alignment. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.3byl4946ogo5/v1>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: September 21, 2024



Last Modified: September 25, 2024

Protocol Integer ID: 108117

Keywords: deepslice automated alignment this protocol, deepslice in the quint workflow, deepslice automated alignment, quint workflow, deepslice, alignment, workflow, quint

Funders Acknowledgements:

Aligning Science Across Parkinson's

Abstract


This protocol explains how to use DeepSlice in the Quint Workflow

Before start

DeepSlice is a deep neural network that automatically aligns mouse histology images through the Allen Brain Atlas coordinate framework. Alignments are viewable and refinable in QuickNII and set up sections to give a good starting point near and around each section's correct plane. DeepSlice's alignment is not completely accurate and further fine tuning in QuickNII is necessary. This is a newly developed tool and not imperative for the workflow but helps speed the process of registration and substitutes the FileBuilder step.

DeepSlice Automated Alignment

- 1 Open **DeepSlice** in your web browser.

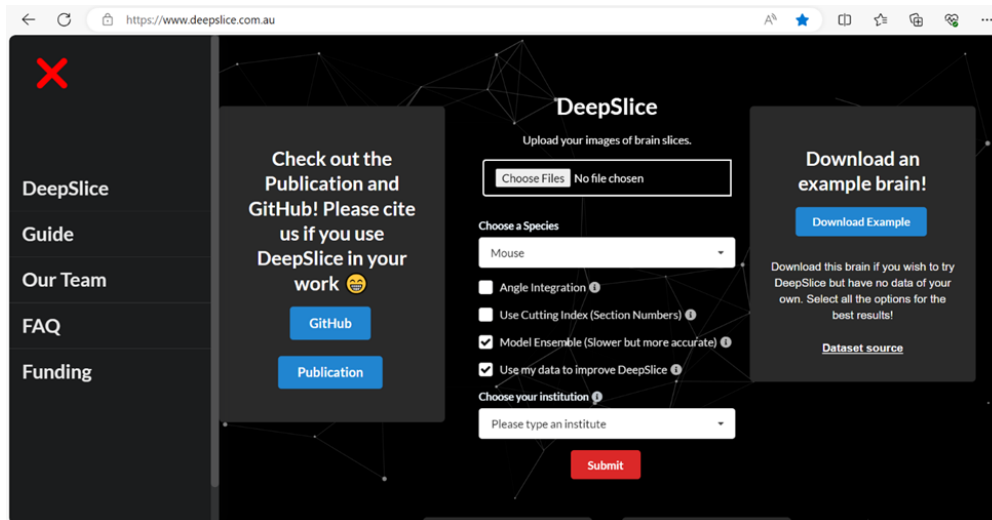

- 2 Select *Choose Files*.
 - a. Upload all images for registration from the **QVN** folder
- 3 Ensure *Mouse* is selected for species, and *Model Ensemble* is checked (uses two DeepSlice versions to optimize alignment).
 - a. Avoid checking *Angle Integration* (this aligns all your brain sections to the same angle, which is inaccurate when blocks have different cutting angles)
 - b. Avoid checking *Using Cutting Index* (this suggests your sections numbers (_s###) correspond to serial section numbers spaced equally apart)
 - c. Optional to allow DeepSlice to use your data to improve the neural network and its predictive accuracy.
- 4 Select **Submit**.
- 5 After all sections/slices are processed, press Download XML.



- 6 Insert the downloaded XML (titled by default 'results.XML') into your QVN folder with all your images for registration.