

May 24, 2024

🌐 Graphing area from an image series using ImageJ: a simple method

DOI

<https://dx.doi.org/10.17504/protocols.io.n2bvj326blk5/v1>


Vanessa Ho¹, Sally P Leys¹

¹University of Alberta



Sally P Leys

University of Alberta

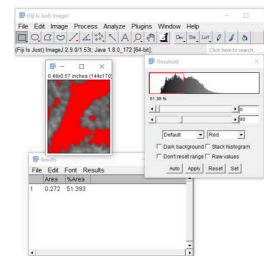

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.n2bvj326blk5/v1>

Protocol Citation: Vanessa Ho, Sally P Leys 2024. Graphing area from an image series using ImageJ: a simple method. protocols.io <https://dx.doi.org/10.17504/protocols.io.n2bvj326blk5/v1>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited



Protocol status: Working

We use this protocol and it's working

Created: October 05, 2023

Last Modified: May 24, 2024

Protocol Integer ID: 88815

Keywords: image analysis, area measurement, behaviour analysis, area from an image series, using imagej, graphing area, used imagej, pixel contrast, image series, image, change in area, good contrast of light, area, dark in the region

Funders Acknowledgements:

NSERC Discovery Grant

Grant ID: Leys 2022-03414

Abstract

This protocol provides a quick way to quantify change in area of an object in a large series of images, by pixel contrast. It's equivalent to rootpainter, except much simpler, but it does rely on good contrast of light and dark in the region you wish to measure.

Others have used ImageJ for this too, but the method and workflow we provide has been tested successfully by many students and therefore offers a more consistent, standardized process.

Image Attribution

Vanessa Ho

Materials

ImageJ Fiji <https://imagej.net/software/fiji/>

Before start

Make sure you have ImageJ downloaded, available here: <https://imagej.net/software/fiji/>

1 **Create a workflow (macro) in ImageJ** with the following steps: This will be run on all your images.

1. Open ImageJ
2. Opening an image,
3. Selecting Plugin>Macro>Record, which opens up a 'recorder' window.

Note

Once you have the recorder window open it will automatically record your actions from this point onward

2 **Recording the macro** for image processing:


1. Set image type to 8-bit grayscale
2. Select the rectangular selection tool and draw a rectangle over the image to create a ROI,
3. Crop the image to that region >Image>Crop,
4. Threshold the image >Image>Adjust>Threshold>Auto>Set>OK,
5. Measure the area fraction, >Analyze>Set Measurements>, check the Area Fraction box (no other measure is needed), >Analyze>Measure. This action sends the measurement to a results window.
6. Save the results to a folder by choosing >File>Save as>[choose your folder and type the file name].
7. Stop the macro recording in the 'recorder' window by clicking >Create.

A Macro.ijm window opens with the macro, click >File>Save as>[choose a folder and file name for the macro].

Close the image and result window now.

Note

The default is to measure the dark areas of the image (coloured in red by ImageJ). To measure the light areas instead check the "dark background" box on the threshold window.

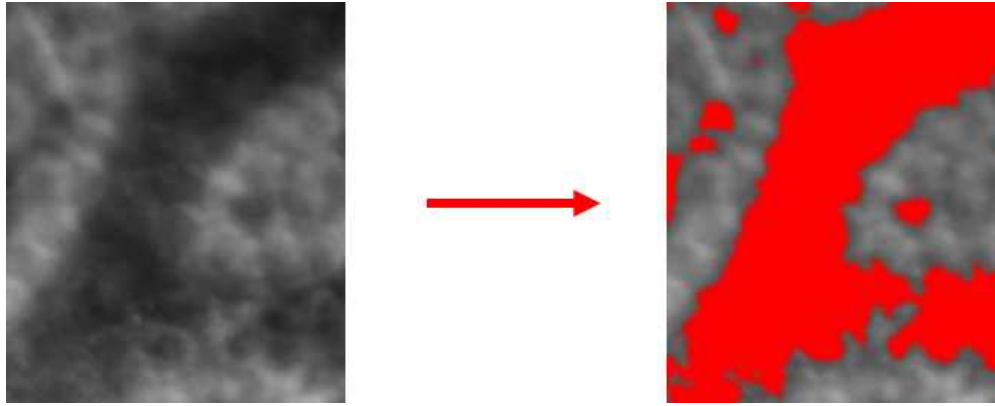


Example of ImageJ thresholding based on pixel brightness/contrast on a selected cropped region

3 **Run the macro** on a set of images:

1. In the main ImageJ toolbar choose >Process>Batch>Macro, to open the batch process window.
2. Select the folder which contains the image set to be analysed, and the folder where the data and the processed images will be placed.
3. Copy and paste the macro.ijm file from above into the window, and click >Process.

The image set will be analysed and the area fraction for each image will be calculated and placed in

a .csv file, which can be opened for graphing or further analysis in Excel or other programs

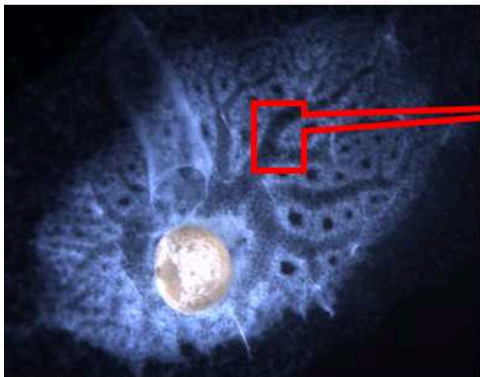

Ephydatia muelleri.