




Apr 30, 2023

Version 1

🌐 Hairy root generation in common bean (*Phaseolus vulgaris* L.) and selection of *Agrobacterium rhizogenes* clones V.1

DOI

<https://dx.doi.org/10.17504/protocols.io.261ge3bpjl47/v1>

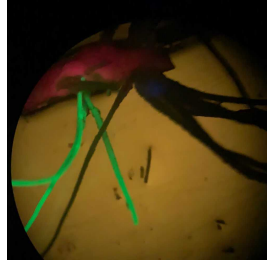

Ronal Pacheco¹, Georgina Estrada-Navarrete¹, Noreide Nava¹, Jorge Solis-Miranda¹, Carmen Quinto¹

¹Instituto de Biotecnología, Universidad Nacional Autónoma de México

Carmen Quinto: Corresponding author.

Ronal's protocols



Ronal Pacheco

Biotechnology Institute, National Autonomous University of M...

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.261ge3bpjl47/v1>



Protocol Citation: Ronal Pacheco, Georgina Estrada-Navarrete, Noreide Nava, Jorge Solis-Miranda, Carmen Quinto 2023. Hairy root generation in common bean (*Phaseolus vulgaris* L.) and selection of *Agrobacterium rhizogenes* clones. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.261ge3bpjl47/v1>

License: This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: April 06, 2023

Last Modified: April 30, 2023

Protocol Integer ID: 80078

Keywords: hairy root generation in common bean, generation of transgenic hairy root, generation of hairy root, hairy root generation, hairy roots in common bean, transgenic hairy root, hairy root, selection of agrobacterium rhizogene, agrobacterium rhizogene, mycorrhizal, legume, rhizobial symbiosis, genetic approach, applying genetic approach, protocol for the generation, common bean

Funders Acknowledgements:

Olivia Santana

Abstract

The common bean (*Phaseolus vulgaris* L.) is one of the legumes used to study the molecular mechanisms that regulate mycorrhizal and rhizobial symbioses. To study these mechanisms, the generation of transgenic hairy roots is a very advantageous method for applying genetic approaches. However, the generation of hairy roots is a difficult task that requires a lot of skill and experience. Here we show an optimized protocol for the generation of hairy roots in common bean.

Citation

Estrada-Navarrete G, Alvarado-Affantranger X, Olivares JE, Guillén G, Díaz-Camino C, Campos F, Quinto C, Gresshoff PM, Sanchez F (2006). Fast, efficient and reproducible genetic transformation of *Phaseolus* spp. by *Agrobacterium rhizogenes*.

<https://doi.org/>

LINK



Protocol materials

✕ Plastic Petri dishes (100×15 mm)

✕ 1.5 mL Eppendorf tubes

✕ 15 mL Falcon tubes

✕ PYREX Glass Rimless Test Tube 12×75mm **The Science Company®** Catalog #NC-0993



Seeds disinfection (when necessary)

1 Immerse the common bean seeds in [IM] 96 % volume ethanol for ⌚ 00:05:00 and wash them three times with sterile water.

5m

2 Immerse the seeds in [IM] 2 % volume sodium hypochlorite for ⌚ 00:05:00 and wash them three times with sterile water.

5m

Note

The concentration of sodium hypochlorite may be higher than indicated, but this depends on the quality of the seeds. When the quality of the seed is not very good, a higher concentration of sodium hypochlorite can damage a large number of them, rendering them useless for germination.

3 Keep the seeds at 🌡 4 °C in a disinfected container, e.g., sterile [X] Plastic Petri dishes (100×15 mm) , or reused Petri dishes previously sterilized with [IM] 96 % volume ethanol.

Seeds germination (1st day)

4 The disinfected seeds are placed, using sterile forceps, in a metal tray on a wet paper towel, previously sterilized in an autoclave, leaving 2 cm between the seeds.



Note

Paper towels must be moistened with deionized and pre-autoclaved water.

- 5 Cover the metal tray with aluminum foil and incubate it in a growth chamber at  28 °C for  46:00:00 to  48:00:00 in the dark.

3d 22h

Note

Position the tray at a slight downward angle to improve seed germination. The hilum should face downward.



Preparation of the inoculum of *A. rhizogenes* K599 (2nd day)

- 6 Spread 150 μL - 200 μL of the inoculum in Petri dishes containing solid LB medium with the appropriate selection antibiotic.

Note

The inoculum consists of a liquid culture of *A. rhizogenes*, transformed with the corresponding vector, and 80% 80 % volume glycerol 50 % (v/v) , stored at -80 °C . It is not recommendable to reuse the inoculum.

- 7 Incubate the Petri dishes inoculated in the previous step, for approximately 30:00:00 at 30 °C .

1d 6h

Note

If after 30 h the *A. rhizogenes* culture has not grown successfully, i.e., the culture layer is dry and very thin, do not use this inoculum for plant transformation and do the following:

- 7.1 Scratch off this thin layer of dried culture with a sterilized yellow tip or something similar.
- 7.2 Transfer this culture to an  1.5 mL Eppendorf tubes and add LB liquid medium.
- 7.3 Centrifuge the 1.5 ml Eppendorf tube  8000 rpm, Room temperature, 00:01:00 and then homogenize the content using a micropipette. The content must be viscous, but liquid enough to be pipetted.

1m



- 7.4 Prepare the inoculum in Eppendorf tubes (preferably 0.6 ml) by mixing an equal volume  50 % (v/v) of the liquid content previously obtained, and  80 % volume glycerol. Mix tubes by inversion and immediately place them in liquid nitrogen; finally, store the inoculum at  -80 °C .

Note

To use this inoculum, spread  150 μL -  200 μL of the inoculum along with an equal volume of sterile LB liquid medium in Petri dishes containing solid LB medium with the appropriate selection antibiotic. Incubate for approximately  30:00:00 at

 30 $^{\circ}\text{C}$.

Seedling transformation by *A. rhizogenes* K599 (3rd day)

- 8 Carefully puncture the hypocotyl area of the seeds several times, using a sterile needle tip (0.4 mm).





Note

Make sure that the hypocotyl tissue is not severely damaged. Damaged hypocotyl reduces seedling development and hairy root generation. The punctured area should not be too close or too far from the cotyledons.

- 9 Apply the inoculum of *A. rhizogenes* on the wounded zone, taken directly from the plates, using an autoclaved micropipette tip

Generation of hairy roots

- 10 Place the infected seedlings on the top of plastic tubes, i.e.,  15 mL Falcon tubes containing B & D medium. <https://link.springer.com/book/10.1007/1-4020-3735-X>

Stock solutions	Solution component	2000x stock concentration	Concentration in final solution
A	CaCl ₂	2.0 M	1.0 mM
B	KH ₂ PO ₄	1.0 M	0.5 mM
C	Fe-citrate ^b	0.02 M	10 µM
D	MgSO ₄	0.5 M	0.25 mM
E	K ₂ SO ₄	0.5 M	0.25 mM
F	MnSO ₄	2 mM	1.0 µM
G	H ₃ BO ₃	4 mM	2.0 µM
H	ZnSO ₄	1 mM	0.5 µM
I	CuSO ₄	0.4 mM	0.2 µM
J	CoSO ₄	0.2 mM	0.1 µM
K	Na ₂ MoO ₄	0.2 mM	0.1 µM

Table 1. Broughton and Dilworth (1971) nutrient solution. ^aAdjust the pH to 6.8 with KOH.

^bDissolve Fe-citrate with heating and keep in dark bottle

- 11 Place the plastic tubes inside glass tubes



PYREX Glass Rimless Test Tube 12×75mm **The Science Company®** Catalog #NC-0993

containing autoclaved deionized water and cover the glass tubes with plastic caps to prevent water evaporation.

- 12 Place glass tubes on racks and incubate in a growth chamber at 28 °C ,



16:00:00 light/ 08:00:00 dark until hairy roots emerge, 10-12 days post-infection (dpi).

1d

Note

When the first pair of leaves reaches the plastic caps (3-5 dpi), remove the caps, and seal the tube hole with parafilm or adhesive plastic (image). During this period, make sure that the level of water and B & D medium contained within the glass tubes and plastics tubes, respectively, is adequate.

- 13 Once the hairy roots have emerged, remove the primary root by cutting the stem 1-2 cm below the hairy roots.
- 14 Transfer the seedlings to autoclaved glass tubes containing sterile B & D medium and seal the tube hole with parafilm or adhesive plastic.



Note

Make sure the level of the B & D medium is below the hairy root calluses, as covering hairy root calli with B & D medium may retard their growth.





- 15 Incubate the seedlings for approximately three days under the same conditions described above to increase the biomass of hairy roots.

Selection of transformed hairy roots

- 16 Observe the fully developed hairy roots (15 to 16 dpi) using an epifluorescence microscope to remove non-fluorescent roots.

Note

Hairy roots must carry a plasmid containing a fluorescent reporter gene e.g., GFP, RFP, or YFP. Commonly, plasmids for RNAi-based gene silencing or gene overexpression carry a fluorescent reporter gene.

Selection of *A. rhizogenes* clones

- 17 Collect the fluorescent hairy roots carrying the RNAi-silencing or the overexpression vector, and the control vector, at the selected sampling time.

Note

If the amount of fluorescent hairy roots collected from a single plant is not enough tissue for RNA extraction, hairy roots from more than one plant will be needed.

- 18 Extract total RNA from hairy roots using an appropriate protocol and perform cDNA synthesis. We recommend using the following protocol

dx.doi.org/10.17504/protocols.io.8epv5jq24l1b/v1

- 19 Quantify transcript levels of the gene of interest by qPCR using two reference genes.

Note

We recommend as reference genes those encoding elongator factor 1 α (*EF1 α* , Phvul.004G075100.1) and β -tubulin (Phvul.009G017300.1). For RNAi silencing-based studies, we recommend choosing *A. rhizogenes* clones with a silencing efficiency of at least 70%.

