

Jun 16, 2019

Manual dissection of the *Schistosoma mansoni* esophagus and back end for proteomic analysis

In 1 collection

DOI

<https://dx.doi.org/10.17504/protocols.io.tq2emye>


Leandro Neves¹, R. Alan Wilson², William de Castro Borges³

¹University of Liverpool; ²University of York; ³Universidade Federal de Ouro Preto



Leandro Xavier Neves

University of Liverpool

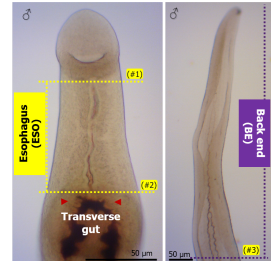

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.tq2emye>

Protocol Citation: Leandro Neves, R. Alan Wilson, William de Castro Borges 2019. Manual dissection of the *Schistosoma mansoni* esophagus and back end for proteomic analysis. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.tq2emye>

**Manuscript citation:**

Neves LX, Wilson RA, Brownridge P, Harman VM, Holman SW, Beynon RJ, Eysers CE, DeMarco R, Castro-Borges W. Quantitative Proteomics of Enriched Esophageal and Gut Tissues from the Human Blood Fluke *Schistosoma mansoni* Pinpoints Secreted Proteins for Vaccine Development. J Proteome Res. 2020 Jan 3;19(1):314-326. doi: 10.1021/acs.jproteome.9b00531. Epub 2019 Dec 2. PMID: 31729880.

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: September 19, 2018

Last Modified: June 16, 2019

Protocol Integer ID: 15866

Keywords: esophageal gland. schistosome esophagus. schistosome gut. gastrodermis. *Schistosoma mansoni*. quantitative proteomics., other schistosome organ, molecular characterisation of other schistosome organ, proteomic analysis schistosome, back end for proteomic analysis schistosome, manual dissection of the schistosoma, schistosome, esophageal gland cell, minor fraction of the whole parasite body, *Schistosoma mansoni*, whole parasite body, schistosoma, esophageal gland, esophagus, posterior esophageal gland, esophagus tube, dissection technique on worm, worms death by starvation, lower parts of the intestine, intestine, parasite, ovary of adult worm, posterior parts of the oesophageal gland, intravenous parasite, esophageal gland constituent, secreted gut protein, oesophageal gland, molecular characterisation of the esophageal gland constituent, secretory cell mass, disease through worm elimination, esophageal function, further digestion nutrient uptake, gut protein, gland secretions as the primary target, worms death

Abstract

Schistosomes are intravenous parasites with ability to survive in the mammalian host for decades, using its blood as a source of nutrients. The feeding process is multistep and takes place along the worm's alimentary tract, which comprises an (i) oral cavity opening to a short (ii) esophagus that is connected to the (iii) the gut caecum. The ultrastructural morphology of *Schistosoma mansoni* and *S. japonicum* has revealed the existence of two secretory cell masses surrounding the esophagus tube, referred to as the anterior and posterior esophageal glands (antESO and postESO, respectively). We recently established that the esophageal glands have a pivotal role in the first steps of blood processing. For instance, erythrocytes and leucocytes are quickly processed along the esophagus before they are propelled to the lower parts of the intestines for further digestion nutrient uptake. We propose that incorrect functioning of alimentary tract is associated with worms death by starvation. This was first observed in the self-cure response of Rhesus macaque (*Macaca mulatta*), one of few known vertebrate hosts capable of combating the disease through worm elimination once the infection is established. Classical immunoproteomics (2D-PAGE and Western blotting) has revealed potential targets in both exposed tegument and secreted gut proteins. Recently, a more detailed investigation using *S. japonicum* in the Rhesus model shed light on the possible operating mechanisms that prevent parasite feeding on blood. In addition, ultrastructural studies and immunocytochemistry on surviving worms indicated the esophageal lumen and the gland secretions as the primary targets of a potent and protective humoral immune response that ultimately disrupts the esophageal functions. Therefore, the molecular characterisation of the esophageal gland constituents is imperative if one intends to emulate the Rhesus self-cure response for therapeutic purposes. However, this is not a trivial task as challenges are multiple. Perhaps, the most important caveat is that both anterior and posterior parts of the oesophageal gland represent a minor fraction of the whole parasite body (or even of its head), meaning that dominant constituents, such as those derived from muscle tissues, suppress the identification of a unique set of gland products. Indeed, our recent investigation on the soluble protein composition of a *S. mansoni* preparation failed to detect the presence of gland products attesting for their low abundance in the whole worm. Although optimized protocols for chemical/enzymatic dissection are reported for isolation of testes and ovary of adult worms, no method has proved feasible for gastrodermal epithelium and the esophageal gland cells. We tackled these challenges by developing a dissection technique on worms preserved in RNAlater solution aided by an essential set of scissors and tweezers that delivers adequate precision during the procedure. The method herein described is compatible with downstream proteomic analyses using mass spectrometry and has permitted a large-scale characterisation of proteins expressed in the *S. mansoni* esophagus and gastrodermis. This methodology can be applied in the molecular characterisation of other schistosome organs and tissues that present a well-defined anatomic location.

Guidelines


This protocol will guide you through the dissection of the esophagus and back end of adult schistosome male worms. Although the recovery of adult worms by portal perfusion of infected animals is obviously indispensable this is not addressed in this protocol. Nevertheless, worm handling and fixation are utterly important and should be performed according to the following instructions:

- Perfuse 45-days infected animals using RPMI-1640 medium (without phenol red) buffered with 10mM HEPES, pH 7.4, containing 4 UI/mL of heparin.
- Wash the parasites in the same medium (pre-warmed at 37°C) until total removal of blood and tissue debris.
OBS.: always use pre-warmed medium to prevent getting worms too curled up as this makes the dissection procedure difficult.
- Carefully transfer the parasites, by pouring them with plenty of pre-warmed medium, into a Corning 50 mL conical tube.
- Discard the medium and add 5-6 volumes of RNAlater for instantaneous worm fixation.
- Keep the worms at 4°C until dissection.

- **CRITICAL POINT** *Parasites respond to the ex vivo environment; thus it is mandatory to avoid a long waiting between the perfusion and fixation in RNAlater. The total procedure time (perfusion, washings and fixation) must not exceed 10 min.*

Procedure overview

Once the parasites are fixed you are able to perform the dissection. In summary, the isolation of the schistosome esophagus of male worms is achieved by the removal of the oral sucker, followed by the detachment of the esophageal region along the line of the transverse gut. The posterior third of the parasite (BE) is also dissected and constitute a second sample for comparative analysis

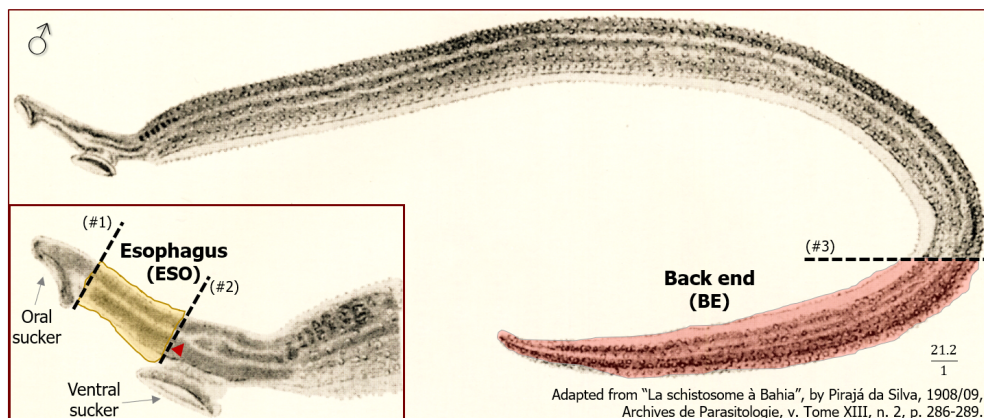

Cut #1 detaches the oral sucker
Cut #2 releases the esophagus
Cut #3 releases the back end

► indicates the bifurcated transverse gut. The presence of hemozoin in this region (a dark pigment) facilitates the visualization of a "limit zone", where incisions can be made with little or no contamination of the ESO with gastrodermis and gut content.



- **CRITICAL POINT** *The procedure must be entirely performed in ice-cold conditions. Keep ESO and BE fragments on ice during the dissection, then store them in the fridge.*

ALTERNATIVE METHOD APPLICATION: Once the set of stereomicroscope, scissors and tweezers are available, different incision points can be used for obtaining alternative body sections. It is important that the target tissue/organ is anatomically well-defined and consider that precision decreases proportionally to the fragment size. In our lab we have successfully dissected the whole head, esophagus, back end, ovary and ootype on either *S. mansoni* or *S. japonicum* male and female worms. Check the protocols available at our collection.

Materials

MATERIALS

 RNA later Merck MilliporeSigma (Sigma-Aldrich) Catalog #R0901-100ML

Safety warnings

- ⚠ Adult schistosomes are not an infective stage thus biological hazard is not evident. Nevertheless, follow the safety instructions of all chemicals used (e.g. RNA later, heparin). Always wear gloves and lab coat.



Before start

Prepare your workplace

- Fill a Styrofoam box with crushed ice. This is where you keep the worms and dissected fragments.
- Wear powder-free nitrile gloves.
- Add approximately 100 microlitres of RNAlater to a 1.5 mL microtube and keep it in the ice. Label this tube 'ESO'.
- Add approximately 500 microlitres of RNAlater to a second 1.5 mL microtube and keep it in the ice. Label this tube 'BE'.
- Add approximately 1 mL of RNAlater to a 2 mL microtube and keep it in the ice. Label this tube 'Dissected body'.
- Fill a 50 mL beaker with deionized water. You may have to immerse the scissor blades and the tweezers into the water periodically to get rid of salt crystals formed during the procedure. Make sure you have easy access to soft and lint-free paper to dry out the tools.
- Fill to half depth a small glass Petri dish (~4.5 cm of diameter) with RNAlater. Let it chill on the ice a few minutes.
- Turn on the stereomicroscope and switch on the bottom light only; this usually gives the best contrast during the dissection.
- Get notepad and pen to make notes about the number of ESO and BE fragments you have dissected.
- Get a laboratory timer.
- Ensure you have a P20 micropipette and disposable tips with you during the procedure.
OBS.: eventually you need to cut the tip to open its orifice. Use ordinary scissors to do that.
- Put few worm pairs in the Petri dish (10-15). Avoid a huge number of parasites in the dish as this complicates the procedure.
- Finally, take the scissors and tweezers and a comfortable chair. Make sure all the items above are reachable while seated.

Place the Petri dish with worms on the microscope

- 1 Take the Petri dish out of the ice bath, touch it against a paper to remove the excess of water, place it under the stereomicroscope and adjust the focus.

Equipment

new equipment

NAME

Stereo microscope EZ4, Leica

BRAND

10447197

SKU

Leica EZ4 educational stereomicroscope.

SPECIFICATIONS



Start a 10 min countdown

- 2 This will help you to control how long the Petri dish is out of the ice.

Zoom in at the head of a male worm

- 3 Look for a male parasite with uncurled head/esophagus. Increase the lenses magnification and use your dissecting tools (tweezers and scissors) to keep the head in the field-of-view.

Equipment

new equipment

NAME

Biology Tweezer 115 mm, Ideal Tak

BRAND

7SG.CX.0

SKU

Biology Tweezer 115 mm (Ideal Tak Chiasso, Switzerland).
Manufacturer Part Number: 7SG.CX.0

SPECIFICATIONS



Equipment

new equipment

NAME

Vannas scissors, John Weiss

BRAND

0103123

SKU


Vannas scissors (John Weiss, Milton Keynes, UK). Straight with sharp tips.

SPECIFICATIONS




Use the scissors to make the first incision at the posterior of the oral sucker

- 4 Hold the parasite carefully with the tweezers, positioning it in an angle that allows you to clearly see the suckers (*i.e.* lateral view, or suckers facing up), then make the first incision (#1) at the posterior of the oral sucker. This releases the oral sucker; this fragment is composed mostly by muscle tissue and usually has limited utility. If you don't need the oral sucker just leave it floating around.

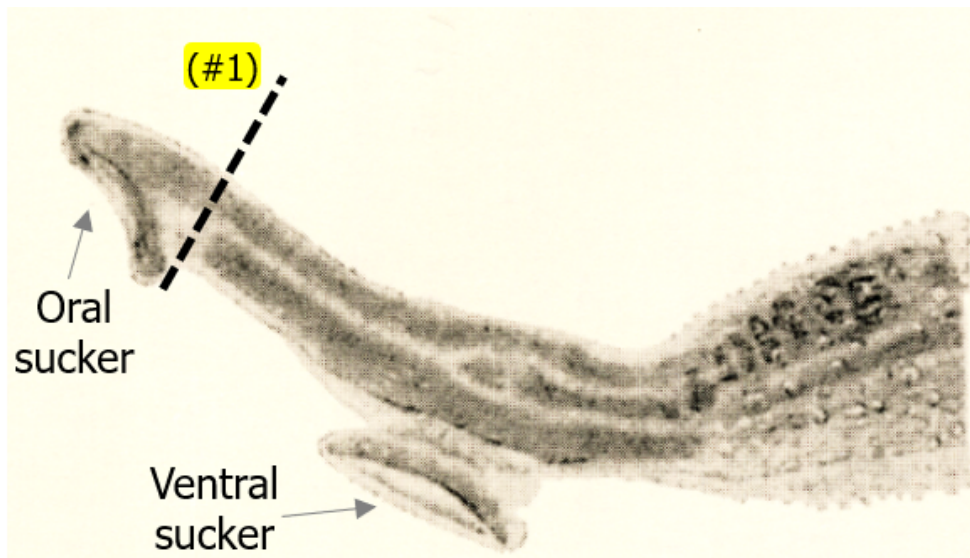

Cut off the esophagus (ESO) along the transverse gut

- 5 The second incision (#2) must be performed along the transverse gut (*i.e.* where the dark pigmentation of hemozoin starts; indicated by the red arrow in the illustration below). Be sure you are not simultaneously cutting the ventral sucker off.

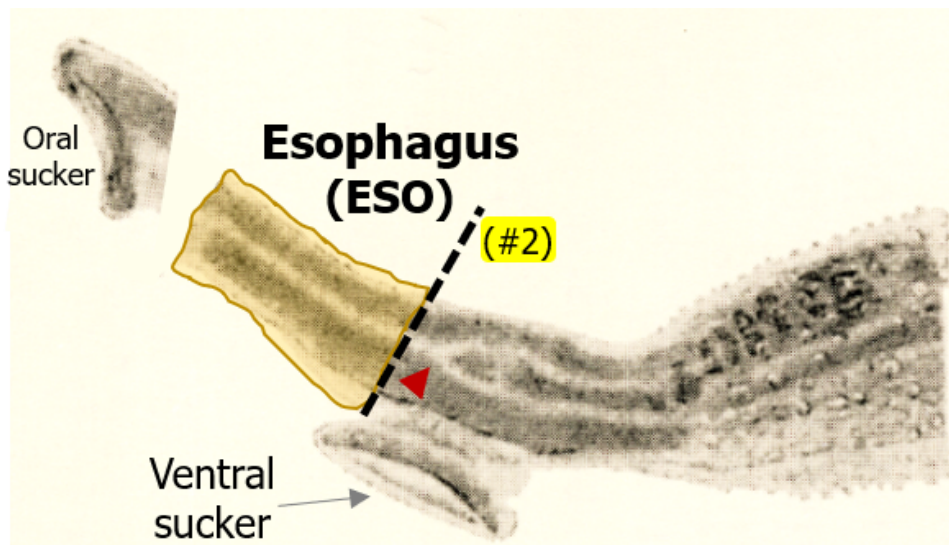

Collect the ESO fragment using a P20 micropipette

- 6 Now, the ESO fragment should be detached and floating around. Aspirate it with the RNAlater solution using the P20 micropipette then dispense it inside the ESO microtube. Always keep the tube on ice.

CRITICAL STEP *do not aspirate any other fragment or particulate material with the ESO.*

Record your progress in the notepad


7

ESO



Cut off the back end

- 8 The final cut (#3) is considerably easier and usually does not require high magnification. Simply cut the posterior third of the parasite body. Collect the BE fragment with the tweezers, transfer it to the BE microtube and keep in the ice box.

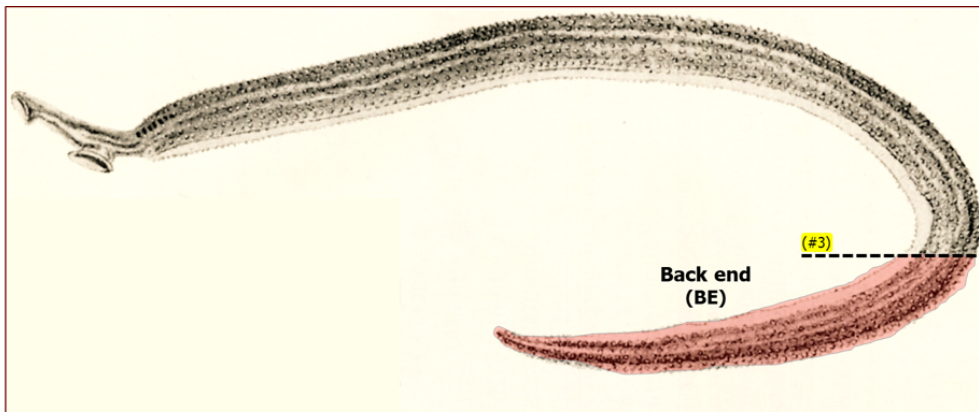

Record your progress in the notepad

9

BE



Find another male parasite and repeat steps 3 to 9

- 10 If needed, you can immerse the tweezers and scissors in the beaker with deionized water to get rid of salt crystals. Dry your tools in a soft and lint-free paper before continuing the dissection.

When the 10 minutes timer runs out place the Petri dish back in the ice

- 11 In Step 2 you started a 10 minutes countdown. You can dissect as many worms as you can in this time, but it is important to chill the Petri dish every 10 minutes to preserve your material.



ATTENTION *This is an opportunity to transfer the dissected worms to the "dissected body" microtube then add new ones to the Petri dish for the next round. In addition, if you notice that salt crystals are accumulating in the RNAlater solution discard it, wash the Petri dish with deionized water, dry it with lint-free paper and add fresh RNAlater solution.*

PAUSE POINT *If you wish, store the dissected fragments in the fridge and continue the procedure another time.*

Start a new round of dissection

- 12 Restart from Step 1.