



Oct 11, 2022

Post-Surgical Dissection of Ovaries

DOI

<https://dx.doi.org/10.17504/protocols.io.j8nlkkzn1l5r/v1>

Stephen Fisher¹, Marielena Grijalva¹, Rong Guo¹, sarahjoh¹, Hieu Nguyen¹, John Renz², Jean G Rosario¹, Steven Rudich², Brian Gregory¹, Junhyong Kim¹, Kate O'Neill¹

¹University of Pennsylvania; ²Gift of Life Donor Program



Stephen Fisher

University of Pennsylvania

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.j8nlkkzn1l5r/v1>

Protocol Citation: Stephen Fisher, Marielena Grijalva, Rong Guo, sarahjoh, Hieu Nguyen, John Renz, Jean G Rosario, Steven Rudich, Brian Gregory, Junhyong Kim, Kate O'Neill 2022. Post-Surgical Dissection of Ovaries. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.j8nlkkzn1l5r/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: August 05, 2022

Last Modified: October 11, 2022

Protocol Integer ID: 68273

Keywords: protocol details dissection of the ovary, aspect of the ovary, ovary, superior to the ovary, attachment to the mesosalpinx, infundibulopelvic ligament, uteroovarian ligament, dissection, preparation for 10x visium, mesosalpinx, marking dye, 10x visium, 10x multiomic, orientation of the specimen, pathology review

Funders Acknowledgements:

NIH

Grant ID: U54HD104392

Abstract

This protocol details dissection of the ovaries in preparation for 10X Visium, 10X Multiomics, pathology review, and biobanking. Care must be taken to systematically divide the ovary and to maintain orientation of the specimens as they are removed using the three points of fixation for the ovary (the uteroovarian ligament, the infundibulopelvic ligament and the mesosalpinx). Once the ovary is isolated it is difficult to discern which aspect of the ovary is superior to the ovary's attachment to the mesosalpinx and which is inferior. Colored marking dyes can be used to help maintain orientation provided they do not interfere with downstream analyses.

Materials

- Ice
- Clean cutting surface
- Disposable scalpels
- Kimwipes
- 3D Mold (pre-chill) (see Figure 2)
- Disposable rulers (with centimeter measurements)
- Ice bucket
- Biohazard bags and container
-  Nuclease-free Water
-  DMEM/F-12 **Thermo Fisher Catalog #11320033**
-  RPMI 1640 Medium **Thermo Fisher Catalog #11875085**
-  MACS® Tissue Storage Solution **Miltenyi Biotec Catalog #130-100-008**
-  DPBS (10X), no calcium, no magnesium **Thermo Fisher Scientific Catalog #14200075** - dilute to 1X with nuclease-free water
-  Marking dye **Cancer Diagnostics Catalog #03000P**

Equipment

VWR® Pour-Boat Weighing Dishes

NAME

flat-bottomed dish

TYPE

VWR International

BRAND

10803-166

SKU

Before start

Dilute phosphate buffered saline (Gibco; 14200-075) to 1X with nuclease-free water.



- 1 Receive excised female reproductive system on ice en bloc. Orient the specimen with the anterior aspect of the uterus facing up based on presence of vesicouterine reflection.
- 2 Using a disposable scalpel, incise the mesovarium. Orient ovary with the uteroovarian ligament on the right, infundibulopelvic ligament on the left, and mesovarium on top for the right ovary (Figures 6).

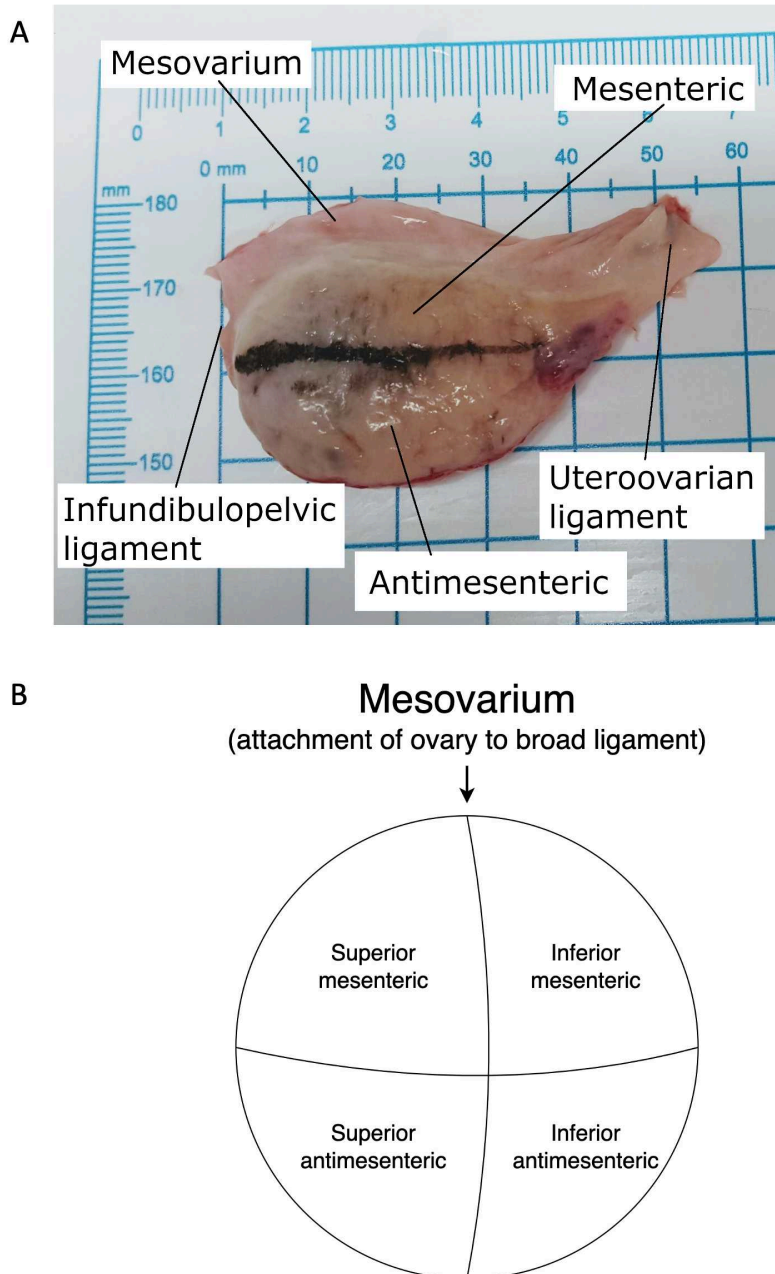


Figure 1. Use of fixed anatomic structures in the ovary. There is laxity in the structures supporting the position of the ovary which results in both intra- and inter-individual variations in the organ's orientation within the pelvis. In order to standardize collection and maintain orientation of tissue sections, a model proposed at a NICHD-sponsored workshop that used fixed anatomic structures to create a regional map of the ovary, was applied. Here the black line was drawn on the inferior half of the ovary when it is divided by its attachment to the mesovarium. The line delineates the mesenteric and antimesenteric portions of the inferior aspect of the right ovary.

- 3 Dry the ovary thoroughly with a Kimwipe.
- 4 Use marking dye to maintain orientation of the ovary with the uteroovarian ligament on the right, infundibulopelvic ligament on the left, and mesovarium on top.
- 5 Take the weight and, using a disposable ruler, take spatial measurement of each ovary.
- 6 Return ovaries to ice.
- 7 Place the first ovary in the pre-chilled mold (Figure 2).

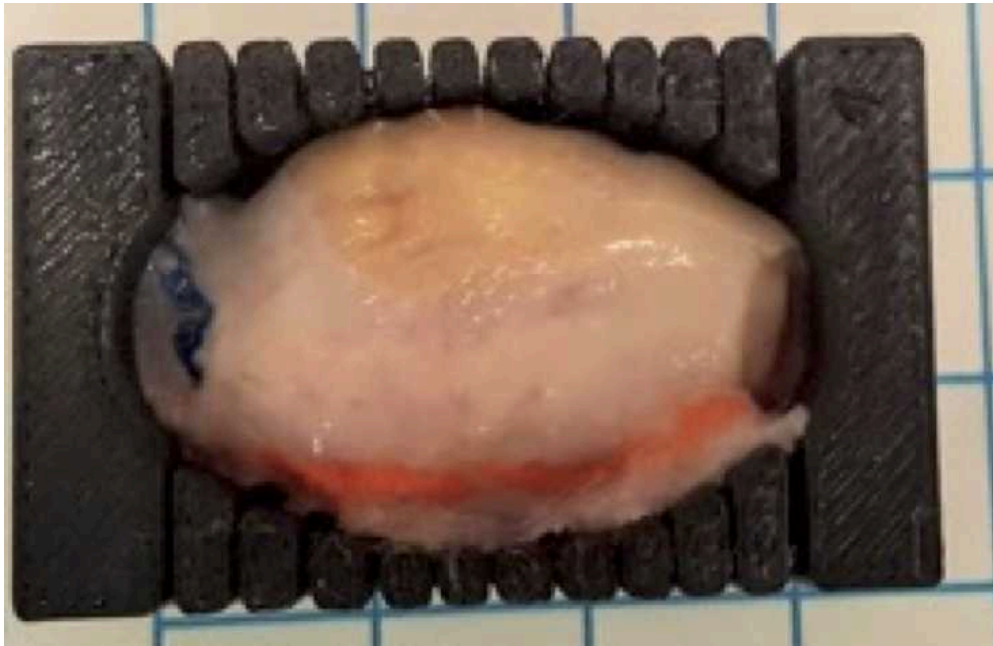


Figure 2. Pre-Chilled Ovary Mold. A 3D mold for long-axis sectioning of the ovary was designed in Autodesk Fusion 360 and printed in polylactic acid (PLA) using a Stratasys F170 Series 3D printer. A razor blade or dissecting knife is passed sequentially through each fenestration which allows for standardization in the width of ovarian slices obtained as well as maintains orientation of individual tissue pieces.

Dr. Alison Pouch, Dr. James Gee, Olivia Sandvold, and Kate da Silva were responsible for the design and execution of the 3D ovarian molds used in Basic Protocol 3. 3D molds were printed courtesy of the University of Pennsylvania Libraries' Biotech Commons.

- 8 Divide the ovary along the long axis of the ovary from the uteroovarian to the infundibulopelvic ligament using the guidelines of the mold (Figure 1a).
- 9 Take each slice and orient mesovarium ridge at 12 o'clock.
- 10 Divide each slice along the line of the mesovarium ridge to divide the ovary into portions superior to and inferior to the mesovarium.
- 11 Take each half and divide it again in half to separate the ovary into the mesenteric and antimesenteric portions (Figure 1b).
- 12 For later 3D anatomical positioning, the quadrant of the tissue piece and the slice number along the long axis should be tracked for each tissue specimen.
- 13 Tissue can be processed with protocols **OCT-Embedded Tissue Preparation**, **Tissue Fixation Preparation**, or **Snap-Frozen Tissue Preparation** depending on the desired downstream processing.
- 14 Repeat steps with the contralateral ovary, using a fresh, ice-cold ovary mold.