

Sep 05, 2024

Version 4

Publishing Generic Organ Scaffold as a SPARC Dataset V.4

DOI

<https://dx.doi.org/10.17504/protocols.io.6qpvr8q72lmk/v4>

Mabelle Lin¹, Hugh Sorby¹

¹University of Auckland



Mabelle Lin

University of Auckland

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.6qpvr8q72lmk/v4>

Protocol Citation: Mabelle Lin, Hugh Sorby 2024. Publishing Generic Organ Scaffold as a SPARC Dataset. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.6qpvr8q72lmk/v4> Version created by **Mabelle Lin**

License: This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: September 05, 2024



Last Modified: September 05, 2024

Protocol Integer ID: 106970

Keywords: Generic organ scaffolds, SPARC dataset, generic organ scaffold for publication, publishing generic organ scaffold, generic organ scaffold, scaffold mapping tool, sparc dataset, sparc, publication, reproducible manner, reusable to the community, mapping

Funders Acknowledgements:

NIH SPARC

Grant ID: OT3OD025347

Abstract

In this protocol, we will demonstrate how to prepare a generic organ scaffold for publication as a SPARC dataset using the Scaffold mapping tools. This protocol will ensure that a generic organ scaffold is made available and reusable to the community in a reproducible manner.

- 1 Use the latest **official** version of the **Scaffold mapping tools** to generate your organ scaffold. **Do not use a developer installation to generate a scaffold.**
- 2 Use the standard generic organ scaffold workflow to generate organ scaffold files. The workflow is installed by default with the mapping tools MAP Client installer and it is also available from **Github**. In this tutorial, as we are using the generic human stomach scaffold as an example, save a copy of the workflow under "sparc-generic-scaffold-stomach-human".

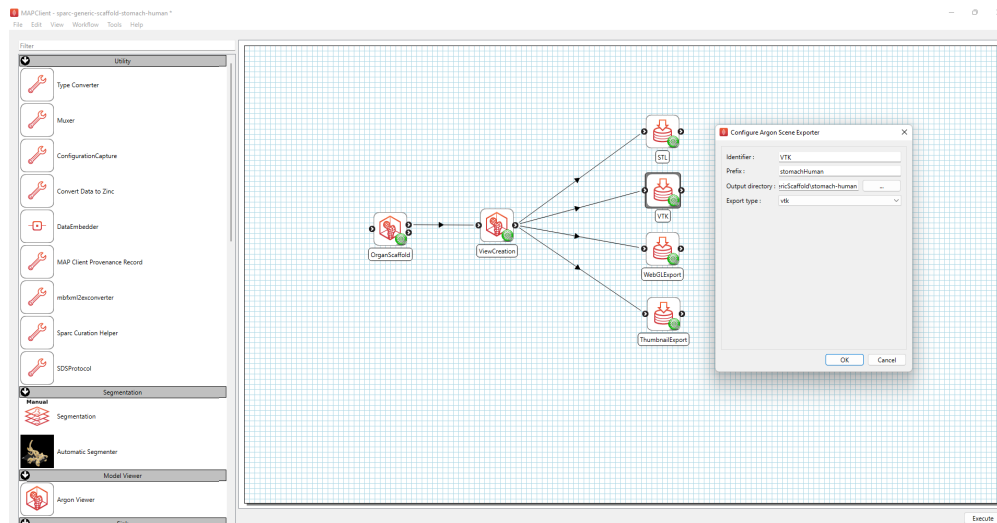


Figure 1. Standard generic organ scaffold workflow.

Ensure that the prefixes in the configuration settings for VTK and STL exporter are the same as the prefix used for the thumbnail exporter. Select the same output directory for each of the webGL, thumbnail, VTK, and STL exports.

After setting up the workflow as shown in Figure 1, save the workflow (Ctrl-s) and execute the workflow by pressing the **Execute** button (or Ctrl-x).

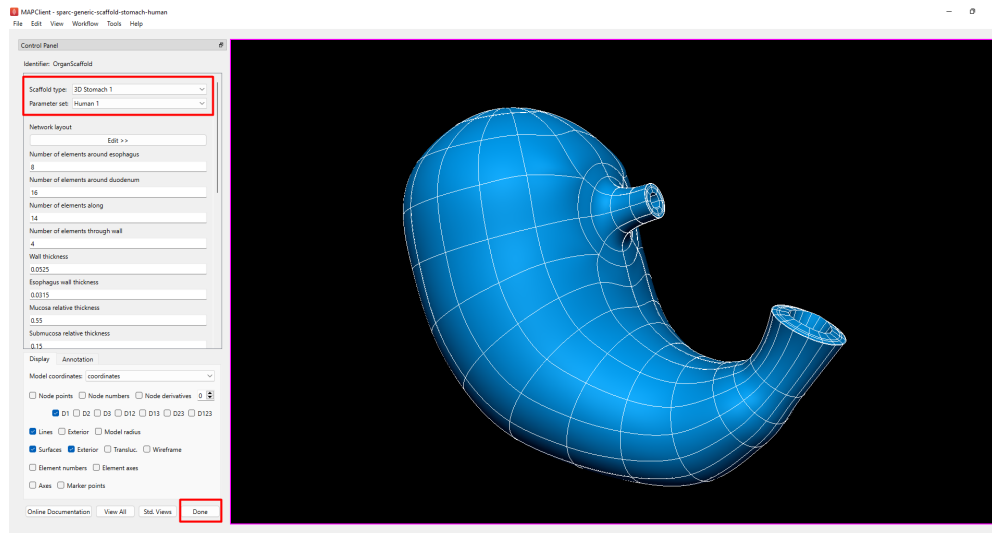


Figure 2. Scaffold creator interface for configuring organ scaffold and species.

In the scaffold creator interface (Figure 2), select the organ scaffold and species you wish to publish as a dataset. Click **Done**.

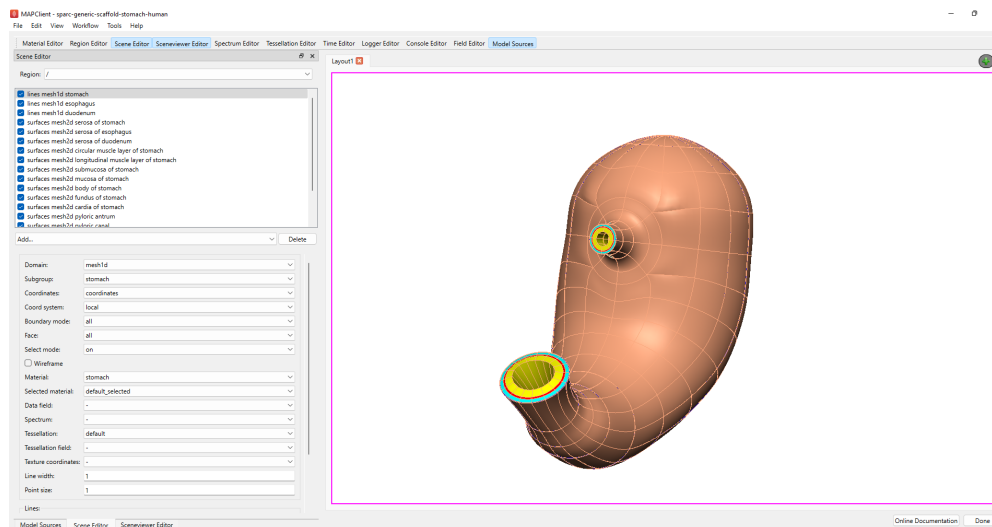


Figure 3. Argon Viewer for configuring graphics for the organ scaffold.

The next step involves configuring Argon Viewer to generate the graphics for the organ scaffold (Figure 3). Add relevant nodes, lines, and surfaces to showcase the annotated markers and groups in the WebGL export. Adjust the visualization viewpoint to an aesthetically pleasing angle and click **Done**. The workflow will create a snapshot of the visualization and use it as the thumbnail for the dataset, for this reason the viewpoint should show the best representation of the organ scaffold.

Upon clicking **Done**, the workflow will generate the WebGL, thumbnail, VTK, and STL files in the output directory that was set in **Step 2**.

- 3 Check that your WebGL visualisation works in ScaffoldVuer. ScaffoldVuer is what the SPARC portal uses to display your scaffolds in the web browser. You can drag and drop the output folder generated by scaffold mapping tool to **ScaffoldVuer** to check that the WebGL files are generating all scaffold features required for your scaffold on the portal. (You can also zip up your files and drop the zipped archive on ScaffoldVuer.)
- 4 Use the **official scaffold SDS creation workflow** to prepare a dataset for submission to **Pennsieve**. The workflow must be cloned from **Github**, *git clone -b v0.1.1 <https://github.com/mapclient-workflows/official-scaffold-sds-creation-workflow>*. The information on the workflow is recorded in the provenance information to aide reuse.

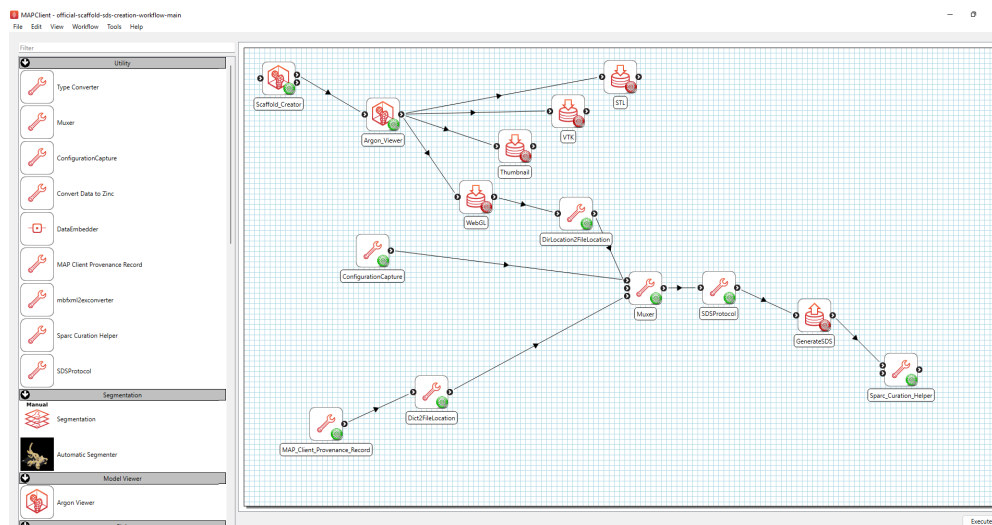


Figure 4. The official scaffold SDS creation workflow.

First, configure the GenerateSDS step (Figure 5). Enter the dataset name in the "Dataset Name" field, select Scaffold under "Dataset type" and set where you would like to save the generated SDS dataset in the "Output directory".

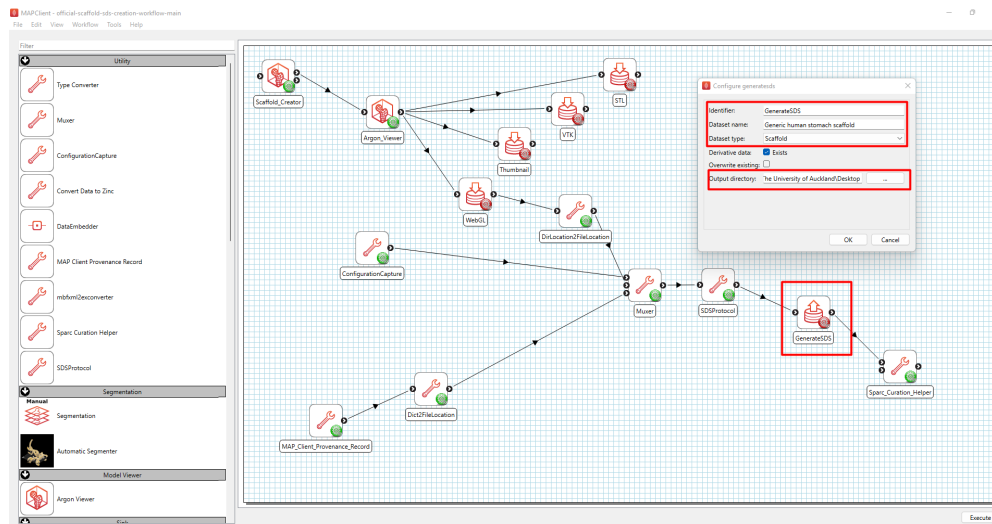


Figure 5. Configuring the GenerateSDS step.

Next, import configuration settings from the standard generic organ scaffold workflow by going to "Import Config" under the "File" tab (Figure 6).

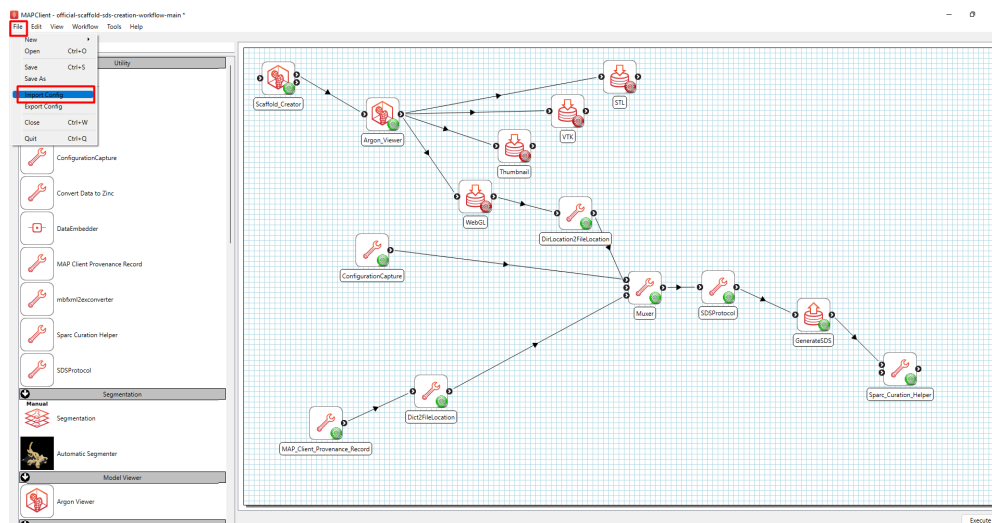


Figure 6. Import configuration settings.

Go to the directory where the standard generic organ scaffold workflow was saved in **Step 2**. Open the filter file type list at the bottom right corner of the dialog and select "MAP Client Project file (map-client-workflow.proj)" so that the workflow description file becomes available for selection (Figure 7).

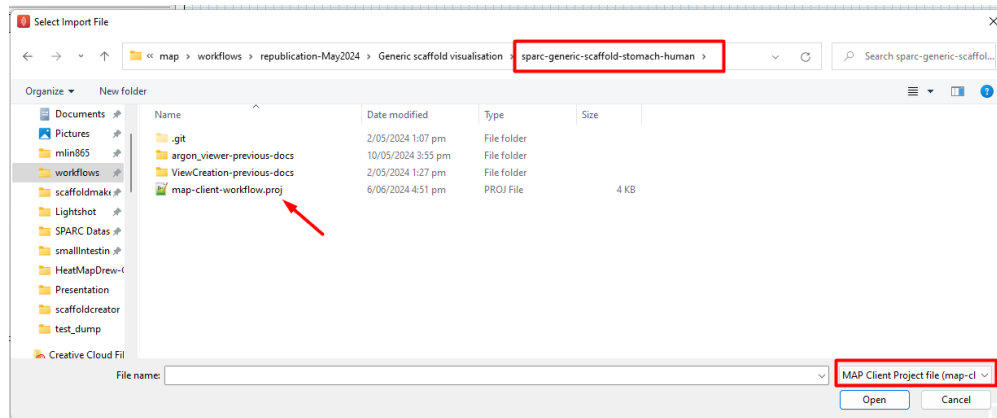


Figure 7. Selection of workflow file for import.

Next, an "Import Workflow Configuration" window will appear. Select the matching configurations as shown in Figure 8 and click **Import**, followed by **Close**. A prompt will appear indicating that the selected configurations have been successfully imported.

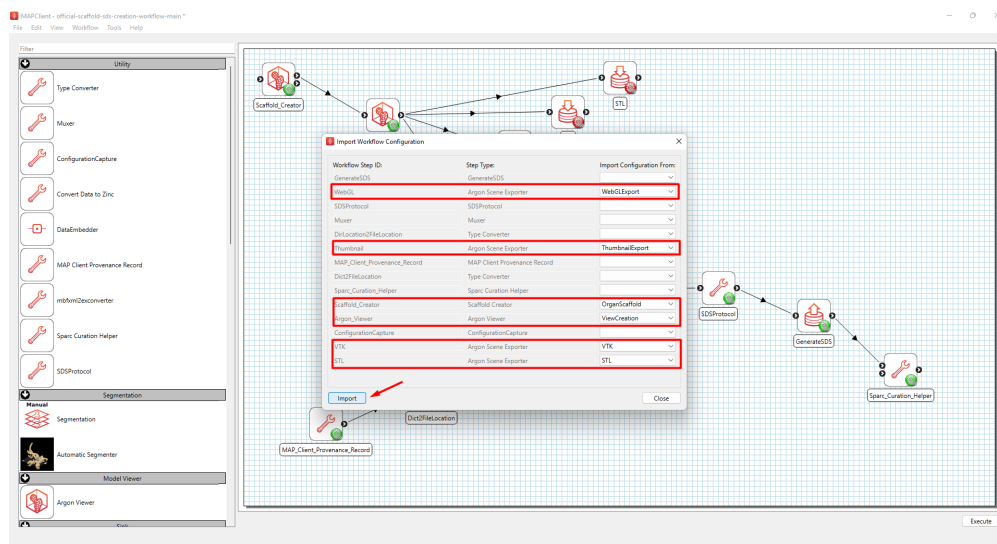


Figure 8. Import workflow configuration settings.

Currently, relative references to the paths or directories can be incorrect when importing configurations. This happens when you import a configuration to a workflow that has a different relative location to the files or directories from the source workflow. If this is the case, when you import configurations from a workflow that is properly configured the corresponding configurations for the importing steps may be/stay improperly configured (an improperly configured step shows a red gear icon). In this situation you will need to reset the relative path for the effected steps.

When following this tutorial, this can be done by clicking on the red gear icon and updating the output directory to point to the output directory set in the "sparc-generic-scaffold-stomach-human" workflow in **Step 2**. The red gear icon should then become

green (once the path to the output directory has been updated correctly). This process needs to be repeated for all four export steps. For convenience, the output directory can be copied and pasted from one step configuration to the next. Once all wheels are green, the workflow is ready to be saved (Ctrl+s) and executed (Ctrl+x).

- 5 You will see the organ scaffold interface. If the scaffold creator settings have been imported correctly, you should see the correct organ and species on the interface. Click **Done** to get to the Argon Viewer interface and check that the graphics are imported correctly and generated as per what were configured in **Step 2**. Click **Done**.
- 6 As shown in Figure 9, the next interface is a list of fields for information required for populating the dataset description, submission and manifest files for the dataset. Note that input history is stored for most of the fields and you can click on the downwards arrow to review the previously saved text and reuse them if appropriate. Otherwise, enter the required information.

MAPClient - official scaffold-ids creation-workflow-main

File Edit View Workflow Tools Help

Basic information

Title: Generic human stomach scaffold

Subtitle: arch scaffold available for registration of segmented neural anatomical-functional mapping of enteric neural circuits

Funding: C730002349

Acknowledgments: Keywords: myenteric plexus, human, stomach

Study information

Study purpose: to create an annotated generic human stomach scaffold for the registration of segmented data obtained by experimental groups

Study data collection: An anatomically-based 3D scaffold of the human stomach is created to map nerve ending pathways. The stomach scaffold is generated with a configurable central path defined from the fundus apex to the duodenum. It uses an annotated fundus, body, antrum, pylorus, and duodenum along the central path and use that to define the respective regions of the stomach. The cross base of the central path provide control of the major and minor radii of the stomach in different sections. The human stomach scaffold is parameterized with data to represent the anatomy as accurate as possible.

Study primary conclusion: None stated

Study organ system: Stomach

Study approach: Computational modelling

Study collection title: computational modelling, anatomical mapping

Study technique: computational modelling, anatomical mapping

Submission information

Consortium data standard: UNMC

Funding consortium: UNMC

Award number: C730002349

Milestone achieved: Documentation, tutorials, and training materials related to MAP functionality

Milestone completion date: August 2024

	Mon	Tue	Wed	Thu	Fri	Sat	Sun
31	29	30	31	1	2	3	4
32	5	6	7	8	9	10	11
33	12	13	14	15	16	17	18
34	19	20	21	22	23	24	25
35	26	27	28	29	30	31	1
36	2	3	4	5	6	7	8

Other information

Identifier description: Protocol for dataset

Relation type: HasProtocol

Identifier: https://doi.org/10.17504/protocols.io.6qpvr8q72lmk/v4

Identifier type: URL

Contributor information

Name: Lin, Mabelle

ORCID: https://orcid.org/0000-0002-4388-2181

Affiliation: Auckland Bioengineering Institute

Role: Corresponding Author Researcher

Figure 9. Information required for generating a dataset description file, submission file and manifest files for the dataset.

Note that for "Keywords" and "Study technique", do not add a list of terms separated by commas. Instead, each term needs to be either selected from the input history or added manually if the term is new. Click on "Add" to add the term. Perform the same steps to add more terms to the field.

Of particular interest to such organ scaffold datasets, under "Other information", make sure that there is an entry for the link to the PMR repository. The details to be entered are as follows:

Identifier description: *Link to repository*



Relation type: *IsVariantFormOf*

Identifier: <https://models.physiomeproject.org/e/8f0> [Should be link to the latest PMR exposure for the organ scaffold]

Identifier type: *URL*

Click **Add** to add a second entry. In this entry, we will include the link to Scaffold Mapping Tools and the details are as follows:

Identifier description: *Link to Scaffold Mapping Tools*

Relation type: *HasSoftware*

Identifier: <https://sparc.science/tools-and-resources/1mv8q3JckdpSYpPK9dvdKx>

Identifier type: *URL*

Click **Add** to add a third entry. In this entry, we will include the link to the protocol for scaffold publication and the details are as follows:

Identifier description: *Protocol for dataset*

Relation type: *HasProtocol*

Identifier: <https://doi.org/10.17504/protocols.io.6qpvr8q72lmk/v4>

Identifier type: *DOI*

Under "Contributor information", make sure that the ORCID is provided as a link to the contributor's ORCID identifier and one of the contributors is assigned as "CorrespondingAuthor" in the "Role" field.

Click **Done** on the bottom left hand corner when all fields have been filled.

- 7 The dataset curation helper tool will guide us through the process for annotating the manifest file for the derivative folder and add context information to the dataset.

Under the "Scaffold Annotation" tab (Figure 10), click on **Fix All Errors** and **Yes** when a confirmation window pops up asking you to "Let this magic tool fix all errors for you?". Perform this repeatedly until there is no more errors reported in the "Errors" box.

If you expand the directory shown under "Annotations", you should see that the VTK, STL and thumbnail files are at the bottom of the list. The thumbnail will also show up in the preview window if you click on the thumbnail path (Figure 10).

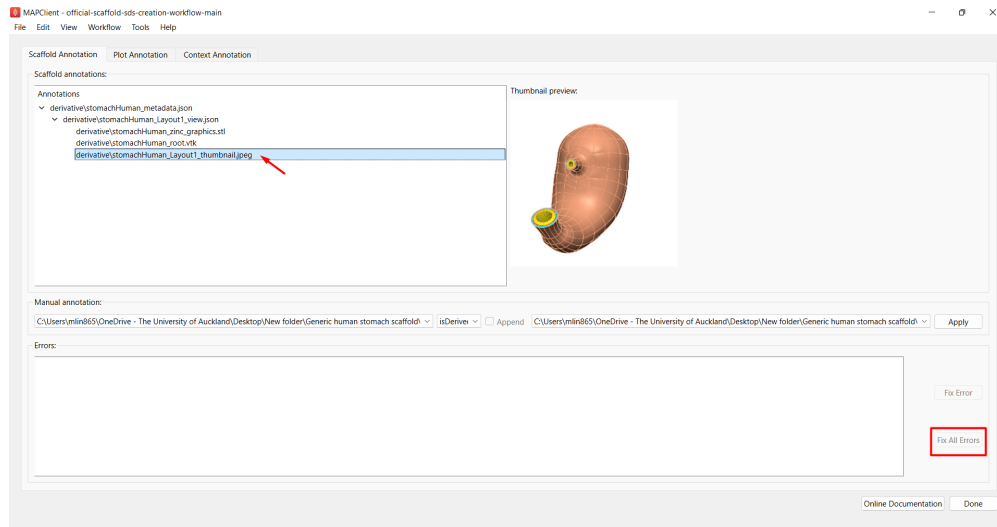


Figure 10. Dataset curation tool for annotating the manifest file for the derivative folder.

Next, click on the **Context Annotation** tab. Here, we are creating a context card which can be used as a sample template for all the other generic organ scaffolds.

First, under "Summary", enter information as shown in Figure 11 but tweak the contents to match the organ species you are working on.

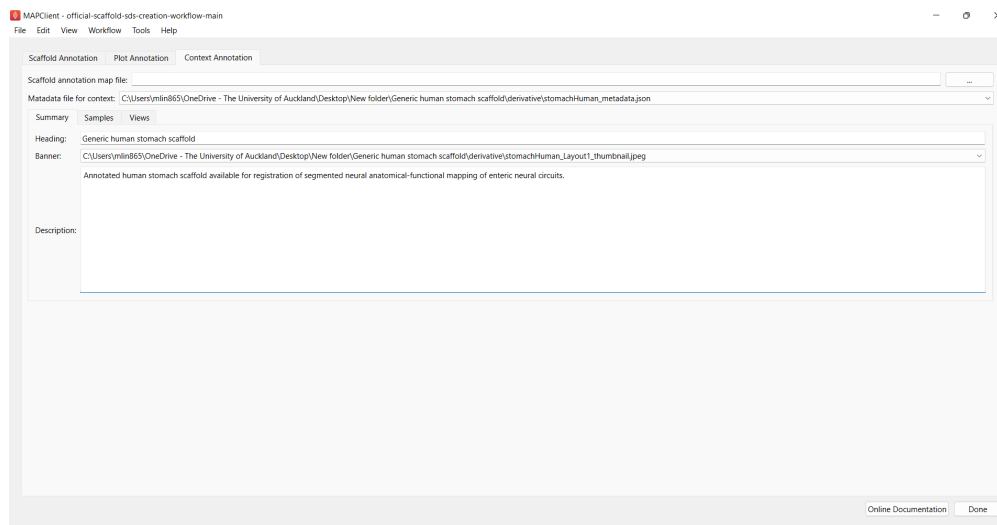


Figure 11. Summary for context card annotation.

Click on the **View** tab next to add views to the context card (Figure 12). Multiple scaffold views can be added to capture specific views. The different view options will be listed on the context card when published and portal users can select a particular view and the scaffold will be displayed in the captured view. This feature is useful for showcasing specific details of the scaffold or its mapped data. In your case, you will only need to add scaffold views if there are specific anatomical features that are useful to highlight so it is

entirely up to you to decide if you want to include views for your generic scaffold. Here we will add a perspective view of the scaffold to document how a view can be added.

Select the **Path** of the view via browsing to the specific view json file that can be found in the derivative folder in the output directory that was set up during the “GenerateSDS” step in **Step 4**. The same goes for the “Thumbnail” field. Lastly, add a description to describe the details of the specific field. Perform the same for other views if you choose to add multiple views of the scaffold.

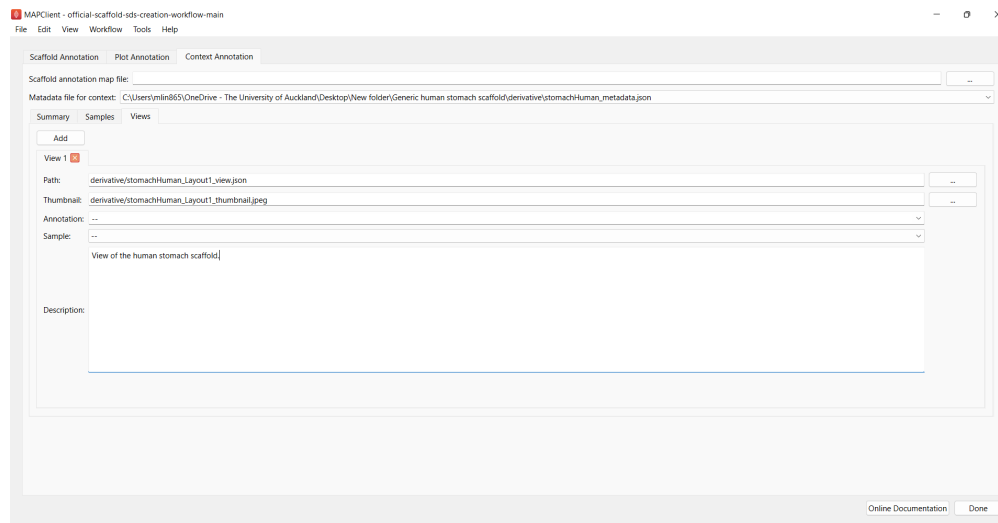


Figure 12. Adding views to the context card.

- 8 After clicking **Done**, go to your generated SDS dataset (output directory set in “GenerateSDS” step in **Step 4**. You should see that the directory contains a folder with the same name as the “Dataset name” set in **Step 4**. There will be the following files:
- derivative folder (manifest.xlsx, webGL files in json, thumbnail file in jpeg, context information file in json, VTK, and STL files)
 - primary folder (see Figure 13 for expected files and folders)
 - dataset_description.xlsx
 - submission.xlsx

Generic human stomach scaffold > primary > Search primary

Name	Status	Date modified	Type	Size
ViewCreation-previous-docs	✓	10/06/2024 9:48 am	File folder	
manifest.xlsx	✓	10/06/2024 12:28 pm	Microsoft Excel W...	6 KB
map-client-workflow.proj	✓	10/06/2024 11:56 am	PROJ File	3 KB
OrganScaffold.conf	✓	10/06/2024 11:56 am	CONF File	1 KB
OrganScaffold-settings.json	✓	10/06/2024 11:56 am	JSON File	13 KB
provenance.json	✓	10/06/2024 11:56 am	JSON File	9 KB
scaffold_info.json	✓	10/06/2024 11:56 am	JSON File	1 KB
STL.conf	✓	10/06/2024 11:56 am	CONF File	1 KB
ThumbnailExport.conf	✓	10/06/2024 11:56 am	CONF File	1 KB
ViewCreation.conf	✓	10/06/2024 11:56 am	CONF File	1 KB
VTK.conf	✓	10/06/2024 11:56 am	CONF File	1 KB
WebGLExport.conf	✓	10/06/2024 11:56 am	CONF File	1 KB

Figure 13. Contents of the primary folder in a generated SDS dataset.

Check that the contents in the generated SDS dataset are correct. For reference, you can use the "Generic human stomach scaffold" dataset as a gold standard for reference as that dataset has been published and meets all curation requirements.

- 9 If the scaffold has not been published before, we will need to create a new PMR workspace and exposure for it. If it is an existing dataset and the exposure has been revised since its last PMR exposure, we will need to update the PMR workspace to include the latest files in the primary folder of the generated SDS folder. Apart from these files, the PMR workspace should also include an updated version (if changed):

- README.rst
- *_thumbnail.jpeg (from derivative folder of generated SDS dataset)

Create a new PMR exposure and retrieve a URL link to the new PMR exposure. You can either update the PMR link in your dataset_description.xlsx file directly or run the saved official scaffold sds creation workflow again to update the exposure link.

- 10 Log in to **Pennsieve**. Create a new dataset if it is a new organ scaffold. If the dataset already exists, go to the existing organ scaffold dataset. Click on **Files** on the left column and delete all files from the dataset. Upload the following files from your generated SDS dataset:

- derivative
- primary
- dataset_description.xlsx
- submission.xlsx

In the "Overview" section, on the top right hand corner of the "Description" section, click on the word "Update". Edit the "Primary vs derivative data" and "Code Availability" sections under "Curator's Notes" to read as follow and save your changes:



Primary vs derivative data: The primary folder contains the mapping tool provenance data file describing the software environment in place when this dataset was created. The primary folder also contains the settings files which in conjunction with the software information in the provenance file, will reproduce the output files stored in the derivative folder. The derivative folder contains JSON files that are used to generate a WebGL visualization of the scaffold on the SPARC portal, as well as the VTK and STL formats of the scaffold.

Code Availability: Scaffold Mapping Tools

Then click on “Update” in “Changelog” to add and save an update describing the changes that have been made to the dataset. For instance,

[Insert date] - Files in primary and derivative folders are replaced by files generated using the sds creation workflow in MAPClient [insert version number]. Dataset description has been updated to the latest SDS format and Curator's Notes have been updated with information about the primary and derivative data and the link in Code Availability has been updated with a link to the Scaffold Mapping Tools. Submission doc is updated with the latest milestone and completion date.

When done, navigate to “Publishing” on the left column and click on the “Request to Publish” button. The dataset will then be reviewed by the curation team. Once curation is complete, you will receive an email confirming that the dataset has been accepted and republished.