

Nov 15, 2023

Version 1

QUINT Workflow Appendix V.1

 In 1 collection

DOI

dx.doi.org/10.17504/protocols.io.36wgq3e6klk5/v1

Michael X. Henderson¹

¹Van Andel Research Institute



Michael Henderson

Van Andel Research Institute

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.36wgq3e6klk5/v1>

Protocol Citation: Michael X. Henderson 2023. QUINT Workflow Appendix. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.36wgq3e6klk5/v1>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working



Created: August 17, 2023

Last Modified: June 01, 2024

Protocol Integer ID: 90666

Keywords: ASAPCRN, quint workflow appendix this protocol, additional details for the modified quint workflow, modified quint workflow, quint workflow, quint workflow appendix, modified quint, quint, protocol

Funders Acknowledgements:

National Institute on Aging

Grant ID: R01-AG077573

Aligning Science Across Parkinson's

Grant ID: ASAP-020616

Abstract

This protocol contains additional details for the modified QUINT workflow.

Attachments



[812-2118.pdf](#)

4.6MB

Troubleshooting

QuPath: Setting Stain Vectors

- 1 Draw a small rectangle over a characteristic DAB or hematoxylin stain.
- 2 On Image tab, double-click the Stain you would like to set. Click Yes on the prompt.
- 3 Give the new settings a unique name (Project ID). Click OK.

Note

***Note:** You could also use "Estimate Stain Vectors."
- 4 This has only saved the parameters for this particular image. To change all images in a project, go to the Workflow tab. This script can also be added to an analysis script and run at the same time.
- 5 "Set color deconvolution stains" should be the last Command. Click on Command, verify that the parameters are correct, and click "Create script."
- 6 In the Script Editor, click Run > Run for Project.
- 7 Move all images from Available to Selected, and click "OK." Modified setting from Aperio Scanner:

Hematoxylin: 0.667, 0.658, 0.349
DAB:0.324, 0.547, 0.772
DAB:0.485, 0.616, 0.620


QuPath: Pixel Classifier Generation

- 8 To create a new pixel classifier for segmentation:

It is easiest if the output segmentation from QuPath is black. Click on the "Annotations" tab,
double click on "Positive" and change the color to black. Click "OK".
- 9 Use the rectangle tool to box an area you can use to test the classifier.

10 Classify > Pixel classification > Create thresholder

11 Example settings are below, but can be optimized by project. Generally higher resolution is better, but slower to run and will create larger files. We generally use "Very high (1.01 um/px)" resolution. Give the classifier a name and click "Save". The classifier will be in your project folder under classifiers and can be moved to another qupath project folder to share with another set of images for consistency.

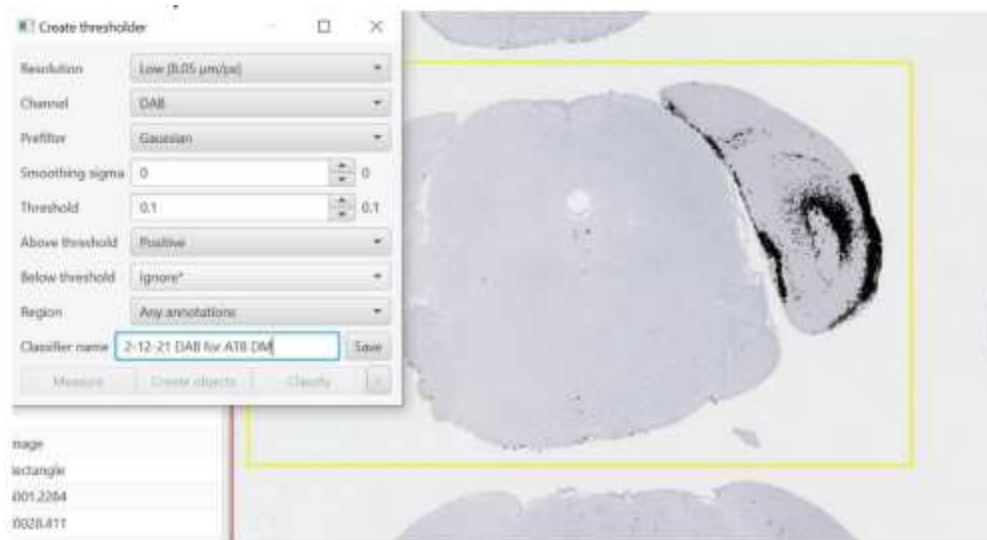

12 Examples of Staining Artifacts:

