

Mar 12, 2019

Serial killer pore investigation: protocol for analyzing fatigue failure pores in additively manufactured Ti6Al4V

DOI

<https://dx.doi.org/10.17504/protocols.io.yebftan>

Muofhe Tshibalanganda¹, Anton du Plessis¹, Stephan le Roux¹

¹CT scanner facility



Muofhe Tshibalanganda

CT scanner facility

Create & collaborate more with a free account

Edit and publish protocols, collaborate in communities, share insights through comments, and track progress with run records.

Create free account

OPEN  ACCESS



DOI: <https://dx.doi.org/10.17504/protocols.io.yebftan>

Protocol Citation: Muofhe Tshibalanganda, Anton du Plessis, Stephan le Roux 2019. Serial killer pore investigation: protocol for analyzing fatigue failure pores in additively manufactured Ti6Al4V. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.yebftan>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited



Protocol status: Working

We use this protocol and it's working

Created: February 21, 2019

Last Modified: March 12, 2019

Protocol Integer ID: 20643

Keywords: dmls ti6al4v specimen, ti6al4v specimen, tensile fatigue specimen, pore investigation, analyzing fatigue failure, ti6al4v, porosity, killer pore, small amount of porosity, eventual fracture surface, mm region near the eventual fracture surface, pores in the vicinity, failure by fatigue test, manufactured metal, fatigue failure, fatigue test, pore, isolated pore, area of failure

Abstract

This protocol outlines the steps followed to determine a potential "killer pore" in the area of failure of DMLS Ti6Al4V specimen.

The methodology can be used more generally for positively identifying and quantifying small isolated pores in additively manufactured metals. In this case the focus was on analysis of a 4 mm region near the eventual fracture surface, to identify all pores in the vicinity of the "killer pore" - the one causing failure by fatigue tests. This protocol was developed for Ti6Al4V with 99.99% density (only small amount of porosity, very small pores only), tensile fatigue specimen round dogbone geometry with 5 mm gauge width.

Guidelines

The samples are scanned before fatigue testing.

Before start

Make sure to scan with the sample number facing up to know which side is up during analysis.

Before Scanning

- 1 Wrap the sample inside a bubble wrap for support. Make sure the marked threaded grip is facing up. This is to know which side is up when analysing.



Sample mounting

- 2

Load the sample on a foam (e.g. florist oasis).



Scanning Stettings

- 3 The sample was scanned using the scanned using Phoenix V|Tome|X L240 (General Electric Sensing and Inspection Technologies Phoenix X-ray, Wunstorff, Germany) microCT system.


	Voltage	160k V
	Current	100 μ A
	X-ray Sensitivity	4
	Voxel size	15 μ m
	Timing	333 m/s
	Images	800
	Averaging	2

Skip

1


Reconstruction

- 4 Reconstruction of the 3D data is done using General Electric Datas 2.2, using specifically beam hardening correction value of 8 - this reduces beam hardening artefacts significantly.

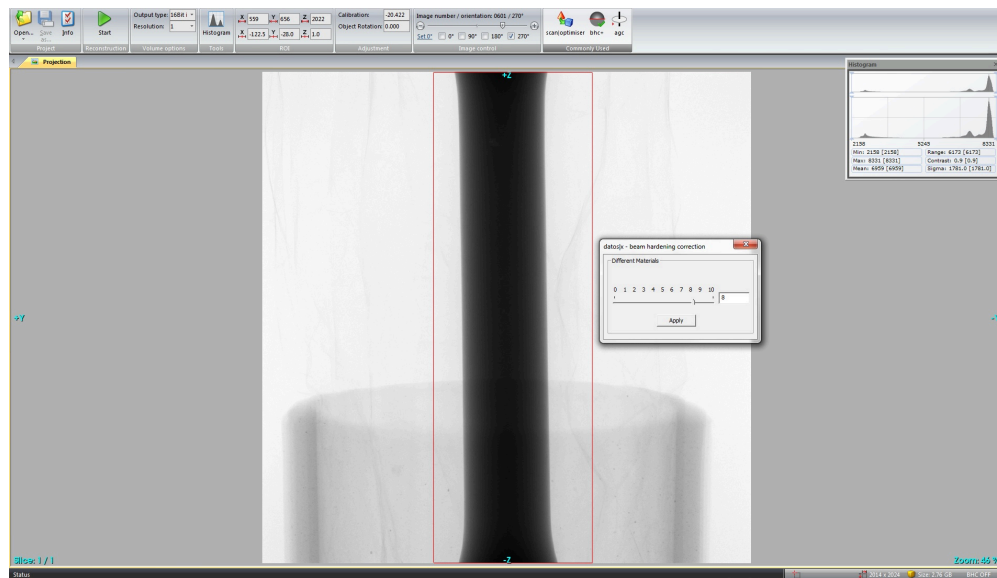


Porosity analysis


- 5 Analysis of 3D data is described here using Volume Graphics VGSTUDIO MAX 3.2 .

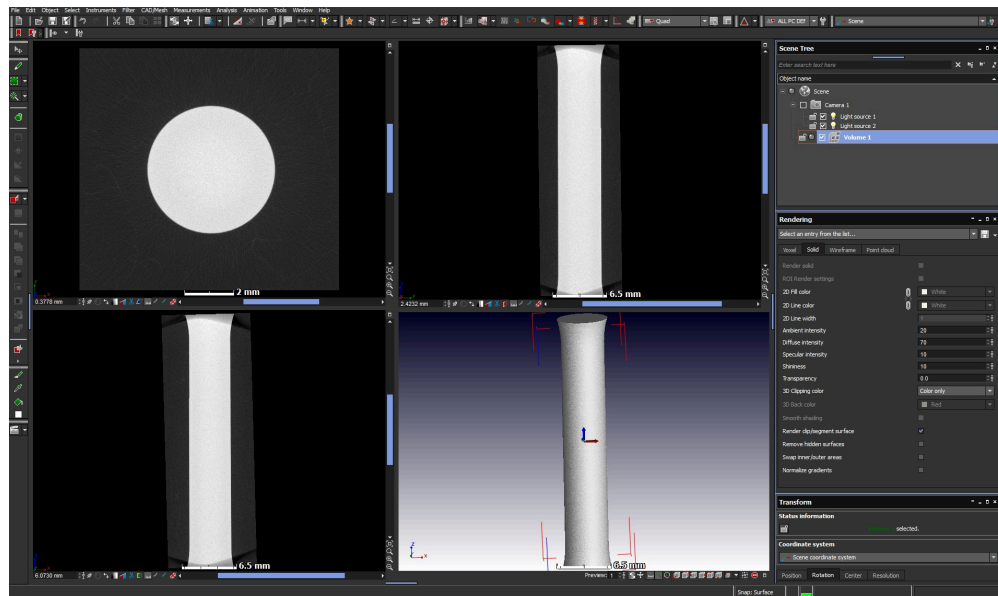


- 6 Register the sample - use simple registration to turn the sample upright in the software windows.

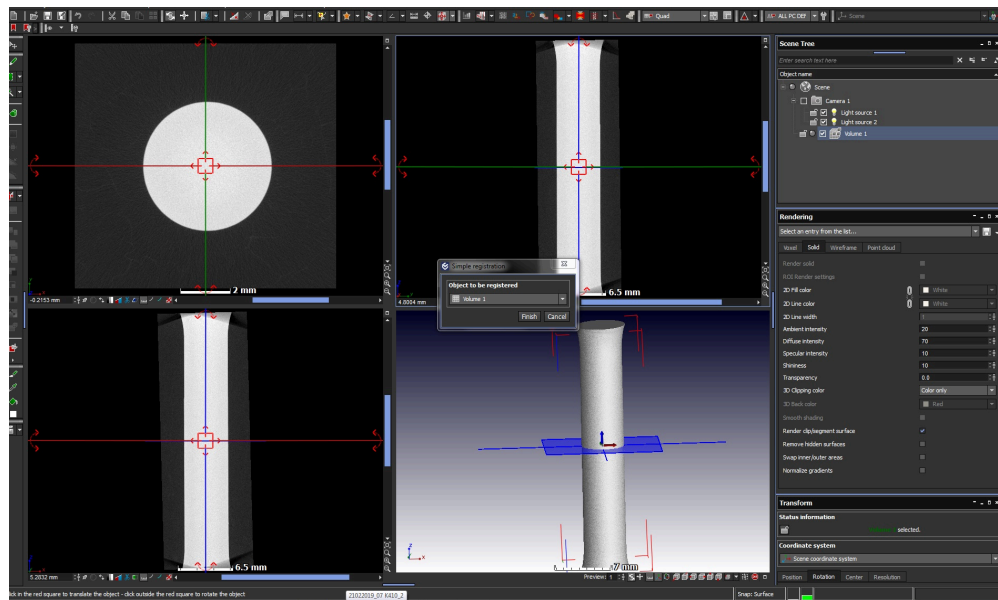


- 7 Use the ellipse Region of interest to crop the sample excluding the area with the cone beam artefact.

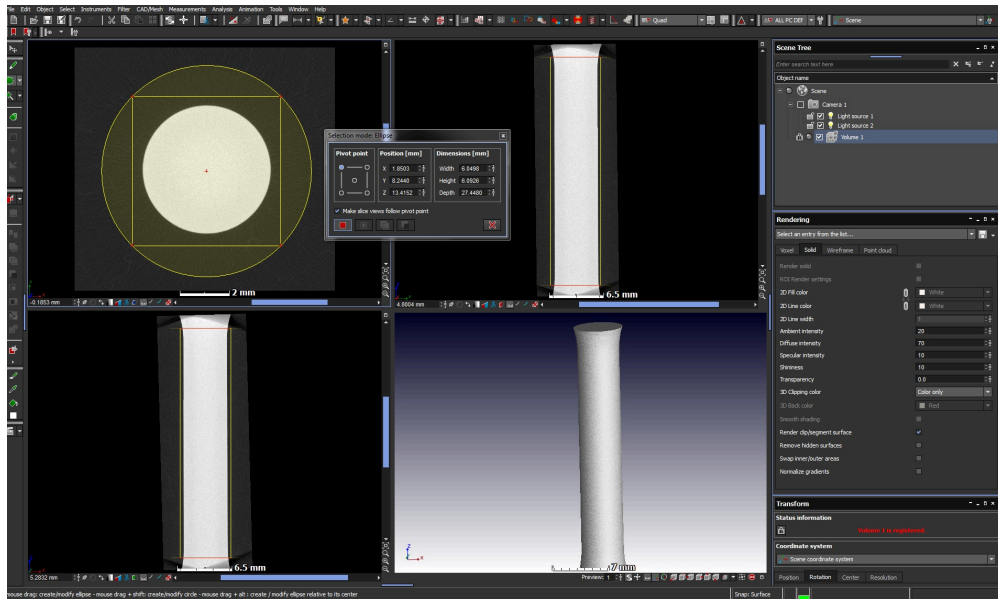


- 8 Extract the region of interest to create new volume.
- 9 De-noising is done using the filters Adaptive Gauss with smoothing factor 0.8, edge threshold=0.1, iteration=1 and create a new volume.

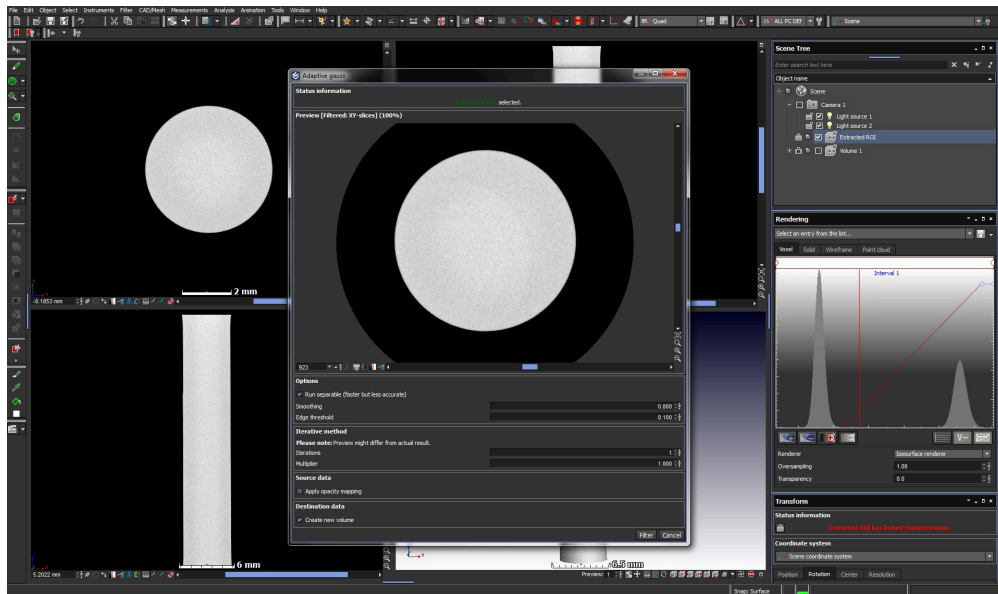


10 Define the material by using the automatic and advanced mode surface determination.

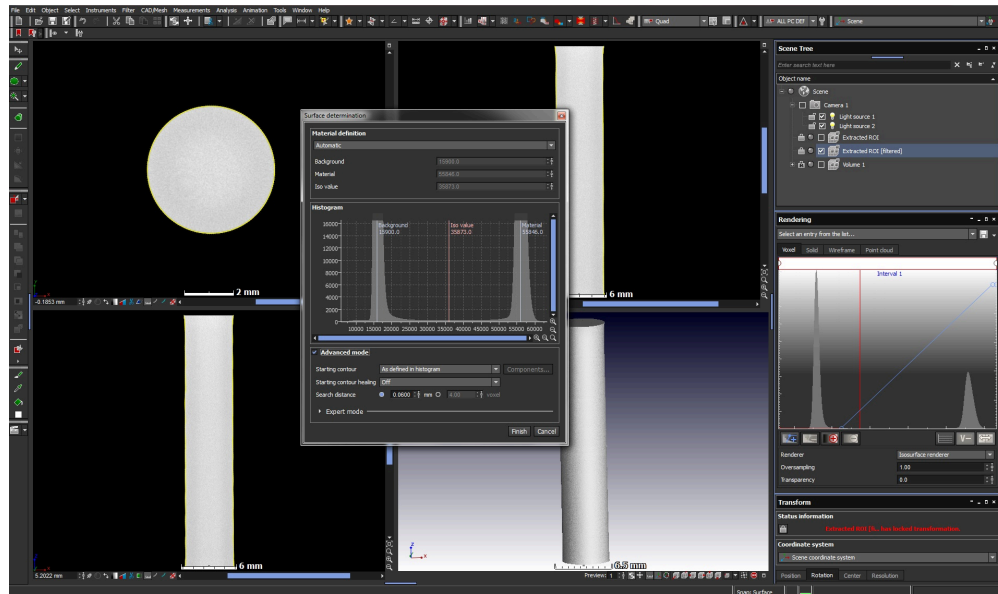


11 In the area of failure crop out 4 mm

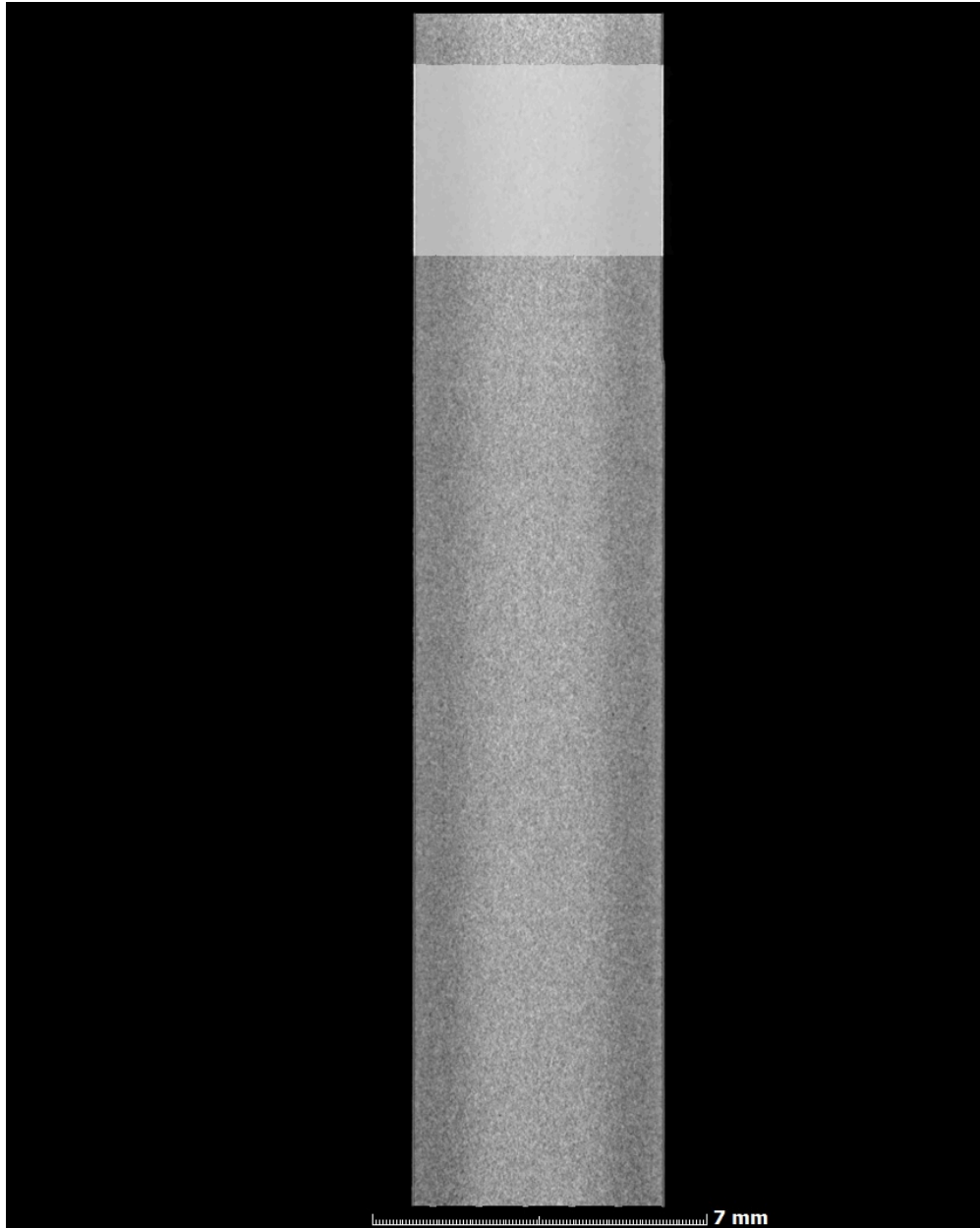



- 12 Denoise using the non local mean filter using the smothing factor 2

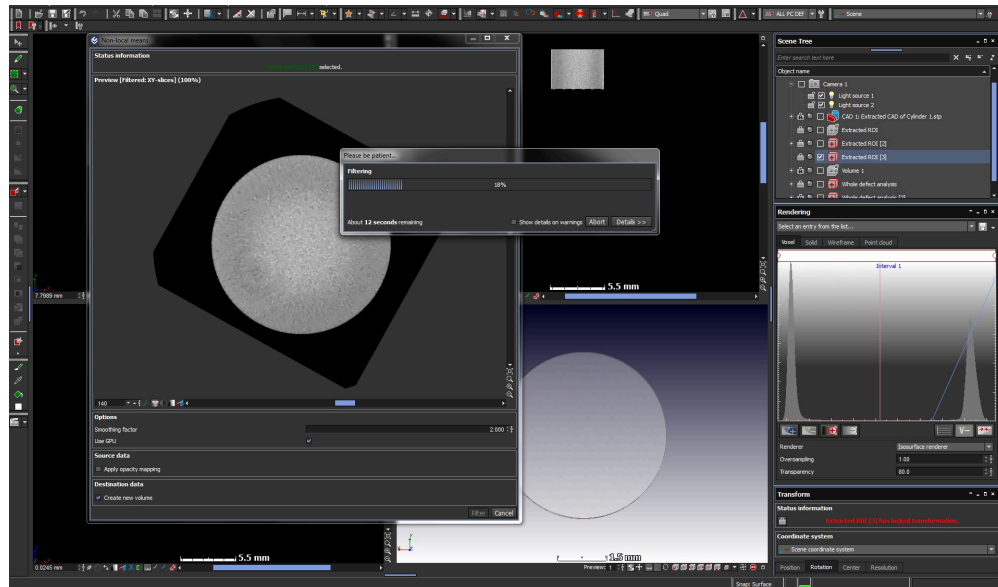

- 13 Define the material by using surface determination. Use the automatic and advanced mode.
- 14 Perform the porosity analysis using
 - Probability 0
 - minimum voxel 8
 - Maximum voxel 400
 - Auto threshold deviation
 - Deviation factor
 - Select the edge distance calculation

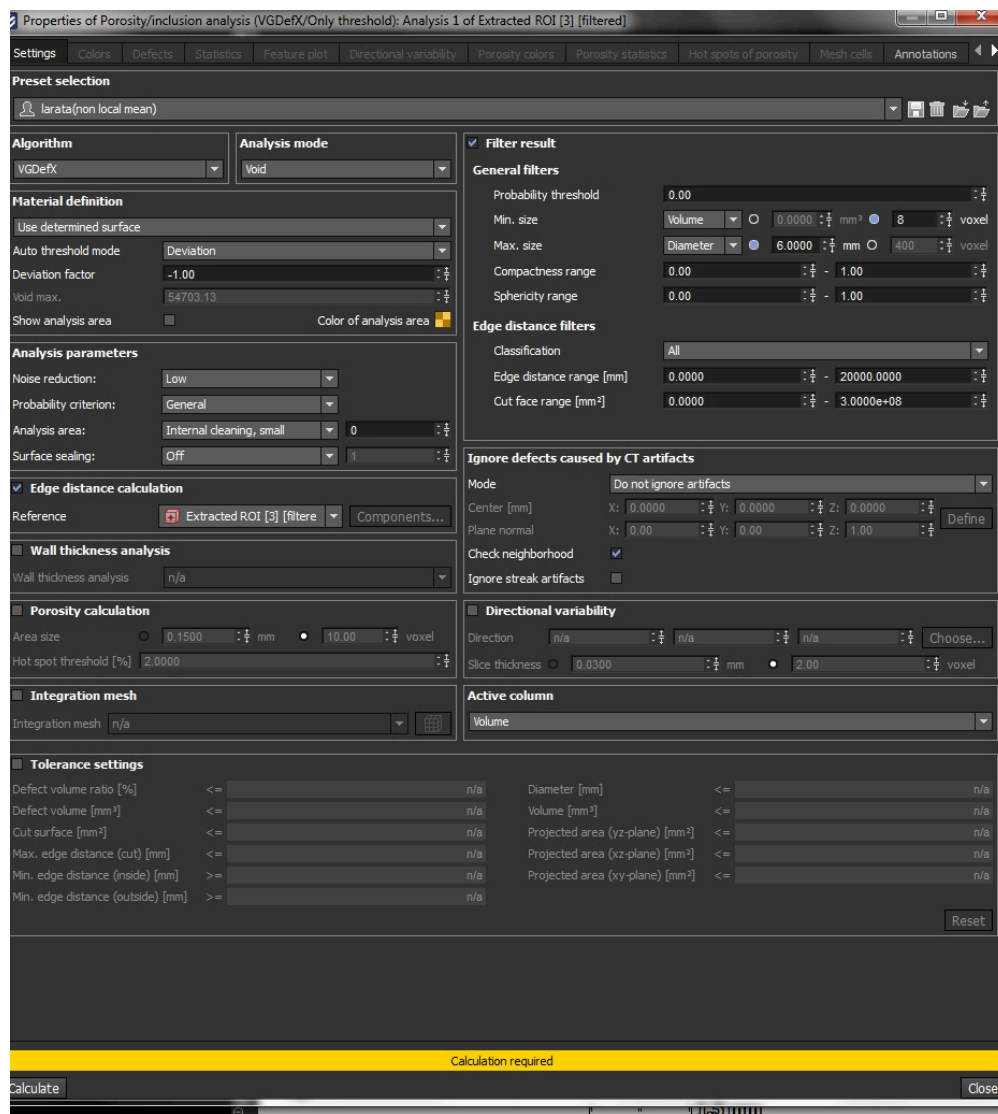

- 15 Sort by minimum edge distance to clearly define the potential crack initiating pores closest to the surface.

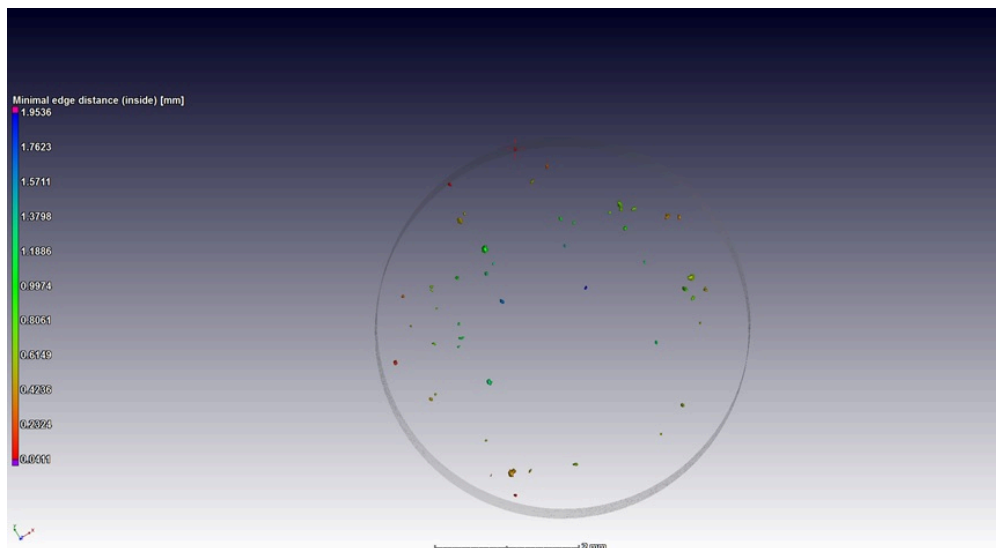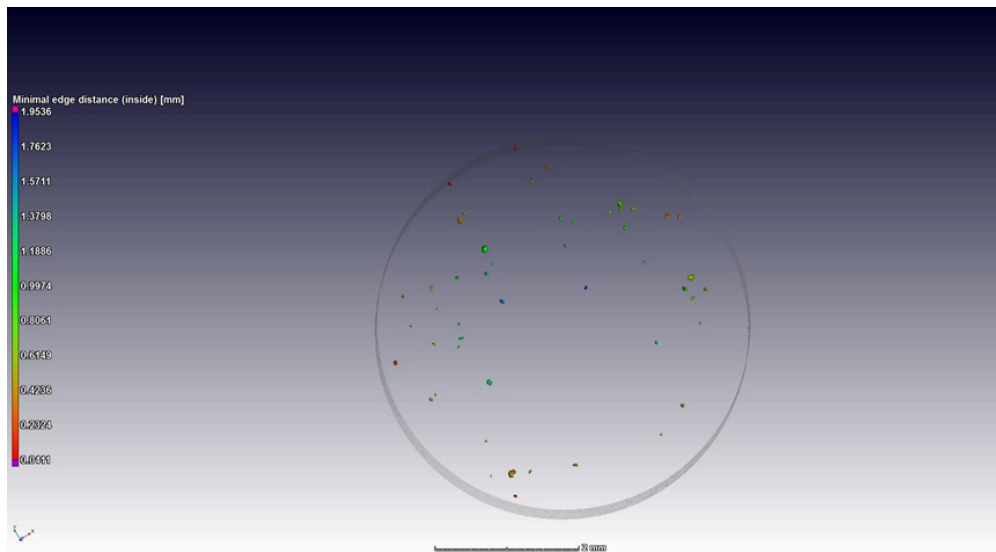

 3D rotation video_1.wmv

 3D rotation video_2.wmv

16 Save the defect analysis as csv file.



- 17 The data is saved including the filtered data set, as a image stack (tiff). These can be viewed manually without special software. The analyzed volume can be opened in free myVGL software.